



Defence Research and  
Development Canada

Recherche et développement  
pour la défense Canada



# **Development of Sampling Methods for Powders and Soil for Detection of Biothreat Agents by Electrochemiluminescence**

H.G. Thompson, R.E. Fulton, and J. Ranches  
DRDC Suffield

Technical Memorandum  
DRDC Suffield TM 2007-172  
June 2007

**DISTRIBUTION STATEMENT A**  
Approved for Public Release  
Distribution Unlimited

Canada

# **Development of Sampling Methods for Powders and Soil for Detection of Biothreat Agents by Electrochemiluminescence**

H.G. Thompson, R.E. Fulton, and J. Ranches  
DRDC Suffield

**Defence R&D Canada – Suffield**

Technical Memorandum  
DRDC Suffield TM 2007-172  
June 2007

*AQ F07-12-13505*

**20071002077**

Principal Author

*Original signed by H.G. Thompson*

---

H.G. Thompson

Approved by

*Original signed by S.J. Armour*

---

S.J. Armour

DHCBDS

Approved for release by

*Original signed by P.A. D'Agostino*

---

P.A. D'Agostino

DRP Chair

© Her Majesty the Queen in Right of Canada, as represented by the Minister of National Defence, 2007

© Sa Majesté la Reine (en droit du Canada), telle que représentée par le ministre de la Défense nationale, 2007

## Abstract

---

The M1M (BioVeris Corporation), an automated electrochemiluminescent (ECL) assay instrument, has been shown to be a technology that is robust and tolerant of samples in a variety of complex matrices; however, methods for sampling of powders and soils in the field for processing on the M1M have not been described. This report describes the development of simple and reproducible methods that may be used for field sampling of powders and soils for subsequent screening and analysis on the M1M. By sampling with polyester swabs using a standard method, it was found that it was possible to collect a reproducible amount (~12 mg) of powder or soil from a variety of surfaces. Using this value, the volume of diluent required to produce a signal <1.2 on the M1M was empirically determined. The powder or soil suspension was filtered and the filtrate assayed on the M1M. The result was a method that could be used consistently to assay powder and soil by ECL M1M with good sensitivity and minimal background signal.

## Résumé

---

Le M1M (BioVeris Corporation), un instrument automatisé de bioessais électrochimiluminescents (ECL) s'avère être une technique robuste et tolérante des échantillons, dans une variété de matrices complexes; on n'avait cependant pas encore décrit les méthodes d'échantillonnage sur le terrain des poudres et des sols devant être traités sur le M1M. Ce rapport décrit la mise au point de méthodes simples et reproductibles pouvant être utilisées pour l'échantillonnage sur le terrain des poudres et sols, ultérieurement criblés et analysés sur le M1M. On a trouvé qu'en utilisant la méthode standard d'échantillonnage avec des porte-coton en polyester, il était possible de recueillir une quantité reproductible (~12 mg) de poudre ou sol provenant d'une variété de surfaces. On a déterminé empiriquement, en utilisant cette valeur, le volume de diluant requis pour produire un signal <1,2 sur le M1M. La poudre ou sol en suspension a été filtré et le filtrat bio-testé sur M1M. On en conclut que la méthode peut être utilisée régulièrement pour bio-tester les poudres et sols sur le M1M ECL avec une bonne sensibilité et un signal de fond minimal.

This page intentionally left blank.

## Executive summary

---

### Development of Sampling Methods for Powders and Soil for Detection of Biothreat Agents by Electrochemiluminescence

H.G. Thompson, R.E. Fulton, and J. Ranches; DRDC Suffield TM 2007-172; Defence R&D Canada – Suffield; June 2007.

**Introduction or background:** The M1M (BioVeris Corporation), an automated electrochemiluminescent (ECL) assay instrument, is currently being considered by several agencies for field identification of biothreat agents. The technology employed has been shown to be robust and tolerant of samples in a variety of complex matrices. However, methods for sampling of powders and soils in the field for processing on the M1M had not been developed previously.

**Results:** A method for collecting reproducible amounts of powder or soil from a variety of surfaces using polyester swabs was devised. Relatively constant sample size permitted a simple, standard method for sample preparation to be developed. The overall protocol allows suspect powders and soil to be assayed by the M1M with good sensitivity and minimal background signal.

**Significance:** The sampling protocol developed here, and the operation of the M1M, are rapid, relatively simple, and easily transferable from a laboratory setting to the field. Their adoption by the Canadian Forces and by first responders would afford a markedly improved ability to detect and identify agents of biological origin.

**Future plans:** The study will be expanded to include other powders and soil types, as well as other biothreat agents.

## Sommaire

---

### Development of Sampling Methods for Powders and Soil for Detection of Biothreat Agents by Electrochemiluminescence

H.G. Thompson, R.E. Fulton, and J. Ranches; DRDC Suffield TM 2007-172; R & D pour la défense Canada – Suffield; juin 2007.

**Introduction ou contexte:** Le MIM (BioVeris Corporation) un instrument automatisé de bioessais électrochimiluminescents (ECL), est actuellement considéré par plusieurs agences pour l'identification, sur le terrain, d'agents représentant une menace biologique. La technologie employée s'avère robuste et tolérante des échantillons dans une variété de matrices complexes. On n'avait cependant pas décrit auparavant les méthodes d'échantillonnage sur le terrain des poudres et des sols devant être traités sur le MIM.

**Résultats:** On a conçu une méthode, visant à recueillir les quantités reproductibles de poudre ou sols provenant d'une variété de surfaces, au moyen de porte-coton en polyester. La taille relativement constante des échantillons a permis de mettre au point une méthode standard de préparation d'échantillons. Le protocole général permet de bio-tester les poudres et sols suspects sur le MIM avec une bonne sensibilité et un signal de fond minimal.

**Portée des résultats:** Le protocole d'échantillonnage mis au point ici ainsi que l'opération du MIM sont rapides et relativement simples et facilement transférables d'un contexte de laboratoire à celui du terrain. L'adoption de ce protocole par les Forces canadiennes et les premiers intervenants permettrait une amélioration marquante de la capacité à détecter et à identifier les agents d'origine biologique.

**Plans futurs:** On prévoit d'étendre cette étude à d'autres types de poudres et sols ainsi qu'à d'autres agents représentant une menace biologique.

# Table of contents

---

Abstract .....	i
Résumé .....	i
Executive summary .....	iii
Sommaire .....	iv
Table of contents .....	v
List of figures .....	vi
List of tables .....	vi
Acknowledgements .....	vii
1 Introduction.....	1
2 Material and methods.....	2
2.1 MIM Analyzer .....	2
2.2 Reagents and Materials.....	2
2.3 Assays.....	3
2.4 Antigens.....	3
2.5 Laboratory Protocol.....	3
2.6 Field Protocol .....	3
2.7 Analysis of Data and Statistics .....	4
3 Results and Discussion .....	5
4 Conclusion .....	11
References .....	12
Annex A Running a Plate on the MIM.....	15
Annex B Powder Swab Protocol .....	17
Annex C Soil Swab Protocol .....	18

## List of figures

---

Figure 1. MIM assay of <i>F. tularensis</i> in powders and soil using the BioVeris suggested matrix:diluent of 1:10.....	6
Figure 2. MIM assay using ricin reagents: increasing ratios of BV-DILUENT™ to matrices.....	7
Figure 3. MIM assay of ricin using the final protocol for the use of swabs to assay powders and soil .....	9

## List of tables

---

Table 1. Amounts of powders on swabs.....	8
Table 2. MIM assay of dairy creamer, flour, and soil from swabs .....	10
Table 3. MIM assays of flour and soil in triplicate.....	10

## Acknowledgements

---

The authors wish to acknowledge Mr. Raphael Cebellero, BioVeris Corp., who provided on-site training in M1M assay performance at DRDC Suffield. Thanks also are extended to Dr. Xiaojiang Dai for her literature search on electrochemiluminescence.

This page intentionally left blank.

# 1 Introduction

---

Incidents involving suspicious powders, especially those involving anthrax and ricin-laced letters in the US postal system [1], have led to substantial effort at DRDC Suffield and elsewhere to develop suitable technologies and sampling procedures for the analysis of such materials.

Some work on the development of methods to detect bacteria in powders has been reported for *Bacillus thuringiensis* and other viable spores [2, 3]. Successful use of a fluorescence-based hand-held assay system, the Rapid Analyte Measurement Platform (RAMP™), has been described for detection of ricin in a variety of common powders [4]. The Association of Analytical Communities (AOAC) International [5] has developed and approved a standard procedure for collecting, packaging, and transport of samples of visible powders suspected of being biological threat agents [6]. A review of current and developing technologies for monitoring biological threat agents includes methods for sample processing of a variety of matrices including powders and soil [7].

Much has been published on the analysis of a variety of substances in soil ranging from inorganic metals [8–10], to pesticides [11], herbicides [12], polycyclic aromatic hydrocarbons [13], and other pollutants [14, 15]. A common thread in most of these publications is the necessity for clean-up step(s) prior to assay, to prevent interfering components of the soil from generating false positives [8–10]. For biological agents, the usual method of detection is to dilute the soil sample and isolate the agent in suitable media [16]. While this technique may be feasible for identification of live agents, it is not suitable for non-living organisms or toxins. Successful use of an artificial marine worm gut fluid has been described for the extraction of *Bacillus thuringiensis* from soil [17]; however, little analytical data is available on the extraction of biothreat agents from soil prior to analysis. In-house experimentation with soil clean-up for subsequent analysis by ELISA has resulted in little success in removing background interferents\*.

Electrochemiluminescence (ECL) is a relatively new technology for conducting immunoassays with improved assay performance [18, 19]. The advantages of ECL over enzyme-linked immunosorbant assay (ELISA) and many other immunoassays include low detection limits (high sensitivity) due to high luminescent signal-to-noise ratios [20, 21], large dynamic range of analyte detection of approximately five orders of magnitude [21, 22], reduced labour intensity, and short assay time (16 minutes, including incubation).

ECL assays have been reported for a number of biological threat agents including *Staphylococcus enterotoxin B*, ricin, *Vibrio cholerae*, Botulinum toxin, and *Bacillus anthracis* [21, 23, 24]. ECL detection of biological agents in both food and environmental samples has also been studied [25].

The M1M (BioVeris Corp., Gaithersburg, MD), an automated electrochemiluminescent (ECL) assay instrument, is currently being considered by several agencies for field identification of biothreat agents. The United States Critical Reagents Program (CRP) has developed ECL M1M assays for nine high priority biothreat agents and these are available, off-the-shelf, to US and allied military forces and first responders. ECL M1M technology is used in defence laboratories in the United States, United Kingdom, and Canada; and has been successfully fielded by several

---

\* Unpublished results.

Theatre Labs of the US Army and Air Force. The MIM technology has been shown to be robust and tolerant of samples in a variety of complex matrices; however, methods for sampling of powders and soils in the field for processing on the MIM have not been described.

This report describes the development of simple and reproducible methods that may be used for field sampling of powders and soils for subsequent screening and analysis on the MIM. To acquire the sample, polyester swabs were chosen, as they were deemed the best option for dealing with powders and soils on a variety of surfaces. A procedure for “picking up” reproducible amounts of powder or soil was devised and standardized. An average amount of powder or soil collected on the swab when using the standard method was determined to be 12 mg. Using this value, the volume of diluent required to produce a signal <1.2 on the MIM was empirically determined. The powder or soil suspension was filtered and the filtrate assayed on the MIM. The result was a method that could be used consistently to assay powder and soil by ECL MIM with good sensitivity and minimal background signal.

## **2 Material and methods**

---

### **2.1 M1M Analyzer**

An M-SERIES<sup>R</sup> MIMR ECL Analyzer was purchased from BioVeris Corp. (Gaithersburg, MD) and operated in the MIM mode for these experiments.

The BioVeris M-SERIES<sup>R</sup> MIM Analyzer is an automated single channel system based on ECL technology. The components of the ECL immunoassay are biotinylated capture antibody bound to streptavidin-coated magnetic beads; the antibody-bead complex binds antigen and is drawn into the ECL-analyzer flow cell; detector antibody, labelled with ruthenium-trisbipyridal, for the emission of light when electrochemically stimulated; and a precursor molecule, tripropylamine, which is activated on an electrode surface, resulting in an electron transfer reaction that initiates excitation of the ruthenium-trisbipyridal, resulting in the emission of light at 620 nm. This excitation/emission cycle repeats many times, resulting in increased signal from a single binding event, resulting in increased sensitivity [18].

The complete steps for running an assay on the MIM analyzer were abbreviated from the BioVeris Operation Manual, and are presented in Annex A.

### **2.2 Reagents and Materials**

The following M-SERIES<sup>R</sup> reagents were purchased from BioVeris Corp: BV-CLEAN<sup>TM</sup> Plus Solution, BV-GLO<sup>TM</sup> Plus Solution, BV-STORE<sup>TM</sup> Solution, BV-DILUENT<sup>TM</sup> Solution, BV-SANITIZE<sup>TM</sup> Solution, M-SERIES<sup>R</sup> Positive Calibrator, and M-SERIES<sup>R</sup> Negative Calibrator.

BV Sample Preparation Filters, 60 µm, catalogue number 110117, were purchased from BioVeris Corp. Dacron<sup>R</sup> Polyester Tipped Swabs, Fisherbrand<sup>R</sup>, catalogue number 14-959-90, were purchased from Fisher Scientific (Ottawa, ON).

The following powders were purchased from the US CRP: baking soda, lot number DPG1DEC03; dairy creamer, lot number DPG13MAY04; and flour, lot number DPGLJ24MAY05. In addition, flour (white, enriched, all purpose, Safeway brand) and dairy creamer (Coffee Mate<sup>R</sup>, Carnation, Nestlé) were purchased from Canada Safeway Ltd. (Medicine Hat, AB). The soil was from the Suffield Experimental Proving Ground (Suffield EPG) and had been characterized by Alberta Environmental Centre (Mr. P. Yeung, Vegreville, Alberta) as sandy clay loam.

## 2.3 Assays

MINITube assays for *Francisella tularensis*, ricin, and Venezuelan equine encephalitis virus (VEEV) were purchased from the US CRP. MINITube assays incorporate lyophilized reagents required to perform the MIM assays, i.e. biotinylated anti-analyte antibody bound to streptavidin-coated magnetic beads and ruthenium-labelled anti-analyte antibody. Positive control medium (PCM) and positive control low (PCL) antigens for these agents (inactivated) were also purchased from the US CRP.

## 2.4 Antigens

The following <sup>60</sup>Co-irradiated antigens were purchased from the US CRP: *Francisella tularensis* (Schu 4), lot number AGD0000541, stock concentration  $3.1 \times 10^9$  cfu/mL; and VEEV (TC83), lot number AGD0000109, stock concentration  $1.1 \times 10^{10}$  pfu/mL. Ricin (RCA<sub>60</sub>), catalogue no. L8508, lot no. 63H4026, was purchased from Sigma-Aldrich Canada (Oakville, ON). Prior to use, the stock ricin was diluted in phosphate buffered saline (PBS), pH 7.4, to a working stock concentration of 390 µg/mL and stored at 4°C.

## 2.5 Laboratory Protocol

A standard procedure recommended by BioVeris, in which powder and soil were prepared at concentrations of 100 mg/mL in BV-DILUENT<sup>TM</sup> was initially evaluated. Baking soda (2 x 100 mg), CRP flour (2 x 100 mg), and soil (2 x 100 mg) were weighed into each of six 5 mL polystyrene tubes (Fisher Scientific). BV-DILUENT<sup>TM</sup> (1 mL) was added to each tube. To one tube each of the baking soda, flour, and soil mixtures was added 1.61 µL *F. tularensis* stock to yield final concentrations of  $5 \times 10^6$  cfu/mL. The two baking soda solutions (*F. tularensis*-spiked and unspiked) were assayed without further sample preparation, according to the general MIM procedure described in Annex A. The flour and soil mixtures (spiked and unspiked) were filtered through BV Sample Preparation Filter tubes and aliquots of the filtrates assayed, as above.

## 2.6 Field Protocol

Samples of CRP dairy creamer (4 x 200 mg), CRP flour (4 x 200 mg), and soil (4 x 1g) were weighed into each of twelve 5 mL polystyrene tubes (Fisher Scientific). Ricin was prepared in nanopure water at a concentration of 833 ng/mL for use with powder (creamers and flour) and 2.08 µg/mL for use with soil. One mL ricin solution at 833 ng/mL was added to each of the eight tubes containing the powders. Two mL ricin solution at 2.08 µg/mL was added to each of the four tubes

containing soil. All tubes were vortexed for 15 seconds and their contents poured into individual plastic weigh boats (1 inch square) (Fisher Scientific). The weigh boats containing the ricin-spiked powders and soil were left to air dry in a Class II Biosafety Cabinet for 72 hours. Each dried material was crushed in a mortar with pestle and placed in a microfuge tube for later use.

In order to determine if a consistent amount of powder or soil could be picked up by a swab from a small sample of matrix, the following protocol was followed. A swab was wet with BV-DILUENT™ and then pressed against the sides of the tube containing the buffer to remove excess liquid from the swab. The swab was then weighed on a five decimal place balance. The wet swab was rolled in the powder or soil to acquire maximum sample and then tapped vigorously to remove excess matrix, then weighed again. Five replicates of each of dairy creamer, flour, and soil were performed by one individual, while five replicates of flour were performed by a second individual who followed a written protocol and was not further coached as to procedure.

Increasing amounts of BV-DILUENT™ were used in assaying the powders and soils coated on swabs until all unspiked matrices yielded signals that were well below the S/B cut-off of 1.2, while positive signals from spiked matrices were maintained. The final assay protocol was verified for reproducibility on the MIM by evaluating triplicate samples of both unspiked and ricin-spiked flour and soil. A shortage of CRP ricin MINITubes precluded the inclusion of dairy creamer in these final assays. Further verification of optimum matrix to diluent ratios was obtained by using VEEV MINITubes and positive control antigens with flour from a different source (Safeway brand) and soil samples. In these experiments, performed in triplicate, VEEV, was manually added to swabbed flour in 1mL BV-DILUENT™ and swabbed soil in 2 mL BV-DILUENT™, to yield a VEEV concentration of  $1 \times 10^5$  pfu/mL. The final protocols for use in powder and soil sampling for subsequent analysis on the MIM are attached in Annexes B and C, respectively.

## 2.7 Analysis of Data and Statistics

Negative (no antigen) controls were tested in replicates of four (two prior to matrix assay and two at conclusion of matrix assay). Positive control antigens (PCL and PCM) were tested in replicates of two tubes/assay. Unless otherwise noted, all other samples were tested singly.

ECL readings were considered positive if the emission reading was significantly greater than the background emission reading. Statistical significance, pre-programmed into the MIM, was established at 20% above background (no antigen control). Therefore a signal was considered positive if the value of the sample signal to background signal (S/B) was 1.2 or greater. An automated alarm sounded if any one of the signals had a S/B greater than 1.2.

Data was manually exported from the MIM to Microsoft Excel 7.0 and was graphed using Grapher 6.0 (Golden Software Inc., Golden, CO).

### 3 Results and Discussion

---

Soil has proven problematic as a sample matrix in immunoassays, particularly ELISAs, due to the presence of background interferents generating false positive signals. Non-specific positive reactions were observed when aqueous solutions of soil (clay, loam, or sand) were used to challenge agent-specific ELISAs as potential assay interferents [26]. ELISAs of *Bacillus anthracis*, *Francisella tularensis*, and ricin performed on soil that had been sonicated and treated with a variety of extracting media including PBS, PBS/2% BSA, PBS/2% BSA/0.1 % Tween 20, PBS/0.5% BSA/0.1% SDS, and PBS/0.5% BSA/0.1% sodium taurocholate [17] all produced high backgrounds\*. The only success in detecting ricin in soil by ELISA was by the use of solid phase extraction methods in the treatment of soil collected from a remedial waste site at DRDC Ottawa. However, this soil clean-up method was time-consuming and complicated and recovery of the ricin was low (15–20%)\*.

The results for the protocol recommended by BioVeris Corp. (100mg powder/1 mL BV diluent) are graphically represented in Figure 1. MIM assays for baking soda, flour, and soil, both unspiked and spiked with *F. tularensis*, were all positive. BioVeris Corp. recommends a ratio of sample matrix to diluent of 1:10 (1 mg/10  $\mu$ L = 100 mg/mL) for the analysis of powders and soils. However, BioVeris Corp. uses a Signal/Background (S/B) cut-off of 2.0 for their MIM assays, while the CRP MINITube kits effect a software adjustment of S/B = 1.2. Therefore, it was not unexpected that with the more stringent S/B cut-off being set at 1.2, false positive signals occurred. The observation that false positives occurred, however, indicated that more diluent would need to be added to decrease the background signal.

Depletion of the *F. tularensis* MINITubes early in the study dictated that MINITubes for a different agent be used for subsequent experiments; hence ricin MINITubes and antigen were selected for ensuing studies. In order to estimate the approximate ratio of powder and soil to diluent necessary to avoid false positives, a series of ratios of creamer, flour, and soil (unspiked) to diluent were assessed on the MIM: 100 mg/1 mL; 100 mg/2 mL; 100 mg/3 mL; 100 mg/4 mL; 100 mg/5 mL. The results of this experiment are graphically represented in Figure 2. Both powders (dairy creamer and flour) produced S/B signals that were below the S/B cut-off of 1.2. The soil, however, produced signals that were above the S/B 1.2 cut-off (false positives) for all five test samples examined, thus indicating that additional diluent was warranted. From this experiment, it was deemed that, while the powders could be used at the powder to diluent ratio tested, the soil samples would require that additional diluent be added.

The results of experiments to determine the average amount of material coated on a swab when following a standardized swabbing procedure are presented in Table 1. From the data, the average amount of powder “picked up” by a swabbing procedure was 12 mg. This weight represents the mid range between the median (10.9 mg) and the mean (12.89 mg) amount of powder deposited on a swab by two different individuals following the same swabbing procedure. An estimate of the amount of soil deposited on a swab when the same standard sampling procedure was followed was impossible to determine, as wet swabs coated with soil actually weighed less than the wet swabs themselves. This could have been due to the wicking of BV-

---

\* Unpublished results.

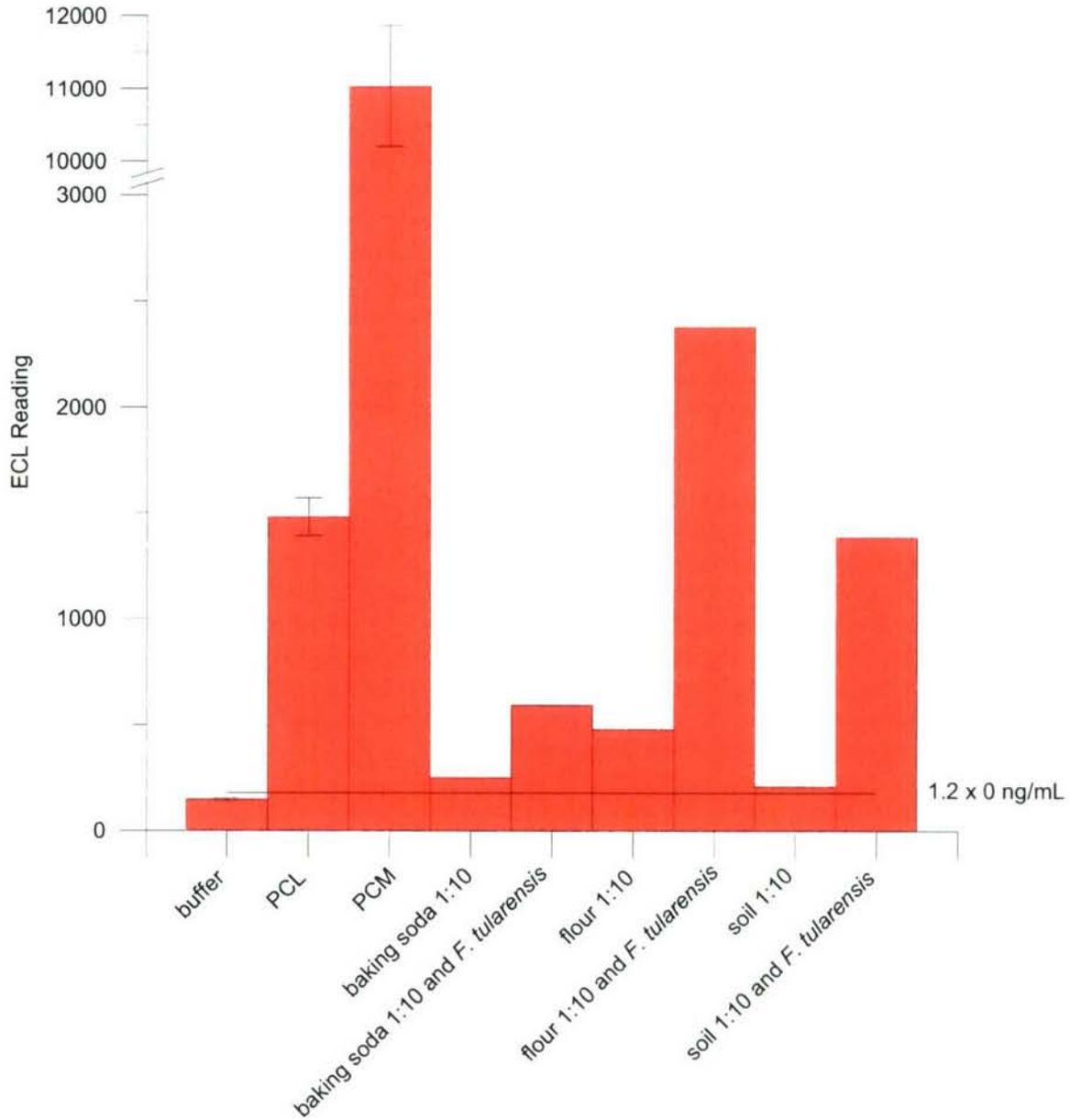


Figure 1. MIM assay of *F. tularensis* in powders and soil using the BioVeris suggested matrix:diluent of 1:10

DILUENT™ from the swab to the dry soil, producing localized muddy sections of soil in the weigh boats. Further experiments are required to estimate the average weight of soil “picked up” by a wet swab, but for the purposes of this study, the average weight for powders (12 mg) was also used for soil ratio calculations.

The next step in this study was to transition procedures from the laboratory method of weighing soils and powder to a method using swabs that would be suitable for use by operators in the field.

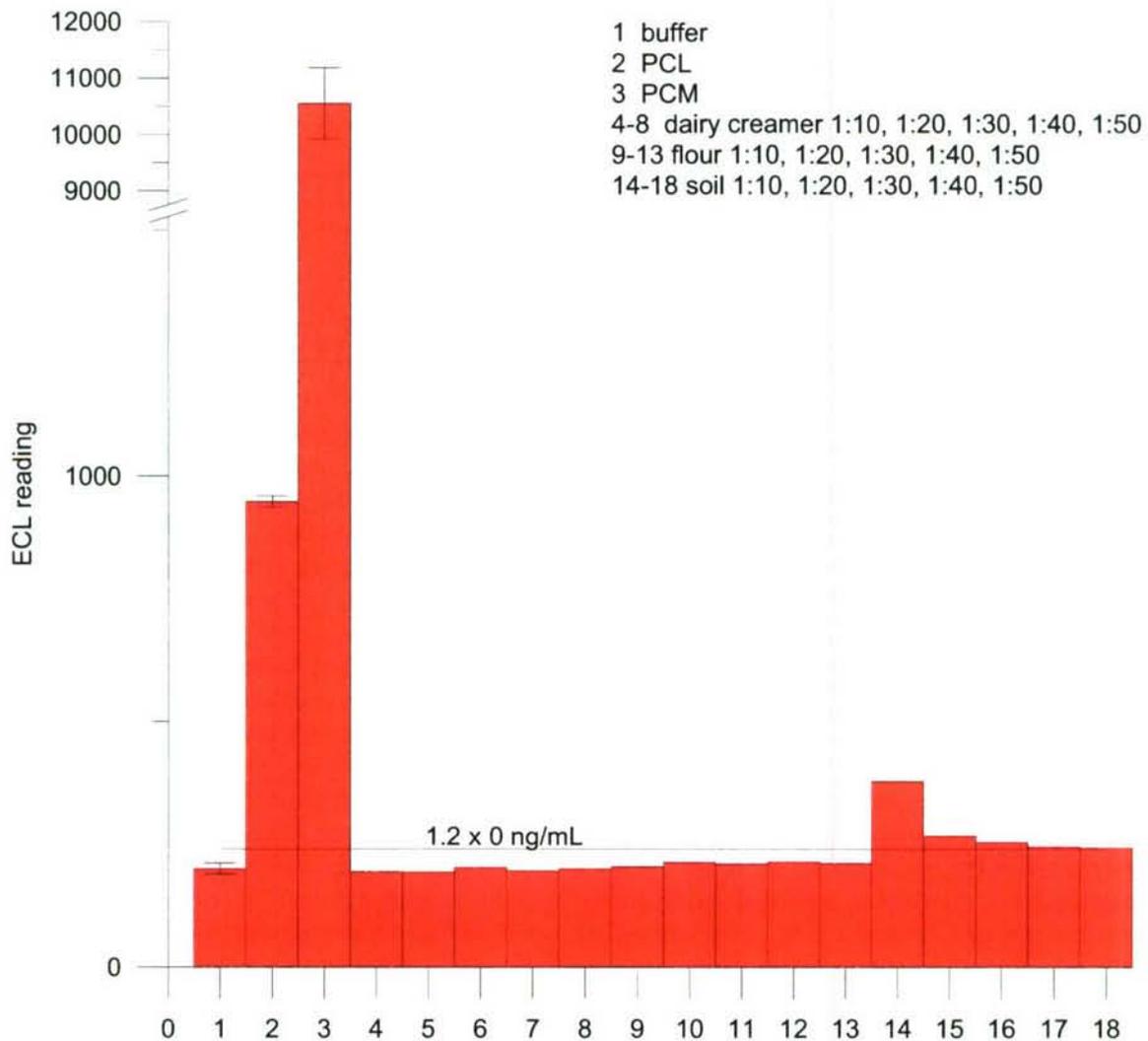


Figure 2. MIM assay using ricin reagents: increasing ratios of BV-DILUENT™ to matrices

Based on the results obtained using weighed amounts of unspiked powder and soil, ratios of ricin-spiked matrix to diluent starting at 1:20 (12 mg/240  $\mu$ L) were used for creamer and flour. Although a ratio of 1 part matrix to 10 parts diluent had not resulted in a false positive when using weighed amounts of matrix (Figure 2), this ratio had resulted in the generation of a false positive in the initial BioVeris protocol experiment using *F. tularensis* reagents (Figure 1). Thus, ratios of 1:20 (12 mg matrix/240  $\mu$ L diluent), 1:30 (12 mg matrix/360  $\mu$ L diluent), and 1:42 (12 mg matrix/ 500  $\mu$ L diluent) were tried for the creamer and flour swabs. No false positives were observed for the dairy creamer at any of the ratios examined (data not shown). Complications in the form of plugged filters were encountered in the filtration of the flour samples at lower matrix to diluent ratios (1:20 and 1:30); therefore additional BV-DILUENT™ was required to be added to the Sample Prep Filter tubes to generate sufficient flour filtrate (50  $\mu$ L) for analysis. In spite of increasing the matrix to diluent ratio to 1:42, the signal for unspiked flour was above the S/B cut-off of 1.2, indicating that even further diluent was necessary.

Table 1. Amounts of powders on swabs

Powder	Individual 1			Individual 2		
	Average weight (mg)	Standard Deviation	CV (%)	Average weight (mg)	Standard Deviation	CV (%)
creamer	15.7	3.4	21	-	-	-
flour	12.5	4.0	32	10.8	2.3	21

When soil was used as matrix, all the soil to diluent ratios examined (1:10–1:50) resulted in false positives (Figure 2). Therefore, a starting ratio of one part soil to 80 parts diluent (1:80), 1:100, and 1:125 were tried. For all of these ratios of soil to diluent, no false positives were generated (data not shown). However, the S/B ratios for all unspiked soil assays were very close to the 1.2 assay cut-off i.e., S/B=1.177, 1.198, and 1.196 for matrix to diluent ratios of 1:80, 1:100, and 1:125, respectively. Hence, the addition of further diluent was considered prudent.

A final protocol for sampling of ricin-spiked powders and soils with polyester swabs was arrived at by selecting a matrix to diluent ratio at which all unspiked matrices produced signals that were acceptably below the 20% above background cut-off (S/B=1.2), while still giving a positive signal for ricin-spiked matrices. Thus, a matrix to diluent ratio of 1:83 (12 mg/1 mL) was selected for powders and 1:167 (12 mg/2 mL) was selected for soils. Results indicated that all unspiked matrices produced signals acceptably below S/B=1.2, while all ricin-spiked matrices were positive (Figure 3; Table 2). Experiments to confirm these findings were conducted using unspiked and ricin-spiked flour and soil matrices in triplicate assays (Table 3). All unspiked matrices resulted in signals well below the 1.2 cut-off, while the ricin-spiked matrices generated positive results. The creamer was not assayed as there were insufficient ricin MINITubes to complete the desired replicates. However, dairy creamer had not posed a problem with false positives from the beginning of the study, and this omission was considered acceptable.

Further verification of the optimum matrix to diluent ratio was obtained using VEEV MINITubes, antigens, and matrices more commonly used at DRDC Suffield. For this verification, replicates of three assays of each of flour (Safeway brand), coffee creamer (Nestlé brand), and soil, both unspiked and spiked, with VEEV ( $5 \times 10^5$  pfu/mL) were tested. Results indicated that no false positives were generated in matrix only assays. Although VEEV-spiked matrices also did not yield positive results, this was attributed to the later finding that VEEV at  $5 \times 10^5$  pfu/mL was considerably below the LOD of the MIM (data not shown). Due to time constraints imposed by the deployment schedule and depletion of available VEEV MINITubes, it was not possible to repeat these assays using spiked VEEV soil at higher concentration. In follow-on experiments, the LOD of VEEV on the MIM was determined to be approximately  $5 \times 10^7$  pfu/mL (data not shown); hence a concentration of  $1 \times 10^9$  pfu/mL was recommended for future use of VEEV.

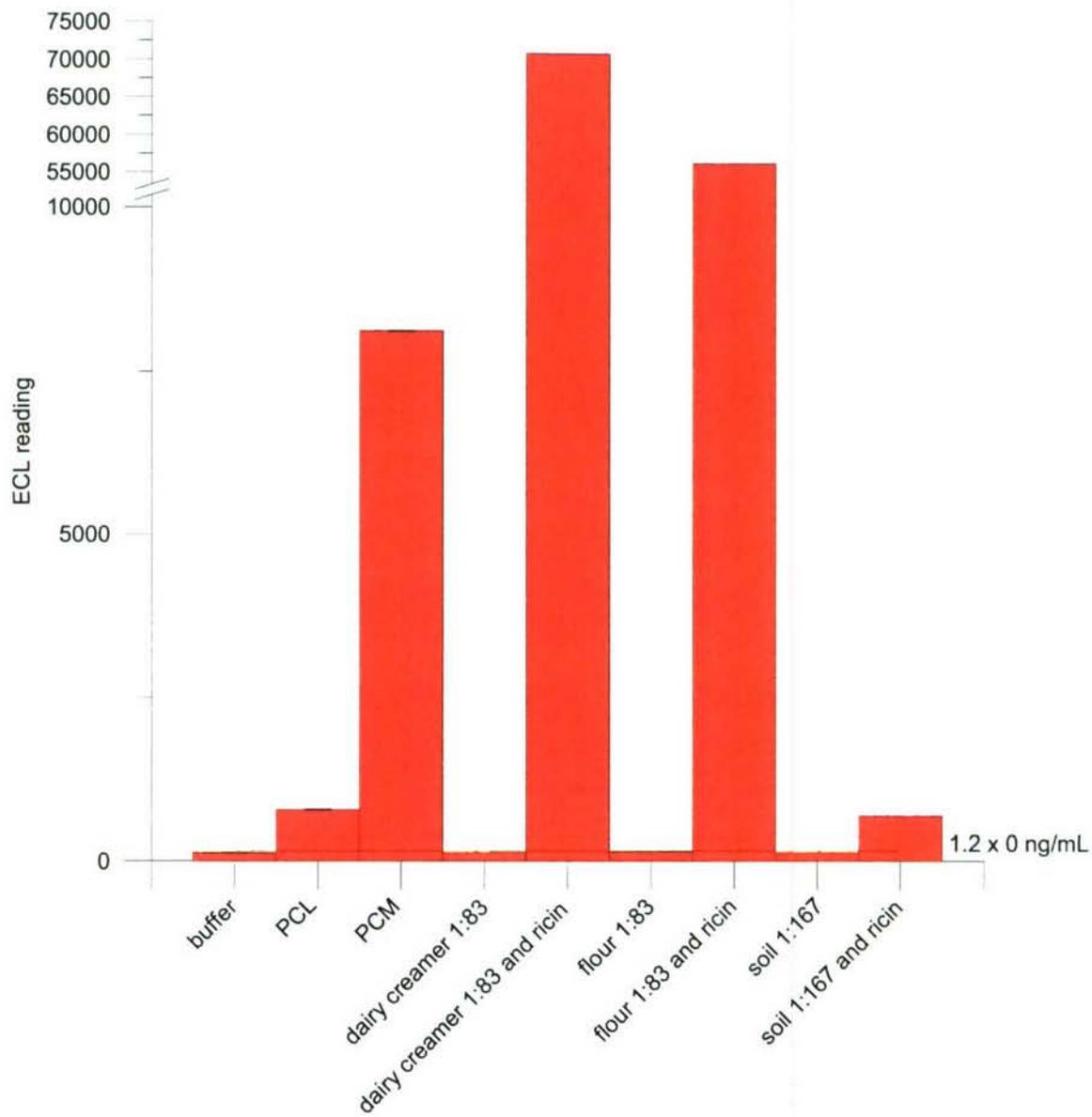


Figure 3. MIM assay of ricin using the final protocol for the use of swabs to assay powders and soil

*Table 2. MIM assay of dairy creamer, flour, and soil from swabs*

<b>Sample</b>	<b>Signal</b>	<b>Standard Deviation</b>	<b>Signal/Background</b>
buffer	132	3.9	-
PCL	780	7.9	5.909
PCM	8123	12.6	61.538
creamer 1:83	141		1.068
creamer and ricin 1:83	70760		536.061
flour 1:83	147		1.114
flour and ricin 1:83	56046		424.591
soil 1:167	141		1.068
soil and ricin 1:167	690		5.227

*Table 3. MIM assays of flour and soil in triplicate*

<b>Sample</b>	<b>Signal</b>	<b>Standard Deviation</b>	<b>Signal/Background</b>
buffer	128	1.8	-
PCL	984	11.1	7.688
PCM	10423	271	81.430
flour 1:83	135	0.58	1.055
flour and ricin 1:83	66447	10327	519.117
soil 1:167	140	7.8	1.094
soil and ricin 1:167	435	106.4	3.399

## 4 Conclusion

---

A procedure for “picking up” reproducible amounts of powder and soil on swabs for subsequent ECL analysis on the MIM platform was devised and standardized for use in the field. The ECL MIM technology is sensitive in detection of biological agents in a variety of complex matrices, such as soil and powders, likely to be encountered by the CF and First Responders. Both the sampling protocol developed here and the operation of the MIM are rapid and relatively simple and easily transferable from a laboratory setting to the field. The study will be expanded to include other powders and soil types as well as other biothreat agents.

## References

---

- [1] Frieden, T. and Cratty, C. (2003-10-23). Toxin Ricin Found at S. Carolina Postal Facility. (Online) CNN. <http://www.cnn.com/2003/US/10/22/ricin.letter> (2004-09-07).
- [2] Min, J.; Lee, J.; and Deininger, R.A. (2006). Simple and Rapid Method for Detection of Bacterial Spores in Powder Useful for First Responders. *J. Environ. Health*, 67(8), 34–37, 44, 46.
- [3] Lee, J. and Deininger, R.A. (2004). A Rapid Screening Method for the Detection of Viable Spores in Powder using Bioluminescence. *Luminescence*, 19(4), 209–211.
- [4] Fulton, R.E. and Thompson, H.G. (2007). Fluorogenic Hand-held Immunoassay for the Identification of Ricin: Rapid Analyte Measurement Platform. *J. Immunoassay and Immunochem.*, in press.
- [5] Anon. (n.d.). AOAC INTERNATIONAL. (Online) AOAC International. <http://www.aoac.org/about/aoac.htm> (2004-11-12).
- [6] Anon. (2006-07-08). Standard Developed For Collection Of Suspicious Powders. (Online) ScienceDaily LCC. <http://www.sciencedaily.com/releases/2006/08/060807121419.htm> (2006-11-23).
- [7] Lim, D. V.; Simpson, J.M.; Kearns, E.A.; and Kramer, M.F. (2005). Current and Developing Technologies for Monitoring Agents of Bioterrorism and Biowarfare. *Clin. Microbio. Rev.*, 18(4), 583–607.
- [8] Yang, J.K.; Barnett, M.O.; Zhuang, J.; Fendorf, S.E.; and Jardine, P.M. (2005). Adsorption, Oxidation, and Bioaccessibility of As(III) in Soils, *Environ. Sci. Technol.*, 39, 7102–7110.
- [9] Szakova, J.; Tlustos, P.; Goessler, W.; Pavlikova, D.; Balik, J.; and Schlagenhaufen, C. (2005). Comparison of Mild Extraction Procedures for Determination of Plant-Available Arsenic Compounds in Soil. *Anal. Bioanal. Chem.*, 382, 142–148.
- [10] Impellitter, C.A.; Lu, Y.; Saxe, J.K.; Allen, H.E.; and Peijnenburg, W.J.G.M. (2002). Correlation of the Partitioning of Dissolved Organic Matter Fractions with the Desorption of Cd, Cu, Ni, Pb and Zn from 18 Dutch Soils. *Environ. Int.*, 28, 401–410.
- [11] Dabrowska, H.; Dabrowski, L.; Biziuk, M.; Gaca, J.; and Namiesnik, J. (2003). Solid-Phase Extraction Clean-Up of Soil and Sediment Extracts for the Determination of Various Types of Pollutants in a Single Run. *J. Chromatog. A*, 1003, 29–42.
- [12] Baranowska, I. and Piesko, C. (2000). Determination of Selected Herbicides and Phenols in Water and Soils by Solid-Phase Extraction and High-Performance Liquid Chromatography. *J. Chromatogr. Sci.*, 38, 211–218.

- [13] Oleszczuk, P. and Baran, S. (2004). Application of Solid-Phase Extraction to Determination of Polycyclic Aromatic Hydrocarbons in Sewage and Sludge Extracts. *J. Haz. Mat. B*, 113, 237–245.
- [14] Dabrowska, H.; Dabrowski, L.; Biziuk, M.; Gaca, J.; and Namiesnik, J. (2003). Solid-Phase Extraction Clean-Up of Soil and Sediment Extracts for the Determination of Various Types of Pollutants in a Single Run. *J. Chromatog. A*, 1003, 29–42.
- [15] Bolygo, E.; Ricketts, D.; Moffatt, F.; Jessop, M.; and Attenborough, S. (2003). Determination of Antifungal Proteins in Soil by Liquid Chromatography. *Anal. Bioanal. Chem.*, 376, 701–705.
- [16] Schoenborn, L.; Yates, P.S.; Grinton, B.E.; Hugenholtz, P.; and Janssen, P.H. (2004). Liquid Serial Dilution is Inferior to Solid Media for Isolation of Cultures Representative of the Phylum-level Diversity of Soil Bacteria. *Appl. Environ. Microbiol.*, 70(7), 4363–4366.
- [17] Shan, G.; Embrey, S.K.; Herman, R.A.; Wolt, J.D.; Weston, D.; and Mayer, L.M. (2005). Biomimetic Extraction of *Bacillus thuringiensis* Insecticidal Crystal Proteins from Soil Based on Invertebrate Gut Fluid Chemistry. *J. Agric. Chem.*, 53, 6630–6634.
- [18] Deaver, D.R. (1995). A New Non-isotopic Detection System for Immunoassays. *Nature*, 377, 758–760.
- [19] Marquette, C.A. and Blum, L.J. (2006). State of the Art and Recent Advances in Immunoanalytical Systems. *Biosens. Bioelectron.*, 21, 1424–1433.
- [20] Gatto-Menking, D.L.; Yu, H.; Bruno, J.G.; Goode, M.T.; Miller, M.; and Zulich, A.W. (1995). Sensitive Detection of Biotoxoids and Bacterial Spores Using an Immunomagnetic Electrochemiluminescence Sensor. *Biosens. Bioelectron.*, 10, 501–507.
- [21] Yu, H.; Raymonda, J.W.; McMahon, T.M.; and Campagnari, A.A. (2000). Detection of Biological Threats by Immunomagnetic Microsphere-based Solid Phase Fluorogenic- and electro-chemiluminescence. *Biosens. Bioelectron.*, 14, 829–840.
- [22] Grimshaw, C.; Gleason, C.; Chojnicki, E.; and Young, J. (1997). Development of an Equilibrium Immunoassay Using Electrochemiluminescent Detection for a Novel Recombinant Protein Product and its Application to Pre-clinical Product Development. *J. Pharm. Biomed. Anal.*, 16, 605–612.
- [23] Kijek, T.M., Rossi, C.A., Moss, D., Parker, R.W., Henchal, E.A. (2000). Rapid and Sensitive Immunomagnetic–electrochemiluminescent Detection of Staphylococcal enterotoxin B. *J. Immunol. Meth.*, 236, 9-17.
- [24] Shelton, D.R. and Karns, J.S. (2001). Quantitative Detection of *Escherichia coli* O517 in Surface Waters by using Immunomagnetic Electrochemiluminescence. *Appl. Environ. Microbiol.*, 67(7), 2908–2915.

- [25] Lyon, W.J. and Niemeyer, D. (n.d.) Electrochemiluminescence Detection of Biological Agents in Foods and Environmental Sample Matrices. (Online) BioVeris Corporation. [http://www.bioveris.com/homelandsecurity/instruments/Lyons\\_WP.pdf](http://www.bioveris.com/homelandsecurity/instruments/Lyons_WP.pdf) (2006-11-24).
- [26] Fulton, R.E. and Stadnyk, L.L. (2001). Enzyme-Linked Immunosorbant Assays for Identification of Biological Agents in Sample Unknowns: NATO SIBCA Exercise III. (DRES Technical Report 2001-175). Defence Research Establishment Suffield.

## Annex A Running a Plate on the M1M

---

1. Turn on the M1M analyzer
2. Turn on the computer and log into Windows. The M-Series M1M analyzer software loads and displays automatically. If it doesn't, double click on the M1M icon on the desktop.
3. Click **New Plate** from the *M-Series M1M Analyzer* window.
4. Enter a title in the *Plate Name* box of the *New Plate Info* box.
5. Click **Next Step** in the bottom right of the window.
6. Click on **Test Kit ID** field. Scan the first Bar Code (A) of the CRP kit package and the box will fill with a code.
7. Click on **Sample ID** field. Enter a sample name such as "unknown no. 1".
8. If doing more than one replicate of each sample, click on **Number of Replicates** and enter number of replicates being tested.
9. Click on **Next Sample**.
10. Enter the second sample name in the **Sample ID** field.
11. When all samples have been entered, click **Next Step** in the bottom right of the window. *Control Lot #* window is displayed.
12. Click **OK**. *Final Template* window is displayed.
13. Check to see if as desired. If not, back up a step. If OK, click **Done**. The *Run a Plate* window is displayed
14. Verify that both the "Plate Summary" and the "Samples & Controls" information are correct. To make corrections, highlight a test in the *Samples & Controls* area and use **Add**, **Edit**, or **Delete**.

### Sample Preparation

- a. Reconstitute PCM and PCL with 250  $\mu$ L of BV-DILUENT to each vial
- b. Open bar-coded package containing tubes of lyophilized material.
- c. Place tubes in Tube Rack according to plate template.
- d. Add 50  $\mu$ L BV-DILUENT to negative control (NC) wells.

- e. Add 50  $\mu$ L of CRP Positive Control Medium (PCM) to specified wells.
  - f. Add 50  $\mu$ L of CRP Positive Control Low (PCL) to specified wells.
  - g. Add 50  $\mu$ L of sample (for example: unknown no. 1) to test wells.
  - h. Add 50  $\mu$ L BV-DILUENT to all wells being used.
15. Place the prepared plate in the plate tray with well A1 towards the back of the MIM analyzer.
  16. Click **Run** on the toolbar to start the plate reading process. As results are calculated, they are displayed in the *Plate Layout* area.
  17. Check the plate status at the bottom of the *Plate Layout* area to track progress of the plate read.
  18. Click **Eject** from the toolbar to eject the plate tray when the plate read is completed

## Annex B Powder Swab Protocol

---

1. Remove swab from wrapping
2. Dip swab in vial or tube of BV-DILUENT and thoroughly wet swab
3. Remove swab from liquid and press swab against sides of tube to remove excess liquid
4. Roll swab in powder to thoroughly coat swab
5. Hold swab between thumb and middle finger and use index finger to firmly tap stem of swab to remove excess powder (as though flicking the ash from a cigarette).
6. Dip coated swab in tube containing 1 mL BV-DILUENT, manually agitating swab in the liquid for about 10 seconds.
7. Remove swab from liquid and press swab against sides of tube to remove excess liquid.
8. Thoroughly mix powder in liquid by shaking or vortexing tube, and pour or pipette contents of tube into the BioVeris Sample Prep Filter.
9. Remove bottom cap of Sample Prep Filter and repeatedly squeeze below filter disk to empty filtered fluid into a clean tube. At least 50  $\mu$ L is needed to run an assay. Filtered liquid may be cloudy.
10. Proceed with assay as described in the manual
11. Soil Swab Protocol assay as described in the manual

## Annex C Soil Swab Protocol

---

1. Remove swab from wrapping
2. Dip swab in vial or tube of BV-DILUENT and thoroughly wet swab
3. Remove swab from liquid and press swab against sides of tube to remove excess liquid
4. Roll swab in soil to thoroughly coat swab
5. Hold swab between thumb and middle finger and use index finger to firmly tap stem of swab to remove excess soil (as though flicking the ash from a cigarette).
6. Dip coated swab in tube containing 2 mL BV-DILUENT, manually agitating swab in the liquid for about 10 seconds.
7. Remove swab from liquid and press swab against sides of tube to remove excess liquid.
8. Thoroughly mix soil in liquid by shaking or vortexing tube, and pour or pipette contents of tube into the BioVeris Sample Prep Filter.
9. Remove bottom cap of Sample Prep Filter and repeatedly squeeze below filter disk to empty filtered fluid into a clean tube. At least 50  $\mu$ L is needed to run an assay. Filtered liquid may be cloudy.
10. Proceed with assay as described in the manual

**DOCUMENT CONTROL DATA**

(Security classification of title, body of abstract and indexing annotation must be entered when the overall document is classified)

<p>1. ORIGINATOR (The name and address of the organization preparing the document. Organizations for whom the document was prepared, e.g. Centre sponsoring a contractor's report, or tasking agency, are entered in section 8.)</p> <p>Defence R&amp;D Canada – Suffield P.O. Box 4000, Station Main Medicine Hat, Alberta T1A 8K6</p>	<p>2. SECURITY CLASSIFICATION (Overall security classification of the document including special warning terms if applicable.)</p> <p align="center">UNCLASSIFIED</p>	
<p>3. TITLE (The complete document title as indicated on the title page. Its classification should be indicated by the appropriate abbreviation (S, C or U) in parentheses after the title.)</p> <p align="center">Development of Sampling Methods for Powders and Soil for Detection of Biothreat Agents by Electrochemiluminescence</p>		
<p>4. AUTHORS (last name, followed by initials – ranks, titles, etc. not to be used)</p> <p align="center">Thompson, H.G.; Fulton, R.E.; and Ranches, J.</p>		
<p>5. DATE OF PUBLICATION (Month and year of publication of document.)</p> <p align="center">June 2007</p>	<p>6a. NO. OF PAGES (Total containing information, including Annexes, Appendices, etc.)</p> <p align="center">30</p>	<p>6b. NO. OF REFS (Total cited in document.)</p> <p align="center">26</p>
<p>7. DESCRIPTIVE NOTES (The category of the document, e.g. technical report, technical note or memorandum. If appropriate, enter the type of report, e.g. interim, progress, summary, annual or final. Give the inclusive dates when a specific reporting period is covered.)</p> <p align="center">Technical Memorandum</p>		
<p>8. SPONSORING ACTIVITY (The name of the department project office or laboratory sponsoring the research and development – include address.)</p> <p align="center">Defence R&amp;D Canada – Suffield P.O. Box 4000, Station Main Medicine Hat, Alberta T1A 8K6</p>		
<p>9a. PROJECT OR GRANT NO. (If appropriate, the applicable research and development project or grant number under which the document was written. Please specify whether project or grant.)</p>	<p>9b. CONTRACT NO. (If appropriate, the applicable number under which the document was written.)</p>	
<p>10a. ORIGINATOR'S DOCUMENT NUMBER (The official document number by which the document is identified by the originating activity. This number must be unique to this document.)</p> <p align="center">DRDC Suffield TM 2007-172</p>	<p>10b. OTHER DOCUMENT NO(s). (Any other numbers which may be assigned this document either by the originator or by the sponsor.)</p>	
<p>11. DOCUMENT AVAILABILITY (Any limitations on further dissemination of the document, other than those imposed by security classification.)</p> <p align="center">Unlimited</p>		
<p>12. DOCUMENT ANNOUNCEMENT (Any limitation to the bibliographic announcement of this document. This will normally correspond to the Document Availability (11). However, where further distribution (beyond the audience specified in (11) is possible, a wider announcement audience may be selected.)</p> <p align="center">Unlimited</p>		

13. **ABSTRACT** (A brief and factual summary of the document. It may also appear elsewhere in the body of the document itself. It is highly desirable that the abstract of classified documents be unclassified. Each paragraph of the abstract shall begin with an indication of the security classification of the information in the paragraph (unless the document itself is unclassified) represented as (S), (C), (R), or (U). It is not necessary to include here abstracts in both official languages unless the text is bilingual.)

The MIM (BioVeris Corporation), an automated electrochemiluminescent (ECL) assay instrument, has been shown to be a technology that is robust and tolerant of samples in a variety of complex matrices; however, methods for sampling of powders and soils in the field for processing on the MIM have not been described. This report describes the development of simple and reproducible methods that may be used for field sampling of powders and soils for subsequent screening and analysis on the MIM. By sampling with polyester swabs using a standard method, it was found that it was possible to collect a reproducible amount (~12 mg) of powder or soil from a variety of surfaces. Using this value, the volume of diluent required to produce a signal <1.2 on the MIM was empirically determined. The powder or soil suspension was filtered and the filtrate assayed on the MIM. The result was a method that could be used consistently to assay powder and soil by ECL MIM with good sensitivity and minimal background signal.

14. **KEYWORDS, DESCRIPTORS or IDENTIFIERS** (Technically meaningful terms or short phrases that characterize a document and could be helpful in cataloguing the document. They should be selected so that no security classification is required. Identifiers, such as equipment model designation, trade name, military project code name, geographic location may also be included. If possible keywords should be selected from a published thesaurus, e.g. Thesaurus of Engineering and Scientific Terms (TEST) and that thesaurus identified. If it is not possible to select indexing terms which are Unclassified, the classification of each should be indicated as with the title.)

electrochemiluminescence, immunoassay, bioassay, ricin, VEE, Francisella tularensis, powder, soil

## **Defence R&D Canada**

Canada's Leader in Defence  
and National Security  
Science and Technology

## **R & D pour la défense Canada**

Chef de file au Canada en matière  
de science et de technologie pour  
la défense et la sécurité nationale



[www.drdc-rddc.gc.ca](http://www.drdc-rddc.gc.ca)