

# UNCLASSIFIED

AD NUMBER
AD835495
NEW LIMITATION CHANGE
TO Approved for public release, distribution unlimited
FROM Distribution authorized to U.S. Gov't. agencies and their contractors; Critical Technology; APR 1968. Other requests shall be referred to Department of the Army, Fort Detrick, Attn: Technical Releases Branch, Technical Information Div., Frederick, MD 21701.
AUTHORITY
Army Biological Defense Research Lab ltr dtd 29 Sep 1971

THIS PAGE IS UNCLASSIFIED

AD

AD 835495

TECHNICAL MANUSCRIPT 445

USE OF PLASTIC T-30  
TISSUE CULTURE FLASKS  
UNDER TIME-LAPSE  
PHASE-CONTRAST MICROSCOPY

Louis E. Schneider  
Robert R. Rosato

APRIL 1968

UNCLASSIFIED

This  
trans-  
mission

controls and each  
national may be

DEPARTMENT OF THE ARMY

Fort Detrick

Frederick, Maryland 21701

Attn: Technical Release Branch  
Technical Information Div

2  
Reproduction of this publication in whole or in part is prohibited except with permission of the Commanding Officer, Fort Detrick, ATTN: Technical Releases Branch, Technical Information Division, Fort Detrick, Frederick, Maryland, 21701. However, DDC is authorized to reproduce the publication for United States Government purposes.

#### DDC AVAILABILITY NOTICES

Qualified requesters may obtain copies of this publication from DDC.

Foreign announcement and dissemination of this publication by DDC is not authorized.

Release or announcement to the public is not authorized.

#### DISPOSITION INSTRUCTIONS

Destroy this publication when it is no longer needed. Do not return it to the originator.

The findings in this publication are not to be construed as an official Department of the Army position, unless so designated by other authorized documents.

DEPARTMENT OF THE ARMY  
Fort Detrick  
Frederick, Maryland 21701

TECHNICAL MANUSCRIPT 445

USE OF PLASTIC T-30 TISSUE CULTURE FLASKS UNDER  
TIME-LAPSE PHASE-CONTRAST MICROSCOPY

Louis E. Schneider

Robert R. Rosato

Process Development Division  
AGENT DEVELOPMENT AND ENGINEERING LABORATORY

Project 1B533001D426

April 1968

USE OF PLASTIC T-30 TISSUE CULTURE FLASKS UNDER  
TIME-LAPSE PHASE-CONTRAST MICROSCOPY

ABSTRACT

With the long-working-distance phase-contrast condensers now available, it is possible to use plastic T-30 tissue culture flasks in phase-contrast time-lapse cinephotomicrographic studies. Data and photographs obtained with such a system are presented, and the advantages of using plastic flasks are discussed.

A continuing need exists to study both the gross effect of viral infection upon mammalian cells and the precise generation times of uninfected mammalian cells. This dual need has been filled by the use of Sage cinephotomicrographic equipment\* that incorporates a Bolex 16-mm motion-picture camera\*\* and a Zeiss-Plankton inverted microscope supplied with a Zeiss Model I-S Zernicke-type phase-contrast condenser. This particular condenser was chosen for its long working distance, which permits the phase-contrast illumination of plastic T-30 tissue culture flasks\*\*\* when a 10X/0.22 Zeiss phase-contrast objective is used. Some data obtained with this system are presented in this study, which further considers both the advantages and difficulties encountered in use of plastic T-30 tissue culture flasks.

Living material studied under phase-contrast illumination reveals structural details otherwise visible only in killed and stained preparations. The ability to use phase-contrast illumination in conjunction with plastic T 30 tissue culture flasks is desirable for (i) its convenience, (ii) the availability of reliably clean, optical-quality growth surfaces, and most important, (iii) the retention of growth conditions identical with those of the normal laboratory propagation of monolayer cultures. Thus, data obtained can be directly related to parallel cultures carried in plastic T-30 flasks. Less valid comparative

---

\* Sage Instruments, Inc., White Plains, New York 10601.

\*\* Model H-16M, Paillard S.A., Ste-Croix, Switzerland.

\*\*\* Falcon Plastics, Los Angeles, California 90045.

relationships can be drawn if the usual perfusion chamber, with its glass cover slips and sealing gasket, is used because of (i) the difference in the growth surfaces, both in cleanliness and in composition, and (ii) the absence of the prominent gas phase, which is characteristic of the normally propagated monolayer cultures.

Perfusion chambers are still needed for studies requiring the greatest resolution. Most microscope objectives are designed for use with glass cover slips that are either 0.17 or 0.18 mm thick (the former value is used by Continental manufacturers, the latter by English and American manufacturers).<sup>1</sup> A recent study<sup>2</sup> describes the incorporation of such cover slips into a plastic\* Cooper dish,<sup>3</sup> which is then sealed at the edges. This system permits observation under high-resolution objectives and also retains a considerable gas phase.

Because high resolution was not required in these studies, the considerable thickness (0.9 mm) of the floor in the plastic tissue culture flask was not objectionable. Figure 1 shows an L-cell monolayer maintained in a plastic T-30 culture flask containing 5 ml of tissue culture medium with 10% bovine serum.\*\* The cells were held under phase-contrast illumination and photographed on a 16-mm motion picture film\*\*\* by means of electronic flash.\*\*\*\* The cells were photographed under a 10X phase-contrast objective and a 5X planar eyepiece and at the 2X magnification of a Zeiss Optovar attachment.

The Optovar is a turret-type magnifier located between the microscope objective and eyepiece. In addition to the four magnification powers, 1X, 1.25X, 1.6X, and 2X, the Optovar also contains a magnifier for aligning the phase rings. Using this attachment, one can easily change the magnification of the field or check the alignment of the phase-contrast illumination during the course of a time-lapse experiment.

---

\* Falcon Plastics, Los Angeles, California 90045.

\*\* Rosensteel, J.F.; Jordan, W.C.; Daniels, W.F. Growth potential of a variation of Eagle's minimal essential media for spinner culture. Unpublished data.

\*\*\* Kodak Ektachrome MS, daylight.

\*\*\*\* Flashtube FT-230, General Electric Co.



Figure 1. L-Cells in a Plastic T-30 Tissue Culture Flask, Phase-Contrast Illumination, Time-Lapse Cinephotomicrography. 200X.

Time-lapse studies for the determination of precise generation times were performed on the greatest possible field of view, achieved by using the Optovar at 1X, with a 10X phase-contrast objective and a 5X planar eyepiece. At this magnification and at a time-lapse speed of one frame every 2 minutes, the contrast proved adequate for the differentiation of cellular cleavage, and even cellular death, between two successive motion picture frames. Figure 2 is a diagram of the growth and progeny of a single mouse L cell maintained and photographed under the conditions just described. The progeny of this cell at 80.2 hours of growth are shown in Figure 3. Two of the cells are dead (the round cells), and two are still alive (the cells with pseudopods). The pseudopods are indicative of the cellular mobility that necessitated a single recentration of the cells to keep them all in view.

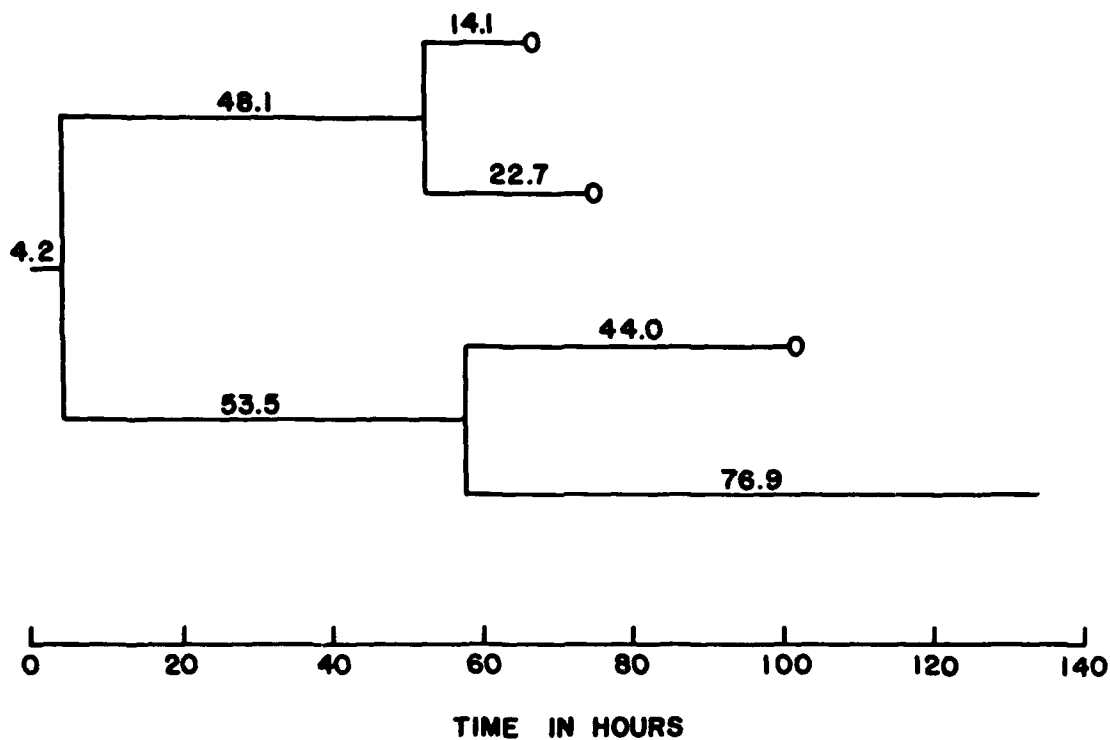


Figure 2. Growth and Progeny of a Single L-Cell from a Culture in a Plastic T-30 Tissue Culture Flask, Phase-Contrast Illumination, Time-Lapse Cinephotomicrography.  $\circ$  = death. Vertical line = cellular cleavage. Horizontal line = elapsed time in hours between cellular cleavages.



Figure 3. Progeny at 80.2 Hours of the Single L-Cell Described in Figure 2. 100X.

An unforeseen difficulty was encountered in the use of plastic T-30 flasks. Although the cellular monolayer was in focus at the start of an experiment, the cells were considerably out of focus less than a day later. After the microscope itself was eliminated as the source of the difficulty, we discovered that the floor and the ceiling of the flasks were not rigid enough to prevent their minute bulging because of gaseous expansion within the sealed flask. Because preliminary temperature equilibration of the flasks and the cultures did not solve the difficulty, we suspected that the cells themselves might be modifying the sealed atmosphere. For this reason, the following simple mechanical expedient was devised. A 1/8-inch-diameter hole was drilled into the center of the tops of spare Bakelite caps without permitting the drill to pierce the white rubber lining. The holes were filled flush with a silicone-rubber compound,\* and the caps were conveniently autoclaved in a glass petri dish. Such a modified cap was aseptically substituted for the regular cap of a monolayer-containing plastic T-30 flask and tightly screwed down. A sterile, disposable hypodermic needle (19-gauge, 1.5-inch) containing sterile, nonabsorbent cotton in the hub was then aseptically pushed through the rubber port in the cap.

The preparation of separate, sterile caps and needles was the most suitable method when infectious material was to be studied. The cap with a port was added in a biological safety cabinet, and the cap was sealed to the neck of the flask with plastic tape. The outside of the flask was decontaminated as it was passed from the safety cabinet to a bench top, where the needle was inserted. When noninfectious material was to be studied, the drilled cap and needle were assembled as a single unit prior to sterilization.

During the preparation of the monolayers, the flasks were carefully examined to see that medium was not escaping through an improperly sealed seam. Leaky flasks were reported by Petrusson and Fogh,<sup>4</sup> and we also discovered some. As an additional precaution, particularly when infectious material is to be studied, a bead of silicone-rubber compound can be run all along the seam of the flask.

The inoculating needle inserted through the cap served as an adequate vent for releasing the pressure from the flask without allowing a loss of carbon dioxide sufficient to change the pH of the medium. Flasks vented in this fashion have been used in 5-day time-lapse experiments involving as many as 100,000 cells per ml (500,000 cells per flask). At the end of each experiment the cells were viable (as evidenced by recent division) and still in focus, and the pH of the medium was identical to that of an unvented control flask.

---

\* Lab RTV, Dow Corning Corp., Midland, Michigan 48640.

LITERATURE CITED

1. Norris, K.P. 1961. Some observations on microscope cover-glasses. J. Roy. Microscop. Soc. 79:287-298.
2. Deitch, A.D.; Godman, G.C. 1966. Modification of plastic culture dishes for high resolution microscopy of living cells. Stain Technol. 41:19-22.
3. Cooper, W.G. 1961. A modified plastic petri dish for cell and tissue cultures. Proc. Soc. Exp. Biol. Med. 106:801-803.
4. Petursson, G.; Fogh, J. 1963. High powered phase contrast equipment with long working distance for observation of culture flasks. Proc. Soc. Exp. Biol. Med. 112:523-524.

PRECEDING  
PAGE BLANK

Unclassified

Security Classification

13

DOCUMENT CONTROL DATA - R & D		
(Security classification of title, body of abstract and indexing annotation must be entered when the overall report is classified)		
1. ORIGINATING ACTIVITY (Corporate author)		2a. REPORT SECURITY CLASSIFICATION
Department of the Army Fort Detrick, Frederick, Maryland, 21701		Unclassified
		2b. GROUP
3. REPORT TITLE		
USE OF PLASTIC T-30 TISSUE CULTURE FLASKS UNDER TIME-LAPSE PHASE-CONTRAST MICROSCOPY		
4. DESCRIPTIVE NOTES (Type of report and inclusive dates)		
5. AUTHOR(S) (First name, middle initial, last name)		
Louis E. Schneider Robert R. Rosato		
6. REPORT DATE	7a. TOTAL NO. OF PAGES	7b. NO. OF REFS
April 1968	13	4
8a. CONTRACT OR GRANT NO.	8b. ORIGINATOR'S REPORT NUMBER(S)	
a. PROJECT NO. 1B533001D426	Technical Manuscript 445	
c.	9b. OTHER REPORT NO(S) (Any other numbers that may be assigned this report)	
d.		
10. DISTRIBUTION STATEMENT		
Qualified requesters may obtain copies of this publication from DDC. Foreign announcement and dissemination of this publication by DDC is not authorized. Release or announcement to the public is not authorized.		
11. SUPPLEMENTARY NOTES		12. SPONSORING MILITARY ACTIVITY
		Department of the Army Fort Detrick, Frederick, Maryland, 21701
13. ABSTRACT		
With the long-working-distance phase-contrast condensers now available, it is possible to use plastic T-30 tissue culture flasks in phase-contrast time-lapse cinephotomicrographic studies. Data and photographs obtained with such a system are presented, and the advantages of using plastic flasks are discussed.		
14. Key Words		
*Plastics Tissue cultures *Flask *Cinephotomicrography *Phase microscope		

DD FORM 1473 REPLACES DD FORM 1473, 1 JAN 66, WHICH IS OBSOLETE FOR ARMY USE.

Unclassified

Security Classification