

AD-A118 334 ARMY MOBILITY EQUIPMENT RESEARCH AND DEVELOPMENT COMM--ETC F/G 7/4
IDENTIFICATION OF DRYING OILS IN PAINTS BY GAS-LIQUID CHROMATOGRAPH--ETC(U)
JUN 82 L. HERRIOTT
UNCLASSIFIED MERADCOM-2360 NL

1 of 1
AD-A
118 334

END
9-82

AD A118334

DTIC FILE COPY



AD

12


Report 2360

IDENTIFICATION OF DRYING OILS IN PAINTS
BY GAS-LIQUID CHROMATOGRAPHY

by
Laurie Herriott

June 1982

Approved for public release; distribution unlimited

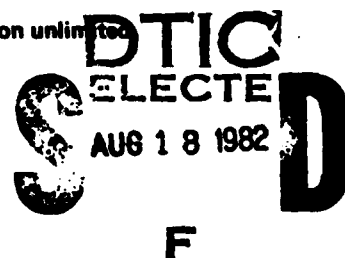

U.S. ARMY MOBILITY EQUIPMENT
RESEARCH AND DEVELOPMENT COMMAND
FORT BELVOIR, VIRGINIA

82 08 17 063

**Destroy this report when it is no longer needed.
Do not return it to the originator.**

**The citation in this report of trade names of
commercially available products does not constitute
official endorsement or approval of the use of such
products.**

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER 2360	2. GOVT ACCESSION NO. AD-A118334	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) IDENTIFICATION OF DRYING OILS IN PAINTS BY GAS-LIQUID CHROMATOGRAPHY		5. TYPE OF REPORT & PERIOD COVERED Interim
		6. PERFORMING ORG. REPORT NUMBER
7. AUTHOR(s) Laurie Herriott		8. CONTRACT OR GRANT NUMBER(s)
9. PERFORMING ORGANIZATION NAME AND ADDRESS Material Technology Lab; ATTN: DRDME-VO; US Army Mobility Equipment Research and Development Command; Fort Belvoir, VA 22060		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS
11. CONTROLLING OFFICE NAME AND ADDRESS		12. REPORT DATE June 1982
		13. NUMBER OF PAGES 22
14. MONITORING AGENCY NAME & ADDRESS (If different from Controlling Office)		15. SECURITY CLASS. (of this report) Unclassified
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE
16. DISTRIBUTION STATEMENT (of this Report) Approved for public release; distribution unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)		
18. SUPPLEMENTARY NOTES		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number)		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number) This gas-liquid chromatographic method is for the identification of drying oils in paint. The identification is made from the distribution of fatty acids as methyl esters shown on the chromatogram.		

DD FORM 1 JAN 73 1473


EDITION OF 1 NOV 65 IS OBSOLETE

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

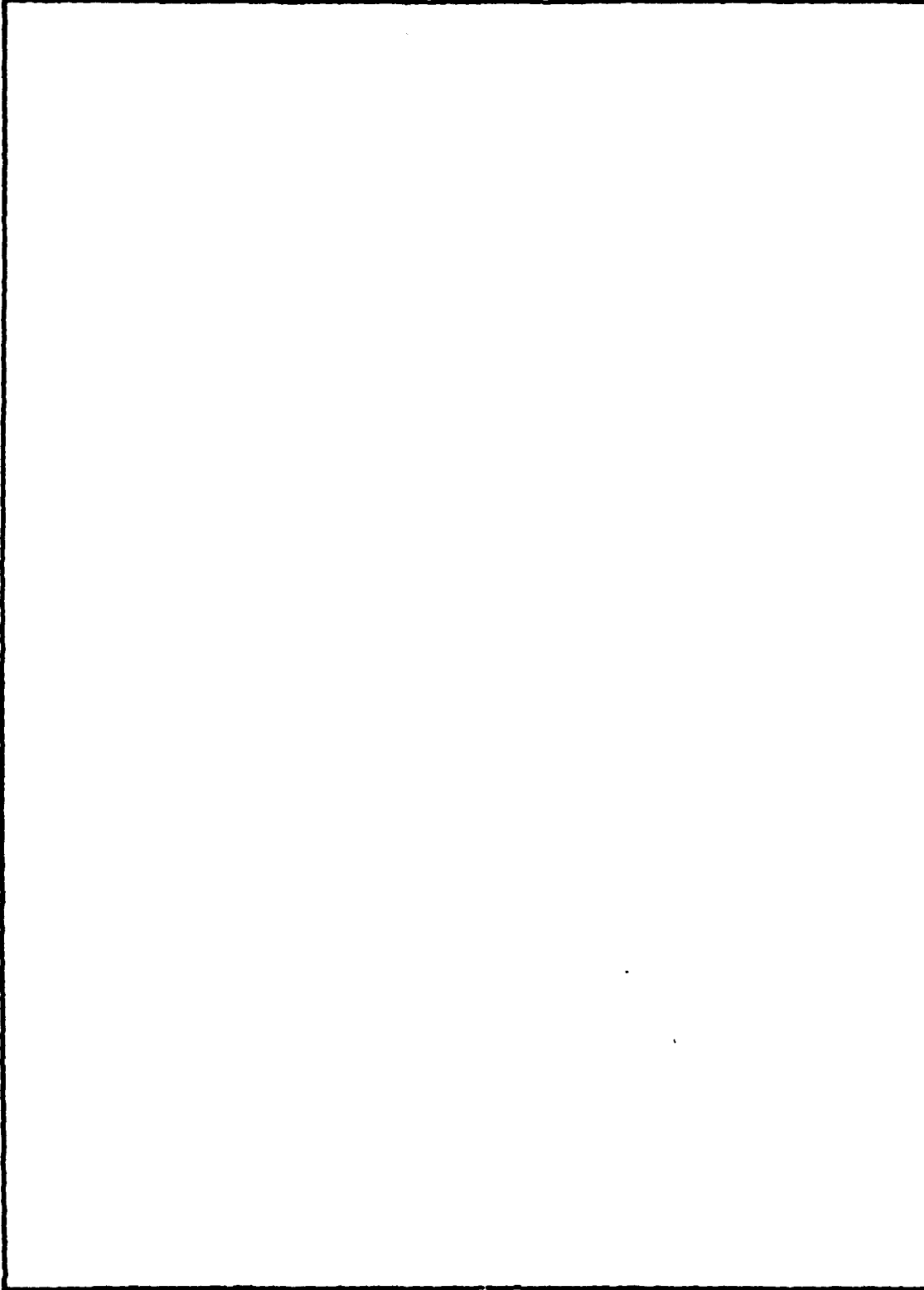

UNCLASSIFIED

ii

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

PREFACE

The work covered by this report was performed by Laurie Herriott and was reviewed by Troy Nichols under the supervision of Emil J. York, Chief, Material Technology Laboratory. The effort constitutes a part of the U.S. Army Materials and Mechanics Research Center's program to provide better procedures in chemical testing techniques to increase inspection efficiency for products procured by the U.S. Army Materiel Development and Readiness Command.



Accession For	
NTIS GRA&I	<input checked="checked" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution/	
Availability Codes	
Dist	Avail and/or Special
A	

CONTENTS

Section	Title	Page
	PREFACE	iii
	ILLUSTRATIONS	v
	TABLES	v
I	INTRODUCTION	
	1. Purpose	1
	2. Background	1
II	DETAILS OF TESTS	
	3. Reagents	1
	4. Gas-Liquid Chromatographic Column	1
	5. Sample Preparation	2
	6. Identification of Ester and Oil	3
III	RESULTS AND DISCUSSION	5
IV	CONCLUSION	5

ILLUSTRATIONS

Figure	Title	Page
1	Chromatogram of Soybean Oil Methyl Ester	6
2	Chromatogram of Castor Oil Methyl Ester	7

TABLES

Table	Title	Page
1	Drying Oils Used in Sample Preparation	2
2	Relative Retention Data of Methyl Esters of Oil Fatty Acids	3
3	Weight Percent Distribution of Fatty Acids in Unmodified Oils	4

IDENTIFICATION OF DRYING OILS IN PAINTS

BY GAS-LIQUID CHROMATOGRAPHY

I. INTRODUCTION

1. Purpose. The purpose of this investigation was to revise and update Federal Test Method Standard (FTMS) 7501, "Identification of Drying Oils (Gas-Liquid Chromatography)."

2. Background. Some of the military specifications for paints specify that certain drying oils are to be used in the formulation. If those drying oils are not used, the paint will prove unsatisfactory in its intended application. Therefore, it is essential for purposes of quality control that an accurate method be available for the identification of these drying oils. In September 1965, Method 7501 of FTMS 141 was approved for this purpose, but it has since been shown to be lacking in necessary flexibility for application to actual samples of specification paints.

II. DETAILS OF TESTS

3. Reagents.

a. Boron Trifluoride solution¹—add 23 ml of boron fluoride methyl ether slowly to 150 ml of chilled absolute methyl alcohol in a 200-ml volumetric flask and dilute the volume with absolute methanol. Stopper tightly and store in a refrigerator. Discard after one month.

b. Petroleum Ether.²

4. Gas-Liquid Chromatographic Column. Prepare a 6-ft by 1/8-in. glass column, packed with 20-percent diethylene glycol succinate on 60-80 mesh chromosorb W acid washed (WAW). Using the thermal conductivity detector (TCD), condition the column at 200° C with helium at a flow rate of 60 ml/min, until the baseline is stabilized.

¹ L. D. Metcalfe and A. A. Schmitz, J. Am. Oil Chem. Soc., 33, 361 (1961).

² L. D. Metcalfe, A. A. Schmitz, and J. R. Felka, Anal. Chem., 38, 514 (1966).

5. Sample Preparation. Isolate the fatty acids from the resin or oil by saponification, acid hydrolysis, and solvent separation as in FTMS 7021 and FTMS 7031. If the fatty acids have become solid, warm them slightly. Then, accurately weigh about 0.2 g in a 6-in. screw cap test tube. Add 6 ml of the Boron Trifluoride solution and immerse in a 65° to 70° C water bath for 10 min. Transfer to a separatory funnel with 20 ml water and 30 ml of petroleum ether. Shake gently. Discard the bottom water layer and filter the remaining layer through fast filter paper into a 50-ml beaker and evaporate in a 65° to 70° C water bath to about 2 ml. Do not evaporate to dryness. Transfer and store in a closed 2-ml vial. Monitoring the specification or statement of composition, inject 1.5 ml of the prepared methyl esters using the optimum conditions for each oil as shown in Table 1.

Table 1. Drying Oils Used in Sample Preparation

Drying Oils	Column Temperature (° C)	Carrier (Helium) Flow (ml/min)
Castor	20 min at 170°; 50 min at 200°; rate 20° C/min	60
Soya	180	60
Cottonseed	180	60
Linseed	180	50
Safflower	190	40

Thermal conductivity detector 250° C.
Injection port temperature 275° C.

6. Identification of Ester and Oil. The individual methyl esters can be identified by their relative retention times as given in Table 2. The drying oil, itself, may be identified in either of two ways. First, by comparison of its chromatogram with those of known oils. Second, identification of the drying oil may be made by matching concentration of the methyl esters with the values given in Table 3. For this purpose, concentration (C_x) of the methyl ester is calculated as follows:

$$C_x = \frac{A_x}{A_T} \times F_x \times 100.$$

Where A_x is the peak area of the methyl ester; A_T is the total area of the methyl ester peaks; and F_x is a correction factor determined from a known sample.

Table 2. Relative Retention Data of Methyl Esters of Oil Fatty Acids

Ester	Relative Retention Time
Caprylic	0.10
Capric	0.19
Lauric	0.32
Myristic	0.56
Palmitic	1.0
Palmitoleic	1.15
Stearic	1.7
Oleic	1.94
Linoleic	2.31
Linolenic	2.97
Licanic	5.2
Eleostearic	5.3
Ricinoleic	10.0

(Ester of Palmitic Acid = 1)

Table 3. Weight Percent Distribution of Fatty Acids in Unmodified Oil

Oil	Caprylic	Capric	Lauric	Myristic	Palmitic	Palmitoleic	Stearic	Oleic	Linoleic	Linolenic	Licanic	Eleostearic	Ricinoleic
Castor					2		1	7	3				87
Coconut	6	6	44	18	11		6	7	2	t			
Cottonseed		t	t	1	29	2	4	24	40				
Linseed					6	t	4	22	16	52			
Oiticica					7		5	6					78
Perilla					7		2	13	14	64			
Safflower				t	8	t	3	13	75	1			
Soybean				t	11		4	25	51	9			
Tall					5		3	46	41	3			
Tung					4		1	8	4	3			
Menhaden		t	7	16	16	2	15	7	2				

t = trace

III. RESULTS AND DISCUSSION

The various drying oils used in paints are composed of specific amounts of fatty acids as given in Table 3. Determination of these fatty acids as methyl esters by gas-liquid chromatography provides a means for confirmation of the required drying oil in the paint.

Because of the wide variation of composition of drying oils it was impractical, if not impossible, to develop a single chromatographic program for satisfactory identification for all of the oil acids. Carrier flow rate and column oven temperature were adjusted for optimum resolution of the methyl esters for each drying oil.

The one complete method for identification of the oil acid consists of the following steps. The acid is saponified and extracted with water as the sodium salt from an ethyl ether despersion of the paint vehicle. On acidification with hydrochloric acid, the fatty acid is then extracted with chloroform (FTMS 7021). Evaporation^{1 2} follows with conversion of the fatty acid to the methyl ester with Boron Trifluoride solution. This well-known reaction converts the acid to a more volatile product which can then be determined by gas-liquid chromatography. An example of this reaction is:

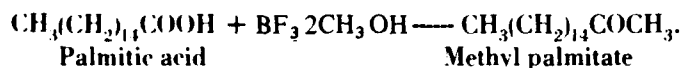


Figure 1 is a chromatogram of the methyl esters of soybean oil obtained with the above procedure. This chromatogram illustrates the use of this method in the identification of drying oils. By its retention time the methyl ester is identified, and by the percentage of each methyl ester, the drying oil is identified.


Figure 2 is the chromatogram of the methyl esters of castor oil. The striking difference between the chromatograms of Figures 1 and 2 illustrate the utilization of this method in identifying drying oils.

V. CONCLUSION

An improved method was developed for the confirmation of the presence of the required drying oil in military specification paints. To a limited extent, the method is useful in the identification of drying oils in nonspecification paints. In some instances, mixtures of drying oils can be determined.


¹ L. D. Metcalfe and A. A. Schmitz, J. Am. Oil Chem. Soc., 33, 361 (1961).

² L. D. Metcalfe, A. A. Schmitz, and J. R. Pelka, Anal. Chem., 38, 514 (1966).

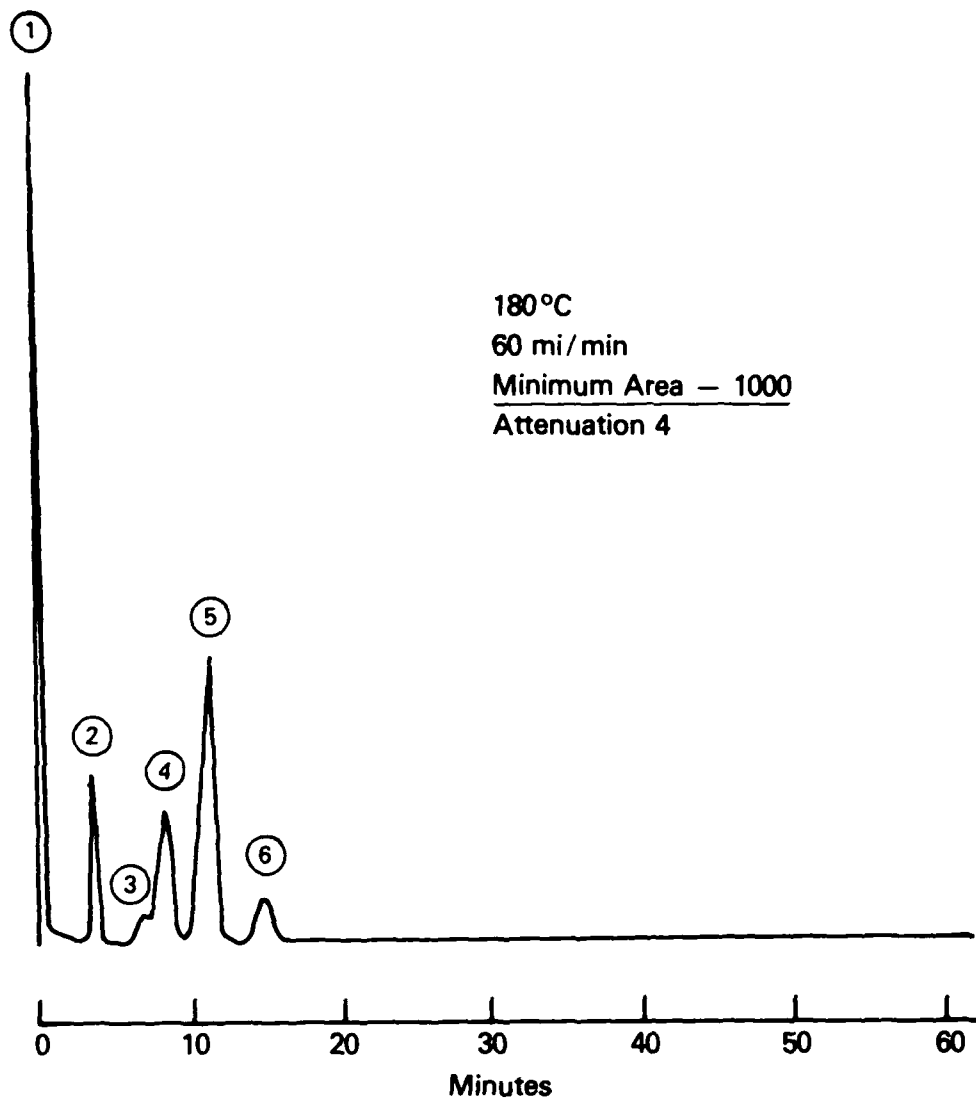

1. Petroleum Ether
2. Palmitic
3. Stearic
4. Oleic
5. Linoleic
6. Linolenic

Figure 1. Chromatogram of soybean oil methyl ester.

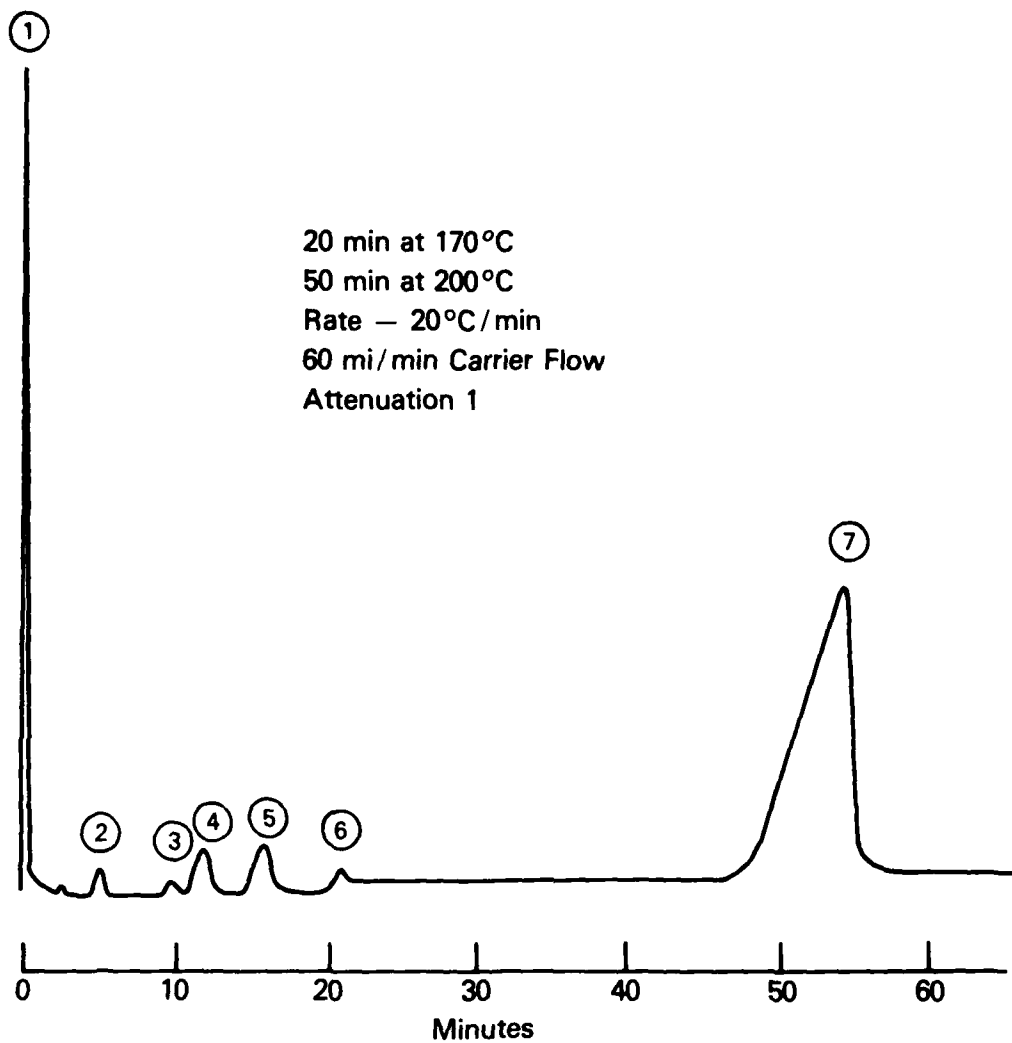

1. Petroleum Ether
2. Palmitic
3. Stearic
4. Oleic
5. Linoleic
6. Impurity
7. Ricinoleic

Figure 2. Chromatogram of castor oil methyl ester.

DISTRIBUTION FOR MERADCOM REPORT 2360

No. Copies	Addressee	No. Copies	Addressee
2	Metals and Ceramics Information Center ATTN: Mr. Harold Mindlin Director Mr. James Lynch Assistant Director 505 King Avenue Columbus, OH 43201	3	Commander US Army Electronics R&D Command ATTN: DRSEL-PA-E Mr. Stan Alster Mr. J. Quinn DELS-D Fort Monmouth, NJ 07703
12	Defense Technical Information Center Cameron Station Alexandria, VA 22314	9	Commander US Army Missile Research and Development Command ATTN: DRDMI-TB, Redstone Scientific Info Ctr (2) DRDMI-TK, Mr. J. Alley DRDMI-M DRDMI-ET, Mr. Robert O. Black DRDMI-QS, Mr. George L. Stewart, Jr. DRDMI-EAT, Mr. R. Talley DRDMI-QP DRSMI-RR Redstone Arsenal, AL 35809
1	Commander US Army Foreign Science and Technology Center ATTN: DRXST-SD3 220 Seventh Street, NE Charlottesville, VA 22901	4	Commander US Army Troop Support and Aviation Materiel Readiness Command ATTN: DRSTS-PL (2) Mr. J. Corwin (1) DRSTS-Q DRSTS-MES (1) DRSTS-M 4300 Goodfellow Blvd St. Louis, MO 63120
2	Office of the Deputy Chief of Staff for Research, Development, and Acquisition ATTN: DAMA-ARZ-E DAMA-CSS Washington, DC 20310	1	Commander US Army Natick Research and Development Command ATTN: DRDNA-EM Natick, MA 01760
2	Commander Army Research Office ATTN: Dr. George Mayer Mr. J. J. Murray P. O. Box 12211 Research Triangle Park, NC 27709	2	Commander US Army Tank-Automotive Materiel Readiness Command ATTN: DRSTA-Q Warren, MI 48090
6	Commander US Army Materiel Development and Readiness Command ATTN: DRCQA-E DRCDE-D DRCDMD-FT DRCLDC DRDCMT DRCMM-M Alexandria, VA 22333		

No. Copies	Addressee	No. Copies	Addressee
6	Commander US Army Armament Materiel Readiness Command ATTN: DRSAR-QA (2) DRSAR-SC DRSAR-RDP DRSAR-EN DRSAE-QAE Rock Island, IL 61299	1	Commander US Army Aberdeen Proving Ground ATTN: STEAP-MT-U (GE Branch) Aberdeen Proving Ground, MD 21005
4	Commander Rock Island Arsenal ATTN: SARRI-EN, Mr. W. M. Kisner SARRI-ENM, Mr. W. D. McHenry SARRI-QA SARRI-LPL Rock Island, IL 61299	1	Director US Army Materiel Systems Analysis Agency ATTN: DRXSY-CM Aberdeen Proving Ground, MD 21005
12	Commander US Army Armament Research and Development Command ATTN: DRDAR-LC, Mr. E. Kelly DRDAR-LCA, Dr. Sharkoff DRDAR-LCE, Dr. Walker DRDAR-QAS, Mr. F. Fitzsimmons (5) DRDAR-SCM, Jr. J. Corrie DRDAR-TSP, Mr. B. Stephans DRDAR-TSS (STINFO)(2) Dover, NJ 07703	1	Director US Army Ballistic Research Lab ATTN: DRDAR-TSD-S (STINFO) Aberdeen Proving Ground, MD 21005
2	Commander Edgewood Arsenal ATTN: DRDAR-CLR, Mr. Montanary DRDAR-QAC, Dr. Maurits Aberdeen Proving Ground, MD 21010	2	Commander Watervliet Arsenal ATTN: DRDAR-LCB Mr. T. Moraczewski SARWV-PPI, Mr. L. Jette Watervliet, NY 12189
1	Technical Library Chemical Systems Laboratory Aberdeen Proving Ground, MD 21010	4	Commander US Army Aviation R&D Command ATTN: DRDAV-EXT DRDAV-QR DRDAV-QP DRDAV-QE St. Louis, MO 63166

No. Copies	Addressee	No. Copies	Addressee
7	<p>Commander US Army Tank-Automotive Research and Development Cnd ATTN: DRDTA-UL, Tech Lib DRDTA-RCKM, Mr. S. Goodman DRDTA-RCKT, Mr. J. Fix DRDTA-RTAS, Mr. S. Catalano DRDTA-TTM, Mr. W. Moncrief DRDTA-ZA, Mr. O. Renius (Only Infrared, Ultrasonic, or Holographic Reports) DRDTA-JA, Mr. C. Kedzior Warren, MI 48090</p>	1	<p>Commander US Army Yuma Proving Ground ATTN: Technical Library Yuma, AZ 85364</p>
		1	<p>Commander US Army Tropic Test Center ATTN: STETC-TD, Drawer 942 Fort Clayton, Canal Zone</p>
		4	<p>Commander Aberdeen Proving Ground ATTN: STEAP-MT STEAP-TL STEAP-MT-M, Mr. J. A. Feroi STEAP-MT-G, Mr. R. L. Huddleston Aberdeen Proving Ground, MD 21005</p>
1	<p>Director US Army Industrial Base Engineering Activity ATTN: DRXIB-MT, Mr. D. Brim Rock Island, IL 61299</p>	1	<p>Commander US Army Cold Region Test Center ATTN: STECR-OP-PM APO Seattle, Washington 98733</p>
1	<p>Commander Harry Diamond Laboratories ATTN: DELHD-EDE Mr. B. F. Willis 2800 Powder Mill Road Adelphi, MD 20783</p>	1	<p>Commander US Army Dugway Proving Ground ATTN: STEDP-MT Dugway, UT 84022</p>
2	<p>Commander US Army Test and Evaluation Command ATTN: DRSTE-TD DRSTE-ME Aberdeen Proving Ground, MD 21005</p>	1	<p>Commander US Army Electronic Proving Ground ATTN: STEEP-MT Fort Huachuca, AZ 85613</p>
4	<p>Commander US Army White Sands Missile Range ATTN: STEWS-AD-L STEWS-ID STEWS-TD-PM ATAA-SL (Tech Lib) White Sands Missile Range, NM 88002</p>	1	<p>Commander Jefferson Proving Ground ATTN: STEJP-TD-I Madison, IN 47250</p>
		1	<p>Commander US Army Aircraft Development Test Activity ATTN: STEBG-TD Fort Rucker, AL 36362</p>

No. Copies	Addressee	No. Copies	Addressee
2	President US Army Armor and Engineer Board ATTN: ATZKOE-TA ATZK-AE-PD-E Fort Knox, KY 40121	1	Commander Sacramento Army Depot ATTN: SDSSA-QA Sacramento, CA 95813
1	President US Army Field Artillery Board ATTN: ATZA-BDOP Fort Sill, OK 73503	1	Commander Savanna Army Depot Activity ATTN: SDSSV-S Savanna, IL 61074
1	Commander Anniston Army Depot ATTN: SDSAN-QA Anniston, AL 36202	1	Commander Seneca Army Depot ATTN: SDSEE-R Romulus, NY 14541
1	Commander Corpus Christi Army Depot ATTN: SDSCC-MEE Mr. Haggerty Mail Stop 55 Corpus Christi, TX 78419	1	Commander Sharpe Army Depot ATTN: SDSSH-QE Lathrop, CA 95330
1	Commander Letterkenny Army Depot ATTN: SDSLE-QA Chambersburg, PA 17201	1	Commander Sierra Army Depot ATTN: SDSSI-DQA Herlong, CA 96113
1	Commander Lexington-Bluegrass Army Depot ATTN: SDSRR-QA Lexington, KY 40507	1	Commander Tobyhanna Army Depot ATTN: SDSTO-Q Tobyhanna, PA 18466
1	Commander New Cumberland Army Depot ATTN: SDSNC-QA New Cumberland, PA 17070	1	Commander Tooele Army Depot ATTN: SDSTE-QA Tooele, UT 84074
2	Commander US Army Depot Activity Pueblo ATTN: SDSTE-PU-Q Pueblo, CO 81001	1	Director DARCOM Ammunition Center ATTN: SARAC-DC Savanna, IL 61704
1	Commander Red River Army Depot ATTN: SDSRR-QA Texarkana, TX 75501	2	Naval Research Laboratory ATTN: Dr. J. M. Krafft, Code 8430 Library, Code 2620 Washington, DC 20375

No. Copies	Addressee	No. Copies	Addressee
3	Air Force Materials Laboratory ATTN: AFML-DO, Library AFML-LTM, Mr. E. Wheeler AFML-LLP, Mr. R. Rowand Wright-Patterson AFB, OH 45433	1	HQDA (DAEN-RDL) Washington, DC 20314
		1	HQDA (DAEN-MPE-T) Washington, DC 20314
9	Director Army Materials and Mechanics Research Center ATTN: DRXMR-P (1) DRXMR-PL (2) DRXMR-M (2) DRXMR-MQ (1) DRXMR-MI, Mr. Darcy (1) DRDMM-L, Dr. Chait (1) DRXMR-RA, Mr. Valente (1) DRXMR-AG-MD (1) DRXMR-X (1) DRXMR-PR (1) DRXMR-T (1) DRXMR-E (1) Watertown, MA 02172	1	Director US Army Engineer Waterways Experiment Station ATTN: Chief, Library Branch Technical Information Center Vicksburg, MS 39180
		2	Director Petrol & Fld Svc Dept US Army Quartermaster School Fort Lee, VA 23801
		1	President US Army Aviation Test Board ATTN: STEBG-PO Fort Rucker, AL 36360
1	Director, Technical Information Defense Advanced Research Projects Agency 1400 Wilson Blvd Arlington, VA 22209	1	US Army Aviation School Library P. O. Drawer O Fort Rucker, AL 36360
		2	HQ, 193D Infantry Brigade (Pan) ATTN: AFZU-FE APO Miami 34004
1	Director Defense Nuclear Agency ATTN: TITL Washington, DC 20305	2	Special Forces Detachment, Europe ATTN: PBO APO New York 09050
1	Commander, HQ TRADOC ATTN: ATEN-ME Fort Monroe, VA 23651	2	Engineer Representative USA Research & Standardization Group (Europe) Box 65 FPO 09510
1	HQDA (DAMA-AOA-M) Washington, DC 20310		
1	HQDA (DALO-TSM) Washington, DC 20310	1	HQDA ODCSLOG DALO-TSE Room 1E588 Pentagon, Washington, DC 20310

No. Copies	Addressee	No. Copies	Addressee
1	Plastics Technical Evaluation Center ARRADCOM, Bldg 3401 ATTN: A. M. Anzalone Dover, NJ 07801		Spec Asst/Matl Asmt, DRDME-ZG Spec Asst/Scs & Tech, DRDME-ZK CIRCULATE
1	Commander Frankford Arsenal ATTN: Library, K2400, B151-2 Philadelphia, PA 19137	1	C, Ctrmine Lab, DRDME-N C, Engy & Wtr Res Lab, DRDME-G C, Elec Pwr Lab, DRDME-E C, Camo & Topo Lab, DRDME-R C, CS/CI Lab, DRDME-M C, Mech & Constr Eqpt Lab, DRDME-H C, Ctr Intrus Lab, DRDME-X C, Matl Tech Lab, DRDME-V Dir, Prod A&T Dir, DRDME-T CIRCULATE
1	Commandant US Army Engineer School ATTN: ATZA-CDD Fort Belvoir, VA 22060		
1	President US Army Airborne, Communications and Electronics ATTN: STEBF-ABTD Fort Bragg, NC 28307	5	Material Technology Lab, DRDME-V
		10	Organic Chem Ctg Rsch Grp, DRDME-VO
		3	Tech Reports Ofc, DRDME-WP
		3	Security Ofc (for liaison officers), DRDME-S
		2	Tech Library, DRDME-WC
1	Commander Headquarters, 39th Engineer Battalion (Cbt) Fort Devens, MA 01433	1	Programs & Anal Dir, DRDME-U
		1	Pub Affairs Ofc, DRDME-I
		1	Office of Chief Counsel, DRDME-L
1	Commander and Director USA FESA ATTN: FESA-TS Fort Belvoir, VA 22060	1	Director, Physics Program (421) Office of Naval Research Arlington, VA 22217
		2	Commander, Naval Facilities Engineering Command Department of the Navy ATTN: Code 032-B, 062 200 Stovall St Alexandria, VA 22332
1	HQ, USAEUR & Seventh Army Deputy Chief of Staff, Engineer ATTN: AEAEN-MT-P APO New York 09403		
1	HQ, USAEUR & Seventh Army Deputy Chief of Staff, Operations ATTN: AEAGC-FMD APO New York 09403	1	US Naval Oceanographic Office Navy Library/NSTL Station Bay St. Louis, MS 39522
	MERADCOM	1	Library (Code L08A) Civil Engineering Laboratory Naval Construction Battalion Ctr Port Hueneme, CA 93043
1	Commander, DRDME-Z Tech Dir, DRDME-ZT Assoc Tech Dir/R&D, DRDME-ZN Assoc Tech Dir/Engrg & Acq, DRDME-ZE		

No. Copies	Addressee
1	Director Earth Physics Program Code 464 Office of Naval Research Arlington, VA 22217
1	Naval Training Equipment Center ATTN: Tech Library Orlando, FL 32813
1	HQ USAF/RDPT ATTN: Mr. Allan Eaffy Washington, DC 20330
1	HQ USAF/LEEEU Chief, Utilities Branch Washington, DC 20330
1	US Air Force HQ Air Force Engineering and Services Center Technical Library FL 7050 Tyndall AFB, FL 32403
1	Chief, Lubrication Br Fuels & Lubrication Div ATTN: AFWAL/POSL Wright-Patterson AFB, OH 45433
1	Department of Transportation Library, FOP 10A, M494-6 800 Independence Ave, SW Washington, DC 20591
1	Professor Raymond R. Fox School of Engineering and Applied Science George Washington University Washington, DC 20052
1	Reliability Analysis Center Rome Air Development Center RADC/RBRAC (I. L. Krulac) Griffiss AFB, NY 13441