



MICROCOPY RESOLUTION TEST CHART NATIONAL BUREAU OF STANDARDS-1963-A

DTIC ELECTE MAY 1 7 1984

18. SUPPLEMENTARY NOTES

19. KEY WORDS (Continue on reverse side if necessary and identify by block number)

Electrical Noise; Noise Measurements; Fast Fourier Transform; Digital Noise Measurements.

20. ABSTRACT (Continue on reverse side if necessary and identify by block number)

The IQS 401 spectrum analyzer, a commercial addition to the Apple II & III personal computers, can be modified for low-cost power spectrum measurements of random noise signals and similar waveforms. The IQS BASIC control program is changed to perform low-frequency sampling, spectral averaging and spectrum scaling. External analog filtering must be provided on the lowest sampling rates; however, no change is made to the original IQS hardware. Resolution of 0.5 dB is obtainable from 10% Hz to 20 kHz.

0.0001

DD 1 JAN 73 1473 EDITION OF 1 NOV 65 IS OBSOLETE

SECURITY CLASSIFICATION OF THIS PAGE (When Date Entered)

84 05 14 067

TE FILE CO.

OFFICE OF NAVAL RESEARCH Contract #NOO014-82-K-0603 TECHNICAL REPORT NO. 4

SPECTRAL AVERAGING AND LOW-FREQUENCY SAMPLING MODIFICATIONS

FOR THE

IQS 401 FFT SPECTRUM ANALYZER

by

Steven W. Smith

Prepared for Publication

in the

Review of Scientific Instruments

Department of Physics University of Utah Salt Lake City, Utah 84112

May, 1983

Reproduction in whole or in part is permitted for any purpose of the United States Government

This document has been approved for public release and sale; its distribution is unlimited

SPECTRAL AVERAGING AND LOW-FREQUENCY SAMPLING MODIFICATIONS FOR THE

IQS 401 FFT SPECTRUM ANALYZER

by

Steven W. Smith
Physics Department
University of Utah
Salt Lake City, Utah 84112

ABSTRACT

The IQS 401 spectrum analyzer, a commercial addition to the Apple II and III personal computers, can be modified for low-cost power spectrum measurements of random noise signals and similar waveforms. The IQS BASIC control program is changed to perform low-frequency sampling, spectral averaging and spectrum scaling. External analog filtering must be provided on the lowest sampling rates; however, no change is made to the original IQS hardware. Resolution of 0.5 dB is obtainable from 10^{-4} Hz to 20 kHz.

i

SPECTRAL AVERAGING AND LOW-TREQUENCY SAMPLING MODIFICATIONS FOR THE IQS 401 FFT SPECTRUM ANALYZER

The IQS 401 FFT spectrum analyzer¹ is a commercial addition to the Apple II and III personal computers. A versatile plug-in hardware card and BASIC control program perform data acquisition, waveform display and spectral analysis. The 401 is well suited for the measurement of repeatable signals such as loudspeaker impulse response testing; however, it has several deficiencies for the analysis of random signals as encountered in electrical noise studies and similar areas. These deficiencies can be corrected to yield an inexpensive power spectrum analyzer capable of measurements extending from 10⁻⁵ Hz to 20 kHz with 0.5 dB resolution. Approximately 30% of the IQS BASIC control program is changed and external analog filtering is provided. No change is made to the original IQS hardware.

These software modifications enhance the 401's capabilities in three areas: spectral averaging, power spectra scaling and low-frequency sampling.

The unmodified 401 estimates power spectra by computing the Fast Fourier Transform using 128, 256, 512 or 1024 points. The uncertainty in this power spectrum estimate is large; however, it can be reduced by smoothing the spectrum with a moving average digital filter provided in the 401 software. This unavoidably removes the fine details of the spectrum

and may not reduce the uncertainty to an acceptable level. Averaging several independent power spectrum estimates reduces the uncertainty without removing the fine details; however, the original 401 has no provision for this spectral averaging. The software modifications make it possible to rms average 1 to 100 separate spectra. A total of 64 averages results in 0.5 dB resolution when measuring Nyquist noise.

The unmodified 401 displays relative power spectra only. The software modifications calculate a scale factor and adjust the display to correspond to the units of Volts squared per Hertz. This procedure accounts for all related system parameters including sample rate, sample length and amplification.

The unmodified 401 samples at a selectable frequency between 200 Hz and 60 kHz, resulting in a lowest analyzable frequency of 0.2 Hz. The software modifications control the IQS hardware to allow for lower sampling rates. The modifications add seven sampling rates between 10 Hz and 0.01 Hz, resulting in a lowest analyzable frequency of 10^{-5} Hz. Each low-frequency sample is displayed on the computer console screen as it is obtained. This provides early detection of out-of-range errors and other low-frequency data acquisition problems.

The low sampling rates require an analog anti-alias filter external to the computer. An external high-pass filter proves useful in removing DC drift with the low sampling rates. The system shown in Figure 1 has been used for noise studies to 10^{-4} Hz. The analog filters are standard

designs² constructed from operational amplifiers. Noise measurements between 0.2 Hz and 20 kHz require only the Princeton Applied Research PAR-113 low noise preamplifier or equivalent.

Figure 2 shows typical power spectra obtained with the system. The input voltage spectral noise density of the PAR-113 and LM-394 preamplifiers is shown together with the observed Nyquist noise of two resistors. This data was obtained by averaging 64 spectra of 256 points at sampling rates of 46 kHz, 1 kHz, 10 Hz and 0.1 Hz. A single sampling rate with 256 points spans about two decades of frequency as shown by the open and filled dated points.

The modified 401's usefulness is limited in some applications by its slowness and limited dynamic range. Approximately 30 minutes is required to calculate a power spectrum from 64 averages of 256 point spectra. The calculation time is limited by the relatively slow 6502 processor used in the Apple II computer. Eight bit analog to digital conversion in the IQS hardware limits the spectral dynamic range to 32 dB.

The ability to adapt this system through software changes makes the 401 an attractive alternative to dedicated power spectrum analyzers costing several times more. Software modifications and filter schematics are available from the author on request. This work was supported in part by the Office of Naval Research.

REFERENCES

- 1. IQS Inc. 5719 Corso Di Napoli, Long Beach, California 90803
- 2. P. Horwitz and W. Hill: "The Art of Electronics", Cambridge University Press, New York, 1980.
- S. W. Smith: "Internal Noise of Low-Frequency Preamplifiers", Review of Scientific Instruments, May 1984.

FIGURE CAPTIONS

- Fig. 1 Typical external filtering for low-frequency sampling.
- Fig. 2 Experimental measurements of amplifier noise and Nyquist noise using modified FFT system.


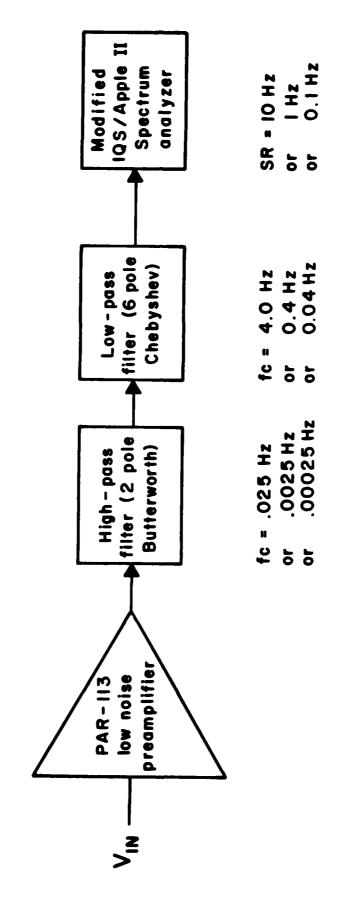


Figure 1. Typical external filtering for low-frequency sampling

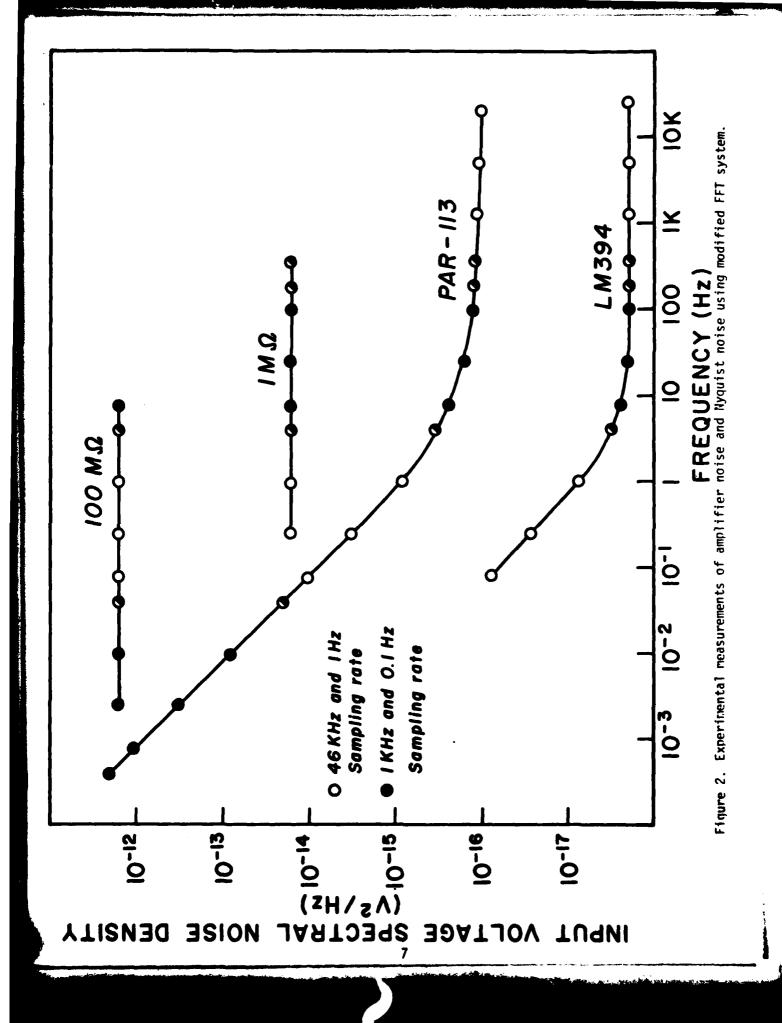

	No. Copies		No. Copies
Office of Naval Research		Naval Ocean Systems Center	
Attn: Code 413		Attn: Mr. Joe McCartney	
800 North Quincy Street		San Diego, California 92152	1
Arlington, Virginia 22217	2	- .	
		Naval Weapons Center	
ONR Pasadena Detachment		Attn: Dr. A. B. Amster,	
Attn: Dr. R. J. Marcus		Chemistry Division	
1030 East Green Street		China Lake, California 93555	, 1
Pasadena, California 91106	1	•	
		Naval Civil Engineering Laboratory	
Commander, Naval Air Systems Command		Attn: Dr. R. W. Drisko	
Attn: Code 310C (H. Rosenwasser)		Port Hueneme, California 93401	1
Department of the Navy		•	
Washington, D.C. 20360	1	Dean William Tolles	
•		Naval Postgraduate School	
Defense Technical Information Center		Monterey, California 93940	1
Building 5, Cameron Station		• :	
Alexandria, Virginia 22314	12	Scientific Advisor	
		Commandant of the Marine Corps	
Dr. Fred Saalfeld		(Code RD-1)	
Chemistry Division, Code 6100		Washington, D.C. 20380	1
Naval Research Laboratory			
-Washington, D.C. 20375	1	Naval Ship Research and Development	
•		Center	
U.S. Army Research Office		Attn: Dr. G. Bosmajian, Applied	
Attn: CRD-AA-IP		Chemistry Division	
P. O. Box 12211		Annapolis, Maryland 21401	1
Research Triangle Park, N.C. 27709	1	•	
•		Mr. John Boyle	
Mr. Vincent Schaper		Materials Branch	
DTNSRDC Code 2803		Naval Ship Engineering Center	
Annapolis, Maryland 21402	1	Philadelphia, Pennsylvania 19112	1
Naval Ocean Systems Center		Mr. A. M. Anzalone	
Attn: Dr. S. Yamamoto		Administrative Librarian	
Marine Sciences Division		PLASTEC/ARRADCOM	
San Diego, California 91232	1	B1dg 3401	
		Dover, New Jersey 07801	1

•	No. Copies	•	No. Copies
Dr. H. Richtol Chemistry Department Rensselser Polytechnic Institute		Dr. R. P. Van Duyne Department of Chemistry	
Troy, New York 12181	• 1	Northwestern University Evanston, Illinois 60201	1
Dr. A. B. Ellis Chemistry Department University of Wisconsin		Dr. B. Stanley Pons Department of Chemistry	
Madison, Wisconsin 53706 Dr. M. Wrighton	1	University of Utah Henry Eyring Building Salt Lake City, Utah 84112	1
Chemistry Department Massachusetts Institute of Technology Cambridge, Massachusetts 02139		Dr. Michael J. Weaver Department of Chemistry Purdue University West Lafayette, Indiana 47097	1 .
Larry E. Plew		Dr. R. David Rauh	•
Naval Weapons Support Center		EIC Corporation	
Code 30736, Building 2906		55 Chapel Street	
Crane, Indiana 47522	1	Newton, Massachusetts 02158	1
S. Ruby		Dr. J. David Margerum	
DOE (STOR)		Research Laboratories Division	
600 E Street		Hughes Aircraft Company	
Providence, Rhode Island 02192	1	3011 Malibu Canyon Road	_
Dr. Aaron Wold	• .	Malibu, California 90265	1
Brown University		Dr. Verste Flateshaue	
Department of Chemistry		Dr. Martin Fleischmann	
Providence, Rhode Island 02192	1 .	Department of Chemistry University of Southampton	
The state of the s	•	Southampton 509 5NH England	1
Dr. R. C. Chudacek		podenampeda 303 3km bagiand	
McGraw-Edison Company		Dr. Janet Osteryoung	
Edison Battery Division		Department of Chemistry	
Post Office Box 28		State University of	
Bloomfield, New Jersey 07003	1	New York at Buffalo	
	•	Buffalo, New York 14214	1
Dr. A. J. Bard			
University of Texas		Dr. R. A. Osteryoung	
Department of Chemistry		Department of Chemistry	
Austin, Texas 78712	. 1	State University of New York at Buffalo	
Dr. M. M. Nicholson		Buffalo, New York 14214	1
Electronics Research Center			
Rockwell International			
3370 Miraloma Avenue			
Anaheim, California 91752	1		•

No. Copies

Dr. Paul Delahay		Dr. P. J. Hendra
Department of Chemistry		Department of Chemistry
New York University		University of Southampton
New York, New York 10003	1	Southampton SOO 5NH
•		United Kingdom
Dr. E. Yeager		•
Department of Chemistry		Dr. Sam Perone
Case Western Reserve University		Chemistry & Materials
Cleveland, Ohio 41106	1	Science Department
	_	Laurence Livermore National Lab.
Dr. D. N. Bennion		Livermore, California 94550
Department of Chemical Engineering		
Brigham Young University		Dr. Royce W. Murray
Provo, Utah 84602	1	Department of Chemistry
	•	University of North Carolina
Dr. R. A. Marcus		Chapel Hill, North Carolina 275
Department of Chemistry		onopol male, more officers and
California Institute of Technology		Naval Ocean Systems Center
Pasadena, California 91125	1	Attn: Technical Library
resevena, valliolina 71125	•	San Diego, California 92152
Dr. J. J. Auborn		Dan Diego, Calliothia 72131
Bell Laboratories		Dr. C. E. Mueller
*Murray Hill, New Jersey 07974	1	The Electrochemistry Branch
nullay nill, new Jetsey 0/7/4	1	Materials Division, Research and
Dr. Adam Heller		Technology Department
Bell Laboratories		Naval Surface Weapons Center
	•	White Oak Laboratory
Murray Hill, New Jersey 07974	1	Silver Spring, Maryland 20910
Dm. M. Vahar		Silver Spring, Maryland 20910
Dr. T. Katan		Dr. G. Goodman
Lockheed Missiles and		_
Space Co., Inc.		Johnson Controls
P. O. Box 504	•	5757 North Green Bay Avenue
Sunnyvale, California 94088	1	Milwaukee, Wisconsin 53201
D		Br. J. Boocklon
Dr. Joseph Singer, Code 302-1		Dr. J. Boechler
NASA-Lewis		Electrochimica Corporation
21000 Brookpark Road	•	Attn: Technical Library
Cleveland, Ohio 44135	1	2485 Charleston Road
	•	Mountain View, California 9404
Dr. B. Brummer		
EIC Incorporated		Dr. P. P. Schmidt
55 Chapel Street	_	Department of Chemistry
Newton, Massachusetts 02158	1	Oakland University
	•	Rochester, Michigan 48063
Library		
P. R. Mallory and Company, Inc.		
Northwest Industrial Park		
Burlington, Massachusetts 01803	1	

	No. Copies		No. Copies
Dr. Donald W. Ernst		Mr. James R. Moden	
Naval Surface Weapons Center		Naval Underwater Systems	
Code R-33		Center	•
White Oak Laboratory		Code 3632	
Silver Spring, Maryland 20910	i	Newport, Rhode Island 02840	1
Dr. R. Nowak		Dr. Bernard Spielvogel	
Naval Research Laboratory		U. S. Army Research Office	
Code 6130 .		P. O. Box 12211	
Washington, D.C. 20375	. 1	Research Triangle Park, NC 27709	1
Dr. John F. Houlihan		Dr. Denton Elliott	
Shenango Valley Campus		Air Force Office of	
Pennsylvania State University		Scientific Research	
Sharon, Pennsylvania 16146	1	Bolling AFB	
À		Washington, D.C. 20332	1
Dr. D. F. Shriver		•	_
Department of Chemistry		Dr. David Aikens	
Northwestern University		Chemistry Department	
Evanston, Illinois 60201	1	Rensselaer Polytechnic Institute Troy, New York 12181	1
Dr. D. H. Whitmore		110y, New 101k 12101	1
Department of Materials Science		Dr. A. P. B. Lever	
Northwestern University		Chemistry Department	
Evanston, Illinois 60201	1	York University	
L'auston, lilinois dozdi		Downsview, Ontario M3J1P3	
Dr. Alan Bewick		Canada	•
Department of Chemistry		Callada	1
The University		No Charlelau Crack	
Southampton, S09 5NH England	1	Dr. Stanislaw Szpak	
Sociampion, 509 Jan England	•	Naval Ocean Systems Center Code 6343	
Dr. A. Himy		San Diego, California 95152	1
NAVSEA-5433		•	
NC #4		Dr. Gregory Farrington	
2541 Jefferson Davis Highway	_	Department of Materials Science	
Arlington, Virginia 20362	1	and Engineering	
No. John Wassadd		University of Pennsylvania	•
Dr. John Kincaid		Philadelphia, Pennsylvania 19104	1
Department of the Navy			
Strategic Systems Project Office		Dr. Bruce Dunn	
Room 901	1	Department of Engineering &	
Washington, D.C. 20376	•	Applied Science	•
		University of California	1
		Los Angeles, California 90024	

	No. Copies		No. Copies
M. L. Robertson		Dr. T. Marks	
Manager, Electrochemical		Department of Chemistry	
and Power Sonices Division		Northwestern University	
Naval Weapons Support Center		Evanston, Illinois 60201	1
Crane, Indiana 47522	1		•
		Dr. D. Cipris	•
Dr. Elton Cairns		Allied Corporation	
Energy & Environment Division		P. O. Box 3000R	
Lawrence Berkeley Laboratory		Morristown, New Jersey 07960	1
University of California			
Berkeley, California 94720	1 .	Dr. M. Philpot	
		IBM Corporation	
Dr. Micha Tomkiewicz		5600 Cottle Road	
Department of Physics		San Jose, California 95193	1
Brooklyn College			
Brooklyn, New York 11210	1	Dr. Donald Sandstrom	
		Washington State University	
Dr. Lesser Blum		Department of Physics	_
Department of Physics		Pullman, Washington 99164	1
University of Puerto Rico	•		
Rio Piedras, Puerto Rico 00931	j	Dr. Carl Kannewurf Northwestern University	
Dr. Joseph Gordon, II		Department of Electrical Engineering	;
IBM Corporation		and Computer Science	
K33/281		Evanston, Illinois 60201	1
5600 Cottle Road			
San Jose, California 95193	1	Dr. Edward Fletcher	
		University of Minnesota	
Dr. Robert Somoano		Department of Mechanical Engineering	
Jet Propulsion Laboratory		Minneapolis, Minnesota 55455	1
California Institute of Technology	_	m	
Pasadena, California 91103	1	Dr. John Fontanella	
De Johann A Josh-Al		U.S. Naval Academy	
Dr. Johann A. Joebstl		Department of Physics	•
USA Mobility Equipment R&D Command DRDME-EC	• •	Annapolis, Maryland 21402	1
Fort Belvior, Virginia 22060	1	Dr. Martha Greenblatt	
		Rutgers University	
Dr. Judith H. Ambrus		Department of Chemistry	_
NASA Headquarters M.S. RTS-6		New Brunswick, New Jersey 08903	1
Washington, D.C. 20546	1	Dr. John Wassib Kings Mountain Specialties	
Dr. Albert R. Landgrebe		P. O. Box 1173	
U.S. Department of Energy		Kings Mountain, North Carolina 2808	6 1
M.S. 6B025 Forrestal Building		•	
Washington, D.C. 20595	1		

	No. Copies		No. Copies
Dr. J. J. Brophy University of Utah Department of Physics Salt Lake City, Utah 84112	· 1		
Dr. Walter Roth Department of Physics State University of New York Albany, New York 12222	1		
Dr. Thomas Davis National Bureau of Standards Polymer Science and Standards Division Washington, D.C. 20234	1	^	
Dr. Charles Martin Department of Chemistry Texas A&M University College Station, Texas 77843	1		
Dr. Anthony Sammells Institute of Gas Technology 3424 South State Street Chicago, Illinois 60616	1		-
Dr. H. Tachikawa Department of Chemistry Jackson State University Jackson, Mississippi 39217	1		
Dr. W. M. Risen Department of Chemistry Brown University Providence, Rhode Island 02912	1		

