

ND-A188 539

QUALITY ASSURANCE/QUALITY CONTROL IN WASTE SITE  
CHARACTERIZATION AND REMEDIAL ACTION(U) OAK RIDGE  
NATIONAL LAB TN H P MASKARINEC ET AL. 1987


1/1


UNCLASSIFIED

OWRL/TH-10600 NIPR-D-6-78-14

F/G 24/3

ML

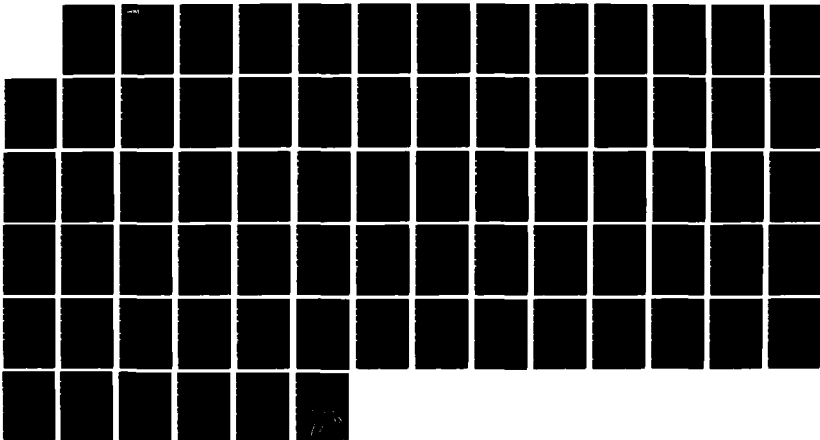

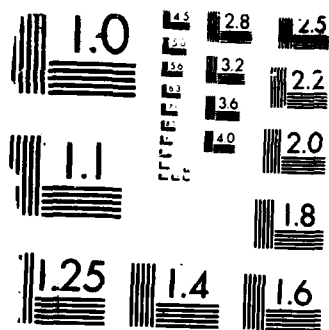

MICROCOPY RESOLUTION TEST CHART

**ornl**


**AD-A188 539**

AD \_\_\_\_\_  
ORNL/TM-10600

**OAK RIDGE  
NATIONAL  
LABORATORY**

**MARTIN MARIETTA**

**Quality Assurance/Quality Control  
in Waste Site Characterization  
and Remedial Action**

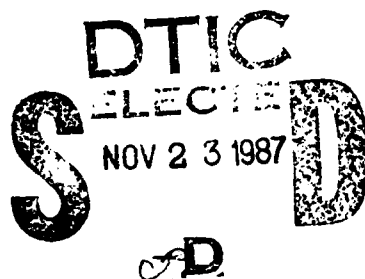

**Final Report**

M. P. Maskarinec  
S. K. Holladay

**Supported by**

U.S. Army Toxic and Hazardous Materials Agency  
Aberdeen Proving Ground, Maryland 21010-5401

**Project Officer: Mary Ann Ryan**

**Approved for public release;  
distribution unlimited**

The findings in this report are not to be construed  
as an official Department of the Army position unless  
so designated by other authorized documents.

OPERATED BY  
MARTIN MARIETTA ENERGY SYSTEMS, INC.  
FOR THE UNITED STATES  
DEPARTMENT OF ENERGY

62 2 1 004

*Approved*

REPORT DOCUMENTATION PAGE				
1a. REPORT SECURITY CLASSIFICATION		1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY		3. DISTRIBUTION/AVAILABILITY OF REPORT Distribution unlimited. Cleared for public release.		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE				
4. PERFORMING ORGANIZATION REPORT NUMBER(S) ORNL/TM-10600		5. MONITORING ORGANIZATION REPORT NUMBER(S)		
6a. NAME OF PERFORMING ORGANIZATION Oak Ridge National Laboratory	6b. OFFICE SYMBOL (if applicable)	7a. NAME OF MONITORING ORGANIZATION USATHAMA		
6c. ADDRESS (City, State, and ZIP Code) P. O. Box X Oak Ridge, TN 37831-6120		7b. ADDRESS (City, State, and ZIP Code) AMXTH-TE-A Aberdeen Proving Ground, MD 21010-5401		
8a. NAME OF FUNDING/SPONSORING ORGANIZATION USATHAMA	8b. OFFICE SYMBOL (if applicable)	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER MIPR #D-6-78-14		
8c. ADDRESS (City, State, and ZIP Code) AMXTH-TE-A Aberdeen Proving Ground, MD 21010-5401		10. SOURCE OF FUNDING NUMBERS		
		PROGRAM ELEMENT NO.	PROJECT NO.	TASK NO.
		WORK UNIT ACCESSION NO.		
11. TITLE (Include Security Classification) Quality Assurance/Quality Control in Waste Site Characterization and Remedial Action				
12. PERSONAL AUTHOR(S) M. P. Maskarinec and S. K. Holladay				
13a. TYPE OF REPORT FINAL	13b. TIME COVERED FROM 12/86 TO 6/87	14. DATE OF REPORT (Year, Month, Day)	15. PAGE COUNT 70	
16. SUPPLEMENTARY NOTATION				
17. COSATI CODES		18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number) IR QA Plan CLP Interagency Task Force, Remedial Action <i>Hazardous Materials</i>		
FIELD	GROUP			
19. ABSTRACT (Continue on reverse if necessary and identify by block number) In this report, efforts to advance the state of the art of Quality Assurance/Quality Control in waste site characterization and remedial action are detailed. The two most widely used and accepted QA/QC programs, the USEPA CLP and the USATHAMA IR QA plan, are compared and recommendations are provided for unification. The findings of a Working Group convened to discuss these issues are reported and the formation of a Task Force on Quality Assurance/Quality Control is announced. <i>Keynote</i>				
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input type="checkbox"/> UNCLASSIFIED/UNLIMITED <input checked="" type="checkbox"/> SAME AS RPT <input type="checkbox"/> DTIC USERS		21. ABSTRACT SECURITY CLASSIFICATION		
22a. NAME OF RESPONSIBLE INDIVIDUAL Mary Ann Ryan		22b. TELEPHONE (Include Area Code) (301) 671-3206	22c. OFFICE SYMBOL	

AD  
ORNL/TM-10600

Quality Assurance/Quality Control in Waste Site  
Characterization and Remedial Action

M. P. Maskarinec and S. K. Holladay

Analytical Chemistry Division  
Oak Ridge National Laboratory  
Oak Ridge, Tennessee 37831-6120

SUPPORTED BY

U.S. ARMY TOXIC AND HAZARDOUS MATERIALS AGENCY  
Aberdeen Proving Ground, Maryland 21010-5401

Project Officer: Mary Ann Ryan

Date Published - October 1987

OAK RIDGE NATIONAL LABORATORY  
Oak Ridge, Tennessee 37831  
operated by  
MARTIN MARIETTA ENERGY SYSTEMS, INC.  
for the  
U.S. DEPARTMENT OF ENERGY  
Under Contract No. DE-AC05-84OR21400

1  
INSPECTED  
EDPA  
DITC

Accession For	
NTIS CRM&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution /	
Availability Codes	
Dist	Avail. and/or Control
A-1	

## EXECUTIVE SUMMARY

This report details efforts to date to advance the state of the art of Quality Assurance/Quality Control (QA/QC) in waste site characterization and remedial action. The report is organized into three sections. The first section provides recommendations on the unification of the two widely used and accepted QA/QC programs: the U.S. Environmental Protection Agency Contract Laboratory Program (USEPA CLP) and the U.S. Army Toxic and Hazardous Materials Agency Installation Restoration Quality Assurance (USATHAMA IR QA) plan. The second section compares the two plans in detail with the Guidelines given by the USEPA. The third section announces the formation of a Task force on Quality Assurance/Quality Control and reports on the findings of a Working Group convened to address these issues.

## TABLE OF CONTENTS

	<u>Page</u>
EXECUTIVE SUMMARY . . . . .	1
UNIFICATION OF THE USATHAMA IR QA PLAN WITH THE USEPA CONTRACT LABORATORY PROGRAM . . . . .	5
REVIEW OF THE USATHAMA QA PROGRAM (MARCH, 1987) AND THE USEPA CONTRACT LABORATORY PROGRAM USING THE USEPA SIXTEEN-POINT QA PROJECT PLAN . . . . .	9
Appendix A . . . . .	39
Appendix B . . . . .	41
Appendix C . . . . .	43
Appendix D . . . . .	54
DEVELOPMENT OF THE QUALITY ASSURANCE TASK FORCE . . . . .	63
DISTRIBUTION LIST . . . . .	71

UNIFICATION OF THE USATHAMA IR QA PLAN WITH THE USEPA  
CONTRACT LABORATORY PROGRAM

There currently exists a widespread agreement on the need for remedial action at past waste disposal sites. The approach usually taken is to study records pertaining to the site (preliminary investigation), to follow the investigation with a survey of contamination, and then to decide on a remedial action plan. Privately used sites are regulated under the Superfund Amendment and Reauthorization Act (SARA) by the USEPA. Sites used by government agencies, such as DOD and DOE, while regulated under SARA, are not generally cleaned up using SARA funds. An important aspect of the entire process is the analysis of large numbers of samples. Because of the increased emphasis on analytical methodologies, and the associated cost, it is crucial to ensure that the data produced be of acceptable quality. Therefore, strict Quality Assurance (QA) measures must be applied. Several different approaches to the various aspects of QA have been developed over the last decade, with perhaps the best known being the approach used by the USEPA under the Contract Laboratory Program (CLP). In addition, the U.S. Army Toxic and Hazardous Materials Agency (USATHAMA) has developed a QA plan to serve the needs of the Army Installation Restoration program. Both systems have been in existence for several years, and substantial experience has been gained. Because of the obvious similarity in the objectives of the two systems, this work was performed in an attempt to draw the two plans together. This was done in the following manner. A workshop was held to bring all interested parties together and discuss differences and similarities. Then, a detailed comparison was made of the two plans. Finally, this report was written to document the results of these efforts.

It must be noted that certain differences exist in the two programs which result from philosophical and logistical considerations beyond the issue of QA. Virtually all data generated in the CLP may eventually be called into court as evidence for prosecution or cost recovery actions. Therefore, it is necessary that the data be of courtroom quality. In the case of the IR plan, use of the data is generally restricted to direction of the remedial action phase. While not removing the requirement for high quality data, this end use does not mandate the degree of documentation used in the CLP. More important in the IR plan is streamlining of the data flow, and rapid identification of QA problems. In addition, the CLP has a relatively large number of participating laboratories compared to the IR system. This, combined with the end use difference, results in a need for more rigid standardization of the entire QA process. Interlaboratory comparability becomes much more crucial to the CLP than to the IR program. A related logistical difference is the fact that in the CLP, samples are collected and sent to a sample management facility, either central or regional. The samples are then distributed to the analytical facilities. In the USATHAMA case, samples are collected by



the prime contractor for the remedial action and either analyzed in house or sent to a subcontracting analytical facility. The CLP can therefore do a more effective job of providing blind QA samples (spikes, splits, or blanks) than can the IR plan. Taken together, these differences result in a marked reliance on external QA in the CLP, and a corresponding reliance on internal QA in the IR plan. Given this difference, it would appear that the two systems are not mutually exclusive, and that reconciliation of the programs might result in an even stronger unified plan.

The first issue to be resolved is the method of assuring initial laboratory proficiency. In the CLP, this is done by analysis of a performance evaluation sample. In the IR plan, the laboratory performs a certification study, which is used to establish the QA parameters for the method. While not specifically required by the CLP, some type of initial certification must certainly be done by the laboratory in order to gain familiarity with the method prior to running the PE sample. It would seem prudent to establish guidelines for the certification process which allow the laboratory to prove competence in the method. The certification procedure used by USATHAMA should be recommended or even required by the CLP.

The use of performance evaluation samples has advantages and disadvantages. PE samples can rarely be provided which are truly blind. There is no assurance that the successful analysis of a PE sample reflects everyday laboratory performance. Furthermore, the analysis of PE samples is restricted to a relatively low frequency (quarterly) so that if problems are identified, large gaps exist in which laboratory performance is in question. Finally, the time required to analyze the PE data further increases the lag time. On the other hand, the PE sample is the only truly external check on laboratory performance, and the only means by which laboratories can be compared and rated. Therefore, it is recommended that USATHAMA adopt the performance evaluation system used by the CLP, and that efforts be made to rapidly evaluate that data and report problems to the laboratory.

Several differences exist between the two plans with respect to sampling and analysis. The only fundamental difference is that the CLP requires the analysis of all samples for compounds on the Hazardous Substances List (HSL) while the Army has contamination from military-specific compounds which do not appear on the HSL. When the IR plan is used for analysis of HSL compounds, the CLP methodology is followed. Therefore, the IR plan is equivalent from an analytical standpoint to the CLP, but includes in addition the QA required for non-CLP methods. Other differences include container cleaning procedures and holding times. Differences of this type can be handled experimentally, by performing an equivalency test. The IR plan does not require chain-of-custody procedures to be followed, unless the data is to be used in litigation. When used, these procedures are functionally identical to

those used in the CLP plan. Chain-of-custody procedures are a necessary part of good laboratory practice, and should always be used. It is recommended that USATHAMA require CLP chain-of-custody procedures be followed for all samples.

In terms of data management and communication, USATHAMA has developed a sophisticated computer-based system. All data is entered by the analyst into a personal computer, checked for completeness, and transferred to a mainframe. The laboratory is required to submit all raw data at the end of the contract. In the CLP, all data generated pertaining to a particular sample is submitted with the results from that sample in a data package. This is clearly an example of differences resulting from the end-use situation mentioned previously. The USATHAMA system is far more workable from the standpoint of remedial action decision making, but the CLP system is required for litigation. However, it must be pointed out that the data in either case is available, and that nothing has been lost. Therefore, it should be possible for the IR software to produce a CLP data package on demand. If this can be done, then there is no practical reason to change either plan. It is recommended that USATHAMA demonstrate the ability to produce a CLP data package.

The software package used in the IR plan has additional features which are quite desirable from a QA standpoint, including the ability to generate QC charts. While it is implied in the CLP that QC charts should be kept, no formal requirement exists and no standardized approach is provided. QC charting has several advantages: rapid identification of out-of-control situations, assurance that performance is consistent on a day-to-day basis, and documentation that the laboratory is performing well on each and every sample. Thus, QC charting can serve as an adjunct to the PE system, and alleviate the drawbacks of PE samples. The question is: what should be charted? Since the surrogates and internal standards used in the CLP are present in every sample, it seems logical to require that the surrogate recoveries and internal standard areas be control charted. It is recommended that the CLP use the USATHAMA software package and require control charts for surrogate recoveries and internal standard areas. It is further recommended that USATHAMA provide USEPA with the software and documentation.

A difference also exists in the area of matrix spiking. The CLP requires a matrix spike and matrix spike duplicate to be run for each matrix. The IR plan uses a standard matrix. The CLP matrix spike does provide additional information on the performance of the methods with respect to individual matrices. However, it can be difficult to determine when one matrix differs from the previous one. On the other hand, the IR method provides a historical record of the performance of the method with time. Given that the surrogates are present in every sample and can be considered matrix spikes, the issue seems to be whether any additional information can be obtained from sample matrix

spikes. Furthermore, the issue of interlaboratory comparability - so important to the CLP program - would be better served by the use of a standard matrix than by use of sample matrices. It is recommended that the CLP drop the requirement for sample matrix spikes and matrix spike duplicates and adopt the standard matrix approach used by USATHAMA.

One of the major problems faced by analytical laboratories doing work in the remedial action area is the audit. Each contracting agency has its own style of auditing, and preparation for the audit depends on the needs and requirements of the auditing agency. In both the CLP and the IR programs, the audit is used as a tool to improve the performance of the laboratory. Because of all of the differences listed previously, the audits take on a different flavor depending on which agency is auditing. However, if the modifications recommended in this document can all be made, the audit could be performed by either USATHAMA or CLP personnel and would suit the needs of both programs. This would result in substantial savings to the agencies involved and would be very convenient for the laboratories.

In summary, the recommendations made here are the result of an objective comparison of the CLP and the USATHAMA IR QA plan. These recommendations are made with the goal of improving quality assurance and quality control in environmental measurements related to waste site characterization and remedial action. An additional goal is the reduction of the cost associated with QA. Two approaches are feasible in this regard. The most easily adopted from the philosophical viewpoint is the declaration of equivalency of the two plans. To this end, a detailed comparison follows of the two plans with the general guidelines set forth in the USEPA sixteen point QA project plan. While this would be expedient, the separate-but-equal approach is far less desirable than the approach of combining the best of both. To that end, continuing communication between the principal agencies and a willingness to cooperate on these issues is mandatory. It is recommended that the USEPA grant equivalency to the USATHAMA IR QA plan, but at the same time strive for unification.

REVIEW OF THE USATHAMA QA PROGRAM (MARCH, 1987) AND THE USEPA CONTRACT  
LABORATORY PROGRAM USING THE USEPA SIXTEEN POINT QA PROJECT PLAN

1. Title Page with Provision for Approval Signatures

2. Table of Contents

3. Project Description

4. Project Organization and Responsibility

A diagram of the lines of communication for USATHAMA IR projects (USATHAMA QA Program, 2nd edition, March, 1987) has been included in Appendix A and a diagram of the program principals of the USEPA (User's Guide to the Contract Laboratory Program, October, 1984) in Appendix B.

5 QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability

Laboratory Certification		EPA-CLP	COMMENTS
USATHAMA			
Contract award	Laboratory selection		USATHAMA specifies that the QA/QC plan must be documented and in practice before samples arrive. EPA-CLP requires a documented QA/QC SOP, but does not specify a timetable.
Development of Project QC Plan	Development of written QA/QC standard operating procedure (SOP) containing these essential elements:		
A statement of adherence or reference to the USATHAMA QA Program	Organization and Personnel QA policy and objectives QA organization Personnel training Document control and revisions		
A detailed account of how the QA Program will be implemented	Facilities and Equipment Procurement and inventory procedures Preventive maintenance		
A description of the organization, responsibilities, and decision-making authorities of the contractor project team	Analytical Methodology Calibration and operating procedures		
A description of sampling team and analyst training in technical skills, standard QC, and essential elements of QA Program	Sample Custody		
Procedures for sampling, preservation, and shipment of samples	Quality Control Quality control procedures Control checks and internal audits Reference material analysis Blank analysis Matrix spike and matrix spike duplicate analysis Internal audits		
Sample inspection and lot sizing			
Instrument calibration			
Logs (field, instrument, sample, QC)			
Analytical reference materials	Data Handling Data handling, reporting, and record-keeping procedures Data validation		
Procedures for verifying and documenting the quality of lab water			
Control charts			
Methods and criteria for determining when sampling or analytical systems are out of control, including holding times			
Actions to be taken to correct out-of-control situations, and how actions will be reported and documented			
A list of personnel responsible for data review and sequence of review prior to submittal			

QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability (cont.)	USATHAMA	EPA CLP	COMMENT
Generation and submission of pre-certification performance data package		Preparation for performance evaluation samples	
Precertification method description (preparation and analysis of standards)			EPA CLP requires the analysis of blind evaluation samples for evaluation of laboratory performance
Standardized method written to be laboratory-specific			USATHAMA requires analysis of certification samples which are not blind
Development of method		Tuning and GC/MS mass calibration	Assumption is made that GC/MS tuning is described as part of the method
Submit documentation for proposed method			
Analytical procedures testing			
Documentation of method in standard format			
Generation of performance data packages			
Review by USATHAMA Analytical Branch			
Assignment of method number after final approval			
Precertification calibration data			
Construction of calibration curve		Construction of Calibration Curves	
Prepare and analyze each standard in duplicate to bracket desired range for certification			
TRL - Target Reporting Limit (designated by USATHAMA)			
Class 1		Pesticides	
Blank, 0, 5, 1, 2, 5, and 10 times the TRL plus expanded range		Established retention time windows	
		Run evaluation standard mix at five concentrations	
		Run standard mix of pesticides	
		Run individual Aroclors	
Class 1A and Class 1B			

5. QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability (cont.)	USATHAMA	EPA-CLP	COMMENTS
Blank, 0.5, 2, and 10 times the TRL + range extension (10% for inorganic and 25% for all others)  Tabulate and graph response vs. concentration Lack of fit (LOF) Zero intercept (ZI)	GCMS: semivolatiles and volatiles require 5 point calibration curve with specified concentrations of 20, 50, 80, 120, and 160 ngs  Tabulation of calibration data GCMS: Relative response factors Relative standard deviation Calibration factors		
Certified calibration check standard - Class 1 and 1B only	GC: % RSD % Breakdown  Verification of performance checks GCMS: System performance check Calibration check		
Class 1	Two calibration check standards should be analyzed, one at the beginning and one at the end of the day - near high end of range	GC: Retention time shifts % Breakdown	
Class 1B	One calibration check standard should be analyzed at the beginning of the day. New high end of range		
Results of identification and purity analyses for all off-the-shelf reference materials	Checklist completed by the QAC		
Approval by USATHAMA of Pre-certification Performance Data Package and Project QA Plan	Generation and submission of Certification Performance Data Package Final USATHAMA-approved copy of the Precertification Performance		

5. QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability (cont.)

USATHAMA	EPA-CLP	COMMENTS
Data Package		
Total method description in USATHAMA format containing approved deviations in the standardized method and laboratory-specific information concerning conduct of the method	Submission of Standard Operating Procedures	
Initial calibration MTR = minimum testing range TRL = target reporting limit * = times		Calibration procedure for semivolatiles and volatiles for EPA-CLP resembles USATHAMA Class 1 more than Class 1A which is reserved for all GC/MS methods.
Class 1  MTR: blank, 0.5, 1, 2, 5, *10, and *10TRL 7 standards + 2 check standards MTR + 1 range extension; 10 standards + 2 check standards (20, 50, 100, 100) MTR + 2 range extensions; 13 standards + 2 check standards (20, 50, 100, 200, 500, 1000, 1000)		However, USATHAMA calibration for pesticides (assuming Class 1) is more stringent than EPA-CLP.
Class 1A  MTR: blank, 0.5, 2, 10 and 10 TRL; 5 standards MTR + 1 range extension; (50, 200, 200); 7 standards MTR + 2 range extension; (50, 200, 500, 2000, 2000); 9 standards Class 1B - same as 1A plus 1 check standard	GCMS: semivolatiles and volatiles require a 5 point initial calibration at specified concentrations	
Class 2 - 6 standards, blank, and 1 triplicate TRL		
Daily calibration  Class 1/Class 1A/Class 1B 2 standards for MTR - *10 and *10 TRL 2 standards for MTR + 1 range	GCMS: 50 total ngs. standard analyzed each 12 hours	Daily calibration for USATHAMA requires analysis of a high standard twice whereas EPA-CLP requires analysis of a lower range standard.



5. QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability (cont.)

USATHAMA	EPA-CLP	COMMENTS
extension - *100 and *100 TRL 2 standards for MTR and 2 range extensions - *1000 and *1000 TRL	Response must be within 25% for organics of mean response of 5 initial calibration standards	With USATHAMA certification samples, the participating laboratory knows immediately whether problems exist in sample preparation and/or analysis. However, this same knowledge is available to EPA-CLP laboratories only if the results of the evaluation samples are returned promptly.
Class 2		
4 standards (MTR) and blank and 1 TRL (duplicate)		
Certification samples (prepared in standard matrix)	Performance Evaluation	
Class 1/Class 18 MTR: 24 Blank, 0.5, 1, 2, 5, and 10 TRL (4 consecutive days)	Samples prepared by EMSL/LV are sent to laboratory	
MTR + 1 range extension: 36 Blank, 0.5, 1, 2, 5, 10, 20, 50, 100 TRL (4 days)		
MTR + 2 range extensions: 48 Blank, 0.5, 1, 2, 5, 10, 20, 50, 100, 200, 500, and 1000 TRL (4 days)		
Class 1A MTR: 8 Blank, 0.5, 2, and 10 TRL (duplicate)		
MTR + 1 range extensions: 12 Blank, 0.5, 2, 10, 50, and 200 TRL (duplicate)		
MTR + 2 range extensions: 16 Blank, 0.5, 2, 10, 50, 200, 500, and 2000 TRL (duplicate)		
Class 2 MTR: 8 Blank, 1 TRL (quadruplicate)		
Statistical Analysis of the Data	Data Package	
Tabulation of found vs. target concentration	Sample Traffic Report Sample Data Summary Package	
LOF and ZI test calculations and	Case narrative	

5. QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability (cont.)		
USATHAMA	EPA-CLP	COMMENTS
<p>results for the pooled data set for found vs. target concentration</p> <p>Linear regression Confidence bounds Reporting limit Accuracy Standard deviation Percent Imprecision (% RSD) Percent Inaccuracy</p> <p>Narrative evaluation of effectiveness of method</p> <p>Checklist completed by QAC</p>	<p>Target compound results-(Form I) Tentatively identified compounds-(Form I)</p> <p>Surrogate spike analysis results-(Form II) Matrix spike/matrix spike duplicate-(Form III) Blank data-(Form IV and Form I)</p>	<p>EPA-CLP surrogate spike is a measure of percent inaccuracy and matrix spike/matrix spike duplicate is a measure of percent imprecision.</p>
	<p>Sample Data Package Case narrative Traffic reports Volatiles data QC summary</p> <p>Surrogate spike results, Form II Matrix spike results, Form III Method blank summary, Form IV GC/MS tuning standard, Form V</p> <p>Sample data TCL results, Form I Ion chromatograms Mass Spectra Library search spectra for TIC Quantitation of TIC Manual work sheets Standards data Initial calibration data, Form VI Continuing calibration data, Form VII Internal standards summary, Form VIII</p> <p>Raw QC data BFB mass spectra Blank data, Form I Ion chromatograms Mass spectra Library search spectra for TIC Quantitation of TIC Manual work sheets Matrix spike results, Form I Matrix spike duplicate results, Form I</p> <p>Semivolatiles data QC summary</p>	

5. QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability (cont.)

USATRAMA

Calibration data, tabulation of concentration vs. response  
(a) Initial  
(b) Daily

Calibration curves (instrument response on ordinate vs. concentration on abscissa)

Data demonstrating that the response for Daily Calibration standards was within required percentage of highest standard

Copies of chromatograms for high and low standards - certification samples

analysis date and time

target concentration

test name

reference to calibration curve

reference to analytical logbook

each peak labeled

Spectra for all target analytes

EPA-CLP

Surrogate spike results, Form II

Matrix spike results, Form III

Method blank summary, Form IV

GC/MS tuning standard, Form V

Sample data

TCL results, Form I

Ion chromatograms

Mass spectra

Library search spectra for TIC

Quantitation of TIC

Manual work sheets

Standards data

Initial calibration data, Form VI

Continuing calibration data, Form VII

Internal standards summary, VIII

Raw QC data

SFB mass spectra

Blank data, Form I

Ion chromatogram

Mass spectra

Library search spectra for TIC

Quantitation of TIC

Manual work sheets

Matrix spike results, Form I

Matrix spike duplicate results, Form I

Pesticides/PCB data

QC summary

Surrogate spike results, Form II

Matrix spike results, Form III

Method blank summary, Form IV

Sample data

TCL results, Form I

Gas chromatograms

Confirmation gas chromatograms

Manual work sheets

GPC chromatograms

GC/MS raw spectra

Standards data

Evaluation standards summary, Form VIII

COMMENTS

5. QA Objectives for Measurement Data in Terms of Precision, Accuracy, Completeness, Representativeness, and Comparability (cont.)

USATHAMA

EPA-CLP

COMMENTS

Standards summary, Form IX  
 Identification, Form X  
 Chromatograms  
 Raw QC data  
 Blank data, Form I  
 Gas chromatograms  
 Matrix spike results, Form I  
 Gas chromatograms and printouts  
 Matrix spike duplicate results, Form I  
 Gas chromatograms and printouts  
 Data evaluated for accuracy by NFO (National Program Office) and audited by EMSL/LV personnel  
 Quality control data goes into EMSL/LV database for trend analyses, etc.  
 On-site laboratory evaluation

# 6. Sampling Procedure

## (A) Sample containers (Appendix C and Appendix D)

Inorganics		COMMENTS
	USATHAMA	
Water	Polyethylene (Exception: glass bottle and top for dissolved oxygen)	EPA-CLP Polyethylene (Medium level requires wide-mouth glass jar)
	Amber glass bottle with Teflon- lined cap	Wide-mouth glass jar
Soil		
Organics		
Water Volatiles	Glass vial with Teflon-lined septum cap	Glass vial
	Amber glass bottle with Teflon-lined cap	Amber glass bottle (wide-mouth glass jar for medium level)
Soil Volatiles	Glass vial with Teflon-lined septum cap	Wide-mouth glass vial
	Amber glass bottle with Teflon-lined cap	Wide-mouth glass jar

6. Sampling Procedure (cont.)

(B) Sample container cleaning procedures (Appendix E and Appendix F)

USATHAMA		EPA-CLP	COMMENTS*
Polyethylene bottles and caps	5% sodium hydroxide deionized water 5% Ultrax nitric acid/water deionized water air dry	Cleaning procedure used by EPA-CLP Sample Bottle Repository not known at this time	These procedures are referenced by companies advertising precleaned bottles
Amber glass bottles or vials	detergent distilled water acetone methylene chloride hexane air dry heat to 200°C		Polyethylene bottles detergent tap water 1:1 nitric acid tap water 1:1 hydrochloric acid tap water distilled water VOA vials detergent tap water distilled water dry at 105°C
Bottle caps	remove paper liners detergent distilled water dry at 40°C		EXTRACTABLE bottles detergent tap water distilled water acetone hexane (pesticide quality) air dry (muffle furnace heating may be substituted for solvent rinses)
Teflon liners	detergent distilled water acetone hexane air dry place liners in cleaned caps heat to 40°C for 2 hours cool		detergent tap water distilled water dry at 105°C for 1 hour

\*EPA 40 CFR 136 "Guidelines for Establishing Test Procedures for the Analysis of Pollutants"

6 Sampling Procedure (cont.)

(C) Sample holding times (Appendix C and Appendix D)

Definitions

USATHAMA - maximum time allowable between sample collection and analysis

EPA-CLP - maximum time allowable between verified time of receipt (VTSR) and analysis

	USATHAMA	EPA-CLP	COMMENTS
Inorganics			
Metals	6 months (chromium VI-24 hours)	6 months	Major difference in definition of holding time probably reflects also the differences noted in the holding times.
Mercury	28 days	26 days	Holding times are selected arbitrarily and by convenience.
Cyanide	14 days	14 days	Neither plan is necessarily correct.
Organics			
Extractables			
Soil	7 days	10 days	
Water	7 days	5 days	
Volatiles			
Soil	14 days	10 days	
Water	14 days (with pH adjustment) 7 days (no pH adjustment)	10 days	

# 6 Sampling Procedure (cont.)

## (D) Sample preservation and storage conditions

### Definition of sample storage termination:

USATHAMA - sample storage shall only be terminated after all analytical results have been validated to Level 3 in the USATHAMA Data Management System

EPA-CLP - sample extracts shall be retained for 365 days after data submission

	USATHAMA	EPA-CLP	COMMENTS
Inorganics			
Metals*			The length of the sample archival depends on the program.
Soil	Cool, 4°C		Sample preservation and storage conditions are basically the same except for volatiles. Preservation in USATHAMA plan is consistent with EPA 40 CFR 136.
Water	HNO <sub>3</sub> to pH<2 (except chromium VI)	HNO <sub>3</sub> to pH<2	
Mercury			
Soil	Cool, 4°C		
Water	HNO <sub>3</sub> to pH<2	HNO <sub>3</sub> to pH<2	
Cyanide			
Soil	Cool, 4°C		
Water	Cool, 4°C NaOH to pH>12 0.6g ascorbic acid	Cool, 4°C NaOH to pH>12 0.6g ascorbic acid	

Organics		
Extractables	Store in dark Cool, 4°C	Store in dark Cool, 4°C
Volatiles	Cool, 4°C 0.008% Na <sub>2</sub> S <sub>2</sub> O <sub>3</sub> if chlorine is present HCl to pH<2 (or analyze within 7 days)	Store in dark Cool, 4°C

\*Samples should be filtered before adding preservative for dissolved metals.



7. Sample Custody

USATHAMA

EPA-CLP

COMMENTS

Field sampling

Sample acquisition as well as distinguishing information recorded in bound logbook

Sample Traffic Report

Sample label

Installation name  
Unique sequential field sample no.  
Sampling date  
Analytes of interest  
Preservative/filtration

Sample tag-information defined by EPA National Enforcement Investigations Center (NEIC)  
CLP case/SAS no(s). CLP sample no.  
Project code Station no.  
Date Time  
Station location Samplers  
Remarks Tag no.  
Lab sample no.

Use of formal Chain-of-Custody procedures implied for litigation

Chain-of-Custody Record

USATHAMA needs to document chain-of-custody procedures.

Laboratory Operation

Sample login

Samples are logged into a project-specific logbook

Samples are grouped into analytical lots, ordered and assigned a USATHAMA sample identification number (QC samples also)

USATHAMA's procedures for sample login and analysis are covered in the Project QC Plan under Laboratory Certification.

Sample analysis

Bound logbooks required for: reference materials; operational activities which occur during sample handling; instrument operation

Standard operating procedures required for: receipt of samples; maintenance of custody; sample storage

Standard operating procedure for tracking the analyses of samples required  
Bound logbooks required for entering all observations and results not on pre-printed data sheets

Document control

All documentation shall be in ink

Errors shall be corrected by crossing a line through the error, entering the correct information, and dating and initialing the change

Standard operating procedure for the assembly of completed data

All documentation shall be in ink

USATHAMA requires a list of personnel responsible for data review and sequence of review in Project QC Plan.

7. Sample Custody (cont.)	USATHAMA	EPA CLP	COMMENTS
Computerized logging systems may not be used for original records	Logbook should be installation-specific	Errors shall be corrected by crossing a line through the error, entering the correct information, and dating and initialing the change	Documentation is cross-checked for consistency Documents are numbered and inventoried

# 8 Calibration Procedures and Frequency

## USATHAMA

## EPA CLP

## COMMENTS

### Initial calibration

#### Frequency

- (a) 1st day of certification analyses
- (b) Instrumental start-up (not daily)
- (c) Analyzing different analytes
- (d) Daily calibration fails

If samples are analyzed on the same day as initial calibration, one standard at the highest concentration must be analyzed after analyses are completed

#### Concentration of standards

Class 1  
MTR: blank, 0.5, 1, 2, 5, \*10, and \*10 TRL, 7 standards + 2 check standards  
MTR + 1 range extension: 10 standards + 2 check standards (20, 50, 100, 100)  
MTR + 2 range extensions: 13 standards + 2 check standards (20, 50, 100, 200, 500, 1000, 1000)  
Class 1A  
MTR: blank, 0.5, 2, 10, & 10 TRL; 5 standards  
MTR + 1 range ext.: (50, 200, 200); 7 standards  
MTR + 2 range ext.: (50, 200, 500, 2000, 2000); 9 standards  
Class 1B - same as 1A plus 1 check standard  
Class 2 - 6 standards, blank, and 1 triplicate TRL

#### Certified check standards

Class 1 - two stds - beginning & end of day  
Class 1B - one std - beginning of day - near high end of range

If acceptability limits are exceeded, immediate reanalysis occurs, followed by a new initial calibration if necessary

#### Frequency

Prior to analysis of samples and if daily calibration fails

Frequency requirements are equivalent

#### Concentration of standards

GC/MS (= Class 1A)  
Volatiles  
20, 50, 100, 150, and 200 µg/L  
The 1 RSD for each calibration check compound must be less than or equal to 30.0%

The minimum acceptable average relative response factor is 0.300, 0.250 for bromoform

#### Semivolatiles

20, 50, 80, 120, and 160 total nanograms

The 1 RSD for each calibration check compound must be less than or equal to 30.0%

The minimum acceptable average relative response factor is 0.05

#### GC (= Class 1)

#### Pesticides

#### Evaluation standard

Mixture of aldrin, endrin, and 4,4'-DDT at concentrations of 20%, 50%, and 100% full-scale  
Individual standard mixes and areolers

Calibration procedure for semivolatiles and volatiles for EPA-CLP resembles USATHAMA Class 1 more than Class 1A which is reserved for all GC/MS methods  
However, USATHAMA calibration for pesticides (assuming Class 1) is more stringent than EPA-CLP.

Certified USATHAMA check standards made from the stock solution used during certification allow a continual check of laboratory performance.

8. Calibration Procedures and Frequency (cont.)

USATHAMA

EPA-CLP

COMMENTS

Daily calibration

Class 1, 1A, 1B zero-intercept  
Highest concentration standard is  
analyzed at beginning and end of day  
Response must be within 10% for  
inorganic and 2% for others of the  
mean response for the same concentra-  
tion as determined for precertifica-  
tion and certification for the 1st  
7 calibrations

After 7 calibrations, response must  
agree within 2 standard deviations

Corrective action

Reanalyze daily standard

Initial calibration repeated

Non-linear or non-zero intercept  
Analyze low, middle, and high  
calibration standards at beginning  
of day and low and high standards at  
the end of the day. (If quadratic,  
four standards)

Responses must fall within 2 std.  
deviations of the mean response

Class 2

One blank and one calibration  
standard at the CRC analyzed at  
beginning and end of sample analysis

The  $\pm$  RSD for evaluation standard  
mix compounds must be  $\leq 10.0\%$ . The  
 $\pm$  breakdown for endrin or 4,4'-DDT  
must not exceed 20.0%. The cali-  
bration factor for each individual  
standard must not exceed a 15.0%  
difference for a quantitation run  
nor exceed a 20.0% difference for  
a confirmation run.

GC/MS ( $\pm$  Class 1A)

Volatiles

50  $\mu\text{g/L}$  standard is analyzed every  
12 hours

The  $\pm$  difference for each calibration  
check compound must be less than or equal  
to 25.0%. The minimum relative response  
factor for the system performance check  
compounds is 0.300 (0.250 for bromoform)

Daily calibration for USATHAMA  
requires analysis of the high  
standard twice whereas EPA CLP  
requires analysis of the lower  
range standard.  
The quality of data should be  
equivalent.

Semivolatile

50 total ngs standard is analyzed  
every 12 hours

The minimum relative response factor  
for the system performance check  
compounds is 0.050

The  $\pm$  difference for each calibration  
check compound must be less than or  
equal to 25.0%

Pesticides

Analyze evaluation standard Mix B and  
individual standard Mix A or B alter-  
nately after every 5 samples

The  $\pm$  difference in retention time for  
the dibutylchloride must not exceed

8. Calibration Procedures and Frequency (cont.)

USATHAMA	EPA-CLP	COMMENTS
<p><u>Reference Materials</u></p> <p>Standard Analytical Reference materials traceable to NBS</p> <p>Interim reference material</p> <p>(a) Central QA Lab</p> <p>(b) EPA</p> <p>(c) NBS</p> <p>Off-the-shelf material</p> <p>(a) Positive identification</p> <p>(b) Estimate of purity</p>	<p>0.3% for capillary or 2.0% for packed column</p> <p>The % breakdown for 4,4'-DDT or endrin must not exceed 20.0%</p>	<p>USATHAMA provides reference materials to prepare all standard solutions.</p> <p>EPA also makes available QC samples intended for periodic (quarterly) use as independent checks on each laboratory's own QC activities.</p> <p>No practical difference</p>
<p>EMSL/LV provides standard materials from its QA Materials Bank for performing initial instrument calibration and as reference standards</p>		

9. Analytical Procedures

USATHAMA uses EPA standardized methods for commonly encountered analytes and USATHAMA-specific methods are used when no EPA comparison is available.

10 Data Reduction, Validation, and Reporting

USATHAMA

EPA-CLP

COMMENTS

Data reporting

CRL = certified reporting limit

Data is not adjusted by any correction factors (such as accuracy, % moisture, and dilution factor), but is reported in the as-received condition

Class 1, 1A, 1B  
All values less than CRL will be reported as <RL

Number of Significant Figures to be used in Reporting Data

Class 1 and 1B

No dilution - 3 significant figures  
Dilution - 2 significant figures  
Noncertified analytes - retention time

Class 1A

No dilution - 2 significant figures  
After dilution - 1 significant figure  
Screening for noncertified - 1 significant figure

Class 2

CRL - 2 significant figures  
Reported as >, <, = CRL

Deliverables

Specific instructions for format, coding, and submission are provided in the IRDMS User's Guide

Soil/sediment data is adjusted to Dry Weight Basis

Values less than quantitation limit are reported with J qualifier

GC/MS:

Report data to 2 significant figures

GC-Pesticides:

Report data to 2 significant figures

Deliverables

Inorganic

- (1) Weekly process reports
- (2) Sample traffic report
- (3) Sample data package
- Tabulated results
- Raw data
- Copies of logbook entries

Organic

- (1) Narrative report
- (2) Sample traffic report
- (3) Quality control summary
- (4) Sample data
- (5) Raw sample data
- (6) Standards package
- (7) QC data package

Note difference in reporting of soil/sediment data. Either report is acceptable if the end user is aware of the difference.

# 11. Internal Quality Control Checks

Types	USATHAMA	EPA-CLP	COMMENTS
<p>Class 1 and Class 1B Method blank Spikes of control analytes in standard matrices</p>		<p><u>Inorganics</u> Preparation blank analysis Interference check sample analysis ICP serial dilution analysis Matrix spike analysis Duplicate sample analysis Furnace AA QC Analysis (Method of Standard Addition may be required under certain conditions) Laboratory quality control sample analysis</p>	<p>USATHAMA does not require matrix spiking (as EPA perceives) for organics. Matrix spikes could easily be added to USATHAMA plan. Frequency should be as in CLP. Matrix spikes are probably not necessary if surrogates are added to each sample, unless surrogate recovery is low.</p>
<p>Class 1A (GC/MS) Method blank/surrogate spikes Surrogates spikes in every field sample</p>			
<p>Class 2 Method blank Spikes of control analytes in standard matrices</p>		<p><u>Organics</u> Method blank analysis Surrogate spike analysis Matrix spike/Matrix spike duplicate analysis</p>	

## Frequency per lot

<p>Class 1 One - standard matrix method blank Three standard matrix spikes 2, 10, &amp; 10 CRL</p>	<p><u>Inorganics</u> Preparation blank - every 20 samples received or with each batch of samples digested whichever is more frequent</p>
<p>Class 1A One - standard matrix method blank spike (surrogate, 10 CRL) All field samples spiked with surrogate - 10 CRL</p>	<p>Interference check sample - analyzed at beginning and end of each analysis run or a minimum of twice per 8 hour working shift</p>
<p>Class 1B One - standard matrix method blank One - standard matrix spikes - 10 CRL</p>	<p>ICP serial dilution - each group of samples of a similar matrix type and concentration for each case of samples or for each 20 samples received, whichever is more frequent</p>
<p>Class 2 (a) One - standard matrix method blank (b) One standard matrix spike - 1 CRL</p>	<p>Spiked sample and duplicate sample - at least one for each group of samples of a similar matrix and concentration for each case of samples or for each 20 samples received, whichever is more frequent</p>



11. Internal Quality Control Checks (cont.)

USATHAMA

EPA-CLP

COMMENTS

Laboratory control sample - one for each group of 20 samples of a similar matrix or for each batch of samples digested whichever is more frequent

Organics

Method blank analysis

Volatiles

For the analysis of volatile TCL compounds, a method blank analysis must be performed once for each 12-hour time period during the analysis of samples from:

- o each case, OR
- o each 14 calendar day period during which samples in a case are received (said period beginning with the receipt of the first sample in that sample delivery group), OR
- o each 20 samples in a case that are of similar matrix (water or soil) or similar concentration (soil only),

whichever is most frequent, on each GC/MS system used to analyze samples

Extractables

For the analysis of extractable TCL compounds, a method blank analysis must be performed once:

- o each case, OR
- o each 14 calendar day period during which samples in a case are received (said period beginning with the receipt of the first sample in that sample delivery group), OR
- o each 20 samples in a case that are of similar matrix (water or soil) or similar concentration (soil only), OR

Method blank requirements are equivalent.

11. Internal Quality Control Checks (cont.)

USATHAMA

KPA (LP)

COMMENTS

- o whenever samples are extracted by the same procedure (separatory funnel or continuous extraction),

whichever is most frequent, on each GC/MS or GC system used to analyze samples

Surrogate spike analysis

All blanks, field samples, matrix spikes, and matrix spike duplicates will be spiked with surrogate compounds

Matrix spike analysis

A matrix spike and matrix spike duplicate must be performed for each group of samples of a similar matrix, once:

- o each case of field samples received, OR
- o each 20 field samples in a case, OR
- o each group of samples of a similar concentration level (soils only), OR
- o each 14 calendar day period during which samples in a case were received (said period beginning with the receipt of the first sample in that sample delivery group).

whichever is most frequent.

See earlier comments on matrix spiking.

11 Internal Quality Control Checks (cont.)

USATHAMA

EPA-CLP

COMMENTS

Preparation

Assigned sample number during  
logging-in process  
Spiked samples (excluding water  
samples) must be allowed to stand  
for one hour before continuing the  
analysis

Data Reporting

Class 1  
Minimum of 3 significant figures  
Method blank: can be corrected -  
reported by concentration  
Control charts

Class 1A  
2 significant figures  
Method blank: can be corrected -  
reported by concentration  
Control charts

Class 1B  
Minimum of 3 significant figures  
Method blank: can be corrected -  
Reported by concentration  
Control charts

Class 2  
Minimum of 2 significant figures  
No control charts

Data Reporting

Reporting of quality control samples  
handled just as samples are

Soil/sediment results are corrected  
for percent moisture and reported  
on a dry weight basis

No corrections are made for method  
blanks

Method blanks can be corrected in  
USATHAMA plan, but cannot be  
corrected in CLP. Blank correction  
is fine, but any time this is done  
the value should be documented.  
CLP should require control charts  
at least for surrogates and  
internal standards. In earlier  
IFB's, control charts were required  
for internal standards.

## 12 Performance and System Audits

### Definition

Performance audit - Evaluation to determine the accuracy of the total measurement system or components thereof

System audit - Evaluation to determine the proper selection and use of the measurement system, or components thereof

### External audits

#### USATHAMA

Reviewer:  
USATHAMA Analytical Branch

#### EPA-CLP

NPO Project Office  
Regional personnel  
EHS/LV personnel  
CLP SMO  
NEIC

#### COMMENTS

No substantial differences  
If Performance Evaluation samples  
were a required part of the  
USATHAMA plan, audits could be  
done simultaneously.

### Frequency:

After review of the project  
QC plan  
After initiation of analyses  
Other visits as deemed necessary

After first performance evaluation  
samples are completed  
Yearly  
Repeat site visit as needed

### Documentation:

Checklist for laboratory adherence

Performance Evaluation sample score  
sheets  
Trend analysis  
Laboratory evaluation checklists

### Circulation of Audit Report:

USATHAMA Project Officer  
Contractor Project Manager  
Analytical Task Manager  
Contractor QAC  
USATHAMA Analytical Branch

### Corrective action:

Serious deficiencies are reported to  
the Contracting Officer at Procure-  
ment for action

Specified in contract under each QC  
section

Evaluated by PO who may initiate a  
site visit, full data audit, or  
analysis of a second PE sample  
Laboratory may be placed on temporary  
hold  
Formulate recovery plan and SHOW CAUSE NOTICE

12 Performance and System Audits (cont )

USATHAMA

EPA-CLP

COMMENTS

Internal audits

Reviewer

Project QC staff

Frequency

Not specified

Should be periodically conducted to evaluate the functioning of the QA SOP and involves an independent check of the laboratory analysts to ensure that procedures are being followed

Documentation:

Verification of maintenance of standards procedures, records, etc.  
Verification of actual practice vs written procedures  
Verification of QA records and results of QC sample analyses

Audit findings must be in a bound logbook

Circulation of Audit Report:

Project Manager  
Analytical Task Leader  
USATHAMA

13 Preventive Maintenance

USATHAMA	EPA-CLP	COMMENTS
<p><b>Schedule:</b>            Must maintain a calibration and maintenance schedule for each instrument as recommended by the manufacturer            Physical or electronic measurements or calibrations must be traceable to NBS</p> <p><b>Supplies:</b>            An adequate supply of critical spare parts must be maintained</p> <p><b>Documentation:</b>            Reports and records must be available for inspection</p>	<p>Not specifically stated in contract, however, the following items are included in the laboratory audit checklist:</p> <ul style="list-style-type: none"> <li>service maintenance</li> <li>in-house replacement parts</li> <li>preventative maintenance</li> <li>permanent service record logbook</li> <li>instrument modifications</li> </ul>	No difference

14	Specific Routine Procedures Used to Assess Data Precision, Accuracy, and Completeness				
		USATHAMA	EPA-CLP	COMMENTS	
	Software provided to assess precision and accuracy during certification	Contract specifies equations to evaluate precision and accuracy of matrix and surrogate spikes			
	USATHAMA maintains a data management system which automates the statistical analyses of the data	Data is manually entered or copied from a floppy diskette into the EMSL/LV database for more extensive statistical review			

USATHAMA procedure is superior and should be implemented if possible in CLP.

15. Corrective Action

USATHAMA

EPA-CLP

COMMENTS

Personnel responsible for initiating action:

Analyst  
QAC  
Analytical Task Manager  
Project Manager

Analyst  
Project Officer

No difference

Action:

Immediate  
Repairing equipment  
Making a new standard  
Long term  
Staff training  
Rescheduling  
Replacing vendors  
Revising of QA system  
Personnel replacement

Specified in contract under each QC section

Evaluated by PO  
Laboratory may be placed on temporary hold  
Formulate recovery plan

Documentation:

Required

SHOW CAUSE NOTICE




16. Quality Assurance Reports to Management

USATHAMA	EPA-CLP	COMMENTS
Precertification and certification performance data package IRDMS submissions Audit reports Control charts - provided to USATHAMA Project Officer weekly	Performance evaluation data package Data package submission Audit reports Quarterly Blind Performance Evaluation samples	QA report differences reflect the differences in other aspects of the plans, such as control charts, and Performance Evaluation samples. The advantages of each could be used.
Final Project QC Data Report		

Appendix A

LINES OF COMMUNICATION FOR USATHAMA IR PROJECTS  
(USATHAMA QA PROGRAM, 2ND EDITION, MARCH, 1987)


Figure 2-1. Lines of Communication for USATHAMA IR Projects

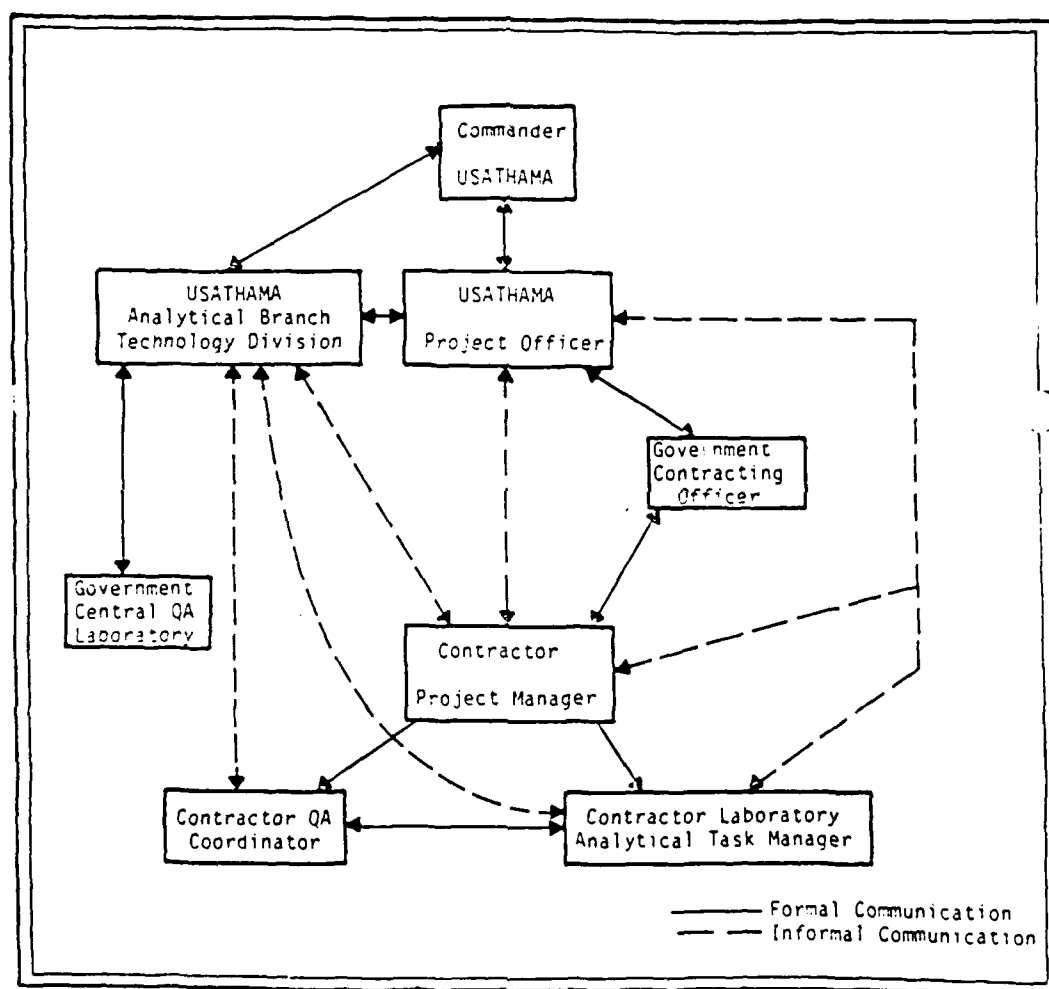

Appendix B

INTERRELATIONSHIP OF PROGRAM PRINCIPALS  
(USER'S GUIDE TO THE CONTRACT LABORATORY PROGRAM, OCTOBER, 1984)

## INTERRELATIONSHIP OF PROGRAM PRINCIPALS

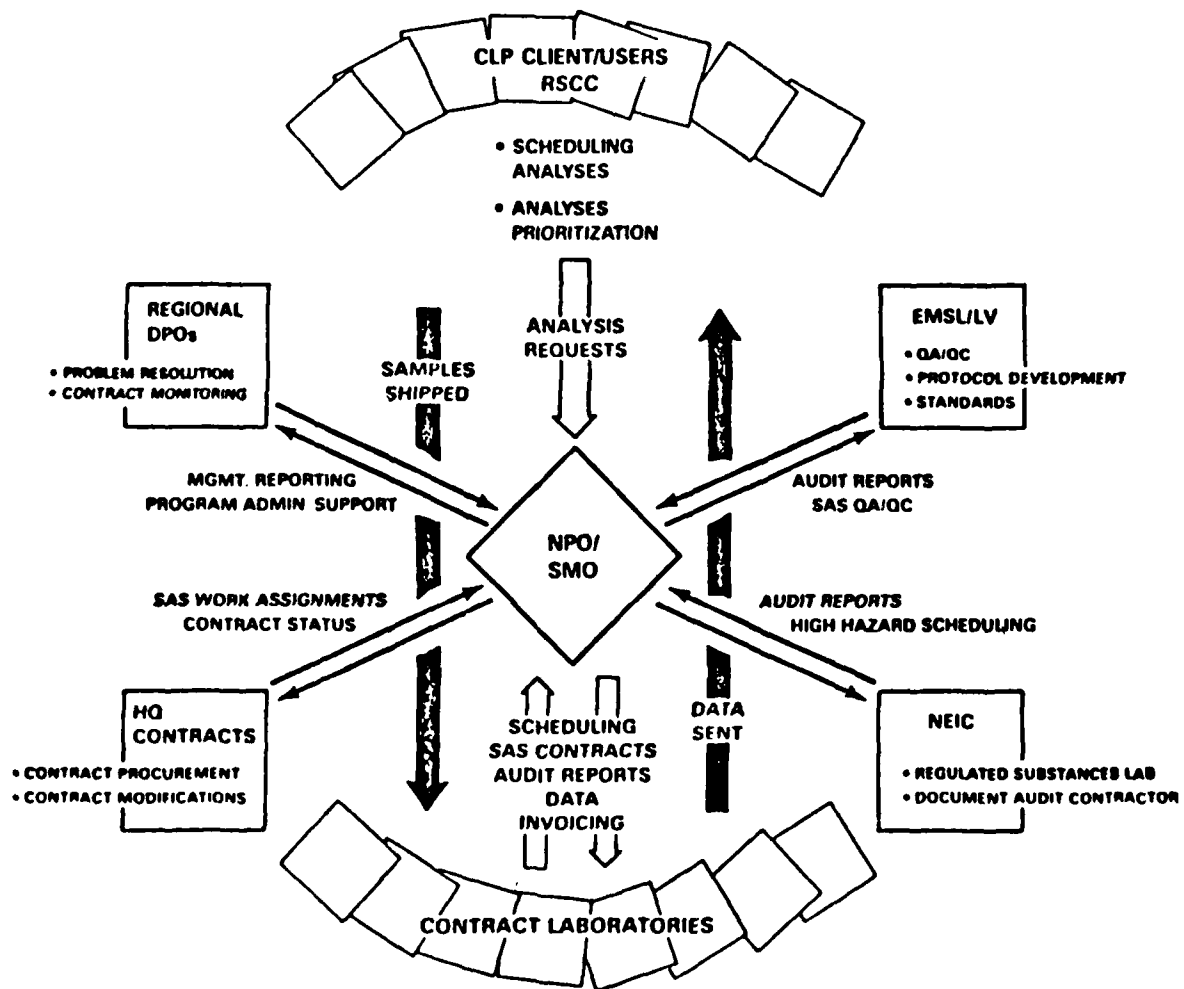

Appendix C

CONTAINERS, PRESERVATION, STORAGE, AND HOLDING TIMES  
AND  
SAMPLE CONTAINER CLEANING PROCEDURES

(USATHAMA QA PROGRAM, 2ND EDITION, MARCH, 1987)

Table H-1. Containers, Preservation, Storage, and Holding Times<sup>a</sup>

Parameter	Container <sup>b</sup>		Preservative <sup>c,d</sup>		Maximum Holding Time for all Matrices <sup>e</sup>
	Water	Soil	Water	Soil	
INORGANIC TESTS					
Acidity	P	G	Cool, 4°C	Cool, 4°C	14 days
Alkalinity	P	G	Cool, 4°C	Cool, 4°C	14 days
Ammonia	P	G	Cool, 4°C H <sub>2</sub> SO <sub>4</sub> to pH < 2	Cool, 4°C	28 days
Asbestos	P	G	Cool, 4°C	Cool, 4°C	48 hours <sup>f</sup>
Bicarbonate	P	G	None Required	None Required	Analyze Immediately
Biochemical Oxygen Demand (BOD) and Carbonaceous BOD	P	G	Cool, 4°C	Cool, 4°C	48 hours
	P	G	None Required	None Required	28 days
Carbonate	P	G	None Required	None Required	Analyze Immediately
Chemical Oxygen Demand (COD)	P	G	Cool, 4°C H <sub>2</sub> SO <sub>4</sub> to pH < 2	Cool, 4°C	28 days
Chloride	P	G	None Required	None Required	28 days
Chlorine, Total Residual	P	N/A	None Required	N/A	Analyze Immediately
Color	P	N/A	Cool, 4°C	N/A	48 hours

Section No. H  
Revision No. 0  
Date December 1985  
Page 3 of 10

Section No. H  
Revision No. 0  
Date December 1985  
Page 4 of 10

Table H-1. (Cont'd.)

Parameter	Container <sup>b</sup>		Preservative <sup>c,d</sup>		Maximum Holding Time for all Matrices <sup>e</sup>
	Water	Soil	Water	Soil	
Cyanide, Total and Amenable to Chlorination	P	G	Cool, 4°C NaOH to pH >12 0.6 g Ascorbic Acid <sup>g</sup>	Cool, 4°C	14 days <sup>h</sup>
Dissolved Oxygen Probe	G Bottle and Top	N/A	None Required	N/A	Analyze Immediately
Winkler	G Bottle and Top	N/A	Fix On Site Store in Dark	N/A	8 hours
Fluoride	P	G	None Required	None Required	28 days
Hardness	P	N/A	HNO <sub>3</sub> or H <sub>2</sub> SO <sub>4</sub> to pH <2	N/A	6 months
Hydrazine	P	G	If not analyzed immediately, collect under acid. Add 90 ml of sample to 10 ml HCl.	Cool, 4°C	7 days
Iodide	P	G	Cool, 4°C	Cool, 4°C	24 hours
Iodine	P	G	None Required	None Required	Analyze Immediately
Kjeldahl and Organic Nitrogen	P	G	Cool, 4°C H <sub>2</sub> SO <sub>4</sub> to pH <2	Cool, 4°C	28 days



Section No. H  
Revision No. 0  
Date December 1985  
Page 5 or 10

Table H-1. (Cont'd.)

Parameter	Container <sup>b</sup>		Preservative <sup>c,d</sup>		Maximum Holding Time for all Matrices <sup>e</sup>
	Water	Soil	Water	Soil	
<b>Metals<sup>1</sup></b>					
Chromium VI	P	G	Cool, 4°C	*Cool, 4°C	24 hours
Mercury	P	G	HNO <sub>3</sub> to pH < 2	Cool, 4°C	28 days
Others	P	G	HNO <sub>3</sub> to pH < 2	Cool, 4°C	6 months
Nitrate	P	G	Cool, 4°C	Cool, 4°C	48 hours
Nitrate plus Nitrite	P	G	Cool, 4°C H <sub>2</sub> SO <sub>4</sub> to pH < 2	Cool, 4°C	28 days
Nitrite	P	G	Cool, 4°C	Cool, 4°C	48 hours
Oil and Grease	G	G	Cool, 4°C H <sub>2</sub> SO <sub>4</sub> to pH < 2	Cool, 4°C	28 days
Orthophosphate	P	G	Filter Immediately Cool, 4°C	Cool, 4°C	48 hours
pH	P	G	None Required	None Required	Analyze Immediately
Phenols	G	G	Cool, 4°C H <sub>2</sub> SO <sub>4</sub> to pH < 2	Cool, 4°C	28 days
Phosphorous, Elemental	G	G	Cool, 4°C	Cool, 4°C	48 hours
Phosphorous, Total	P,G	G	Cool, 4°C H <sub>2</sub> SO <sub>4</sub> to pH < 2	Cool, 4°C	28 days
Silica, Dissolved or Total	P	G	Cool, 4°C	Cool, 4°C	28 days

Section No. H  
Revision No. 0  
Date December 1985  
Page 6 of 10

Table H-1. (Cont'd.)

Parameter	Container <sup>b</sup>		Preservative <sup>c,d</sup>		Maximum Holding Time for all Matrices <sup>e</sup>
	Water	Soil	Water	Soil	
Residue					
Filterable	P	N/A	Cool, 4°C	N/A	7 days
Settleable	P	N/A	Cool, 4°C	N/A	48 hours
Nonfilterable (TSS)	P	N/A	Cool, 4°C	N/A	7 days
Total	P	N/A	Cool, 4°C	N/A	7 days
Volatile	P	N/A	Cool, 4°C	N/A	7 days
Specific Conductance	P	G	Cool, 4°C	Cool, 4°C	28 days
Sulfate	P	G	Cool, 4°C	Cool, 4°C	28 days
Sulfide	P	G	Cool, 4°C Add Zinc Acetate plus NaOH to pH > 9	Cool, 4°C	7 days
Sulfite	P	G	None Required	None Required	Analyze Immediately
Surfactants	P	G	Cool, 4°C	Cool, 4°C	48 hours
Temperature	P	G	None Required	None Required	Analyze Immediately
Turbidity	P	N/A	Cool, 4°C	N/A	48 hours
<u>ORGANIC TESTS<sup>j</sup></u>					
Acrolein and Acrylonitrile	S	S	Cool, 4°C 0.008% Na <sub>2</sub> S <sub>2</sub> O <sub>3</sub> <sup>g</sup> Adjust pH to 4-5 <sup>k</sup>	Cool, 4°C	14 days <sup>k</sup>

Section No. H  
Revision No. 0  
Date December 1985  
Page 7 of 10

Table H-1. (Cont'd.)

Parameter	Container <sup>b</sup>		Preservative <sup>c,d</sup>		Maximum Holding Time for all Matrices <sup>e</sup>
	Water	Soil	Water	Soil	
Benzidines <sup>1</sup>	G	G	Cool, 4°C <sup>m</sup> 0.008% Na <sub>2</sub> S <sub>2</sub> O <sub>3</sub> pH 2-7	Cool, 4°C	7 days until extraction <sup>n</sup>
Chlorinated Hydrocarbons <sup>1</sup>	G	G	Cool, 4°C	Cool, 4°C	7 days until extraction 40 days after extraction
Haloethers <sup>1</sup>	G	G	Cool, 4°C 0.008% Na <sub>2</sub> S <sub>2</sub> O <sub>3</sub>	Cool, 4°C	7 days until extraction 40 days after extraction
Nitroaromatics and Isophorone	G	G	Cool, 4°C Store in Dark	Cool, 4°C Store in Dark	7 days until extraction 40 days after extraction
Nitrosamines <sup>1,0</sup>	G	G	Cool, 4°C Store in Dark 0.008% Na <sub>2</sub> S <sub>2</sub> O <sub>3</sub>	Cool, 4°C Store in Dark	7 days until extraction 40 days after extraction
PCBs	G	G	Cool, 4°C	Cool, 4°C	7 days until extraction 40 days after extraction
Pesticides <sup>1</sup>	G	G	Cool, 4°C pH 5-9 <sup>n</sup>	Cool, 4°C	7 days until extraction 40 days after extraction
Phenols <sup>1</sup>	G	G	Cool, 4°C 0.008% Na <sub>2</sub> S <sub>2</sub> O <sub>3</sub>	Cool, 4°C	7 days until extraction 40 days after extraction
Phthalate Esters <sup>1</sup>	G	G	Cool, 4°C	Cool, 4°C	7 days until extraction 40 days after extraction

Section No. H  
Revision No. 0  
Date December 1985  
Page 8 or 10

Table H-1. (Cont'd.)

Parameter	Container <sup>b</sup>		Preservative <sup>c,d</sup>		Maximum Holding Time for all Matrices <sup>e</sup>
	Water	Soil	Water	Soil	
Polynuclear Aromatic Hydrocarbons	G	G	Cool, 4°C 0.008% Na <sub>2</sub> SO <sub>3</sub> <sup>g</sup> Store in dark	Cool, 4°C Store in Dark	7 days until extraction 40 days after extraction
Purgeable Aromatic Hydrocarbons	S	S	Cool, 4°C 0.008% Na <sub>2</sub> SO <sub>3</sub> <sup>g</sup> HCl to pH < 2 <sup>h</sup>	Cool, 4°C	14 days <sup>q</sup>
Purgeable Halocarbons	S	S	Cool, 4°C 0.008% Na <sub>2</sub> SO <sub>3</sub> <sup>g</sup>	Cool, 4°C	14 days
TCDD <sup>l</sup>	G	G	Cool, 4°C 0.008% Na <sub>2</sub> SO <sub>3</sub> <sup>g</sup>	Cool, 4°C	7 days until extraction 40 days after extraction
Total Organic Carbon	G	G	Cool, 4°C HCl or H <sub>2</sub> SO <sub>4</sub> to pH < 2 <sup>h</sup>	Cool, 4°C	28 days
Total Organic Halogen	G	G	Cool, 4°C 1 ml of 0.1 M sodium sulfite	Cool, 4°C	7 days

Analyses not listed should be preserved at 4°C and held not longer than 7 days.

<sup>a</sup>Preservatives and holding times are from Federal Register, Vol. 49, No. 209, Friday, October 26, 1984, Page 43260 and Characterization of Hazardous Waste Sites: A Methods Manual -- Volume II, Sampling Methods, Second Edition, EPA-600/4-84-076. Container requirements are consistent with these references.

<sup>b</sup>P = Polyethylene

G = Amber Glass with Teflon-lined cap

S = Glass Vial with Teflon-lined septum cap

c<sub>1</sub> Sample preservation should be performed immediately upon sample collection. For composite samples, each aliquot should be preserved at the time of collection. When use of an automatic sampler makes it impossible to preserve each aliquot, samples may be preserved by maintaining at 4°C until compositing and sample splitting is completed.

d<sub>1</sub> When any sample is to be shipped by common carrier or sent through the U.S. Mail, it must comply with the Department of Transportation Hazardous Materials Regulations (49 CFR Part 172). The person offering such material for transportation is responsible for ensuring such compliance. For the preservation requirements in this table, the Office of Hazardous Materials, Materials Transportation Bureau, Department of Transportation, has determined that the Hazardous Materials Regulations do not apply to the following materials: hydrochloric acid (HCl) in water solutions at concentrations of 0.04% by weight or less (pH about 1.96 or greater); Nitric acid (HNO<sub>3</sub>) in water solutions at concentrations of 0.15% by weight or less (pH about 1.62 or greater); Sulfuric acid (H<sub>2</sub>SO<sub>4</sub>) in water solutions at concentrations of 0.3% by weight or less (pH about 1.15 or greater); and Sodium hydroxide (NaOH) in water solutions at concentrations of 0.030% by weight or less (pH about 12.3 or less).

e<sub>1</sub> Samples should be analyzed as soon as possible after collection. The times listed are the maximum times that samples may be held before analysis and still be considered valid.

A Laboratory is obligated to hold the sample for a shorter time if knowledge exists to show this is necessary to maintain sample integrity.

f<sub>1</sub> If samples cannot be filtered within 48 hours, add 1 ml of a 2.71% solution of mercuric chloride to inhibit bacterial growth.

g<sub>1</sub> Should only be used in the presence of residual chlorine.

h<sub>1</sub> Maximum holding time is 24 hours when sulfide is present. Optionally, all samples may be tested with lead acetate paper before pH adjustment in order to determine if sulfide is present. If sulfide is present, it can be removed by addition of cadmium nitrate powder until a negative spot test is obtained. The sample is filtered and then NaOH is added to pH 12.

i<sub>1</sub> For dissolved metals, filter immediately on site before adding preservative.

j Guidance applies to samples to be analyzed by GC, LC, or GC/MS for specific compounds.

k The pH adjustment is not required if acrolein will not be measured. Samples for acrolein receiving no pH adjustment must be analyzed within three days of sampling.

l When the extractable analytes of concern fall within a single chemical category, the specified preservative and maximum holding times must be observed for optimum safeguard of sample integrity. When the analytes of concern fall within two or more chemical categories, the sample may be preserved by cooling to 4°C, reducing residual chlorine with 0.008% sodium thiosulfate, storing in the dark, and adjusting pH to 6-9; samples preserved in this manner may be held for 7 days before extraction and 40 days after extraction. Exceptions to this optimal preservation and holding time procedure are noted in footnotes g, m, and n.

m If 1,2-diphenylhydrazine is likely to be present, adjust the pH of the sample to  $4.0 \pm 0.2$  to prevent rearrangement to benzidine.

n Extracts may be stored up to 7 days before analysis if storage is conducted under an inert (oxidant-free) atmosphere.

o For the analysis of diphenylnitrosamine, add 0.008%  $\text{Na}_2\text{S}_2\text{O}_3$  and adjust pH to 7-10 with NaOH within 24 hours of sampling.

p The pH adjustment may be performed upon receipt at the laboratory and may be omitted if the samples are extracted within 72 hours of collection. For the analysis of aldrin, add 0.008%  $\text{Na}_2\text{S}_2\text{O}_3$ .

q Sample receiving no pH adjustment must be analyzed within 7 days of sampling.

Section No. G  
Revision No. 1  
Date March 1987  
Page 2 of 3

## APPENDIX G

### SAMPLE CONTAINER CLEANING PROCEDURES

To ensure the integrity of aqueous and solid samples, steps must be taken to minimize contamination from the containers in which they are stored. If the analyte(s) to be determined are organic in nature, the container should be made of amber glass. If the analyte(s) are inorganic, the container should be polyethylene. When both organic and inorganic substances are expected to be present, separate samples should be taken. New sample bottles must be cleaned according to the procedure presented below; reuse of sample containers is expressly prohibited. Commercially cleaned containers may be utilized if cleaning procedures comply with those provided in this appendix and prior USATHAMA approval is obtained. The procedures for cleaning the glass and polyethylene containers and their caps are as follows:

- Polyethylene Bottles and Polyethylene Caps
  - (1) Rinse bottles and lids with 5% sodium hydroxide.
  - (2) Rinse with deionized water.
  - (3) Rinse with 5% Ultrex (or equivalent) nitric acid in deionized water.
  - (4) Rinse with deionized water.
  - (5) Drain and air dry.
- Amber-Glass Bottles or 40-ml Vials
  - (1) Scrub and wash bottles in detergent.
  - (2) Rinse with copious amounts of distilled water.
  - (3) Rinse with acetone.
  - (4) Rinse with methylene chloride (Nanograde or equivalent).
  - (5) Rinse with hexane (Nanograde or equivalent).
  - (6) Air dry.
  - (7) Heat to 200°C.
  - (8) Allow to cool.
  - (9) Cap with clean caps with Teflon liners.

Section No. G  
Revision No. 0  
Date December 1985  
Page 3 of 3

- Bottle Caps
  - (1) Remove paper liners from caps.
  - (2) Wash with detergent.
  - (3) Rinse with distilled water.
  - (4) Dry at 40°C.
- Teflon Liners (avoid contact with fingers)
  - (1) Wash with detergent.
  - (2) Rinse with distilled water.
  - (3) Rinse with acetone.
  - (4) Rinse with hexane (Nanograde or equivalent).
  - (5) Air dry.
  - (6) Place liners in cleaned caps.
  - (7) Heat to 40°C for 2 hours.
  - (8) Allow to cool.
  - (9) Use to cap cleaned bottles.



Appendix D

ORGANIC SAMPLE COLLECTION REQUIREMENTS  
AND  
REQUIRED CONTAINERS, PRESERVATION TECHNIQUES, AND HOLDING TIMES



(USER'S GUIDE TO THE CONTRACT LABORATORY PROGRAM, OCTOBER, 1984)

AND

(CONTRACT LABORATORY PROGRAM STATEMENT OF WORK  
FOR INORGANIC ANALYSIS, OCTOBER, 1986 REV.)

# ORGANIC SAMPLE COLLECTION REQUIREMENTS



<u>WATER SAMPLES</u>	<u>REQUIRED VOLUME</u>	<u>CONTAINER TYPE</u>
EXTRACTABLE ANALYSIS (LOW LEVEL)	1 GALLON	2 X 80-OZ. AMBER GLASS BOTTLES OR 4 X 1-LITER AMBER GLASS BOTTLES
EXTRACTABLE ANALYSIS (MEDIUM LEVEL*)	1 GALLON	4 X 32-OZ. WIDE-MOUTH GLASS JARS
VOLATILE ANALYSIS (LOW OR MEDIUM LEVEL*)	80 ML	2 X 40-ML GLASS VIALS

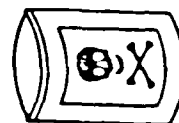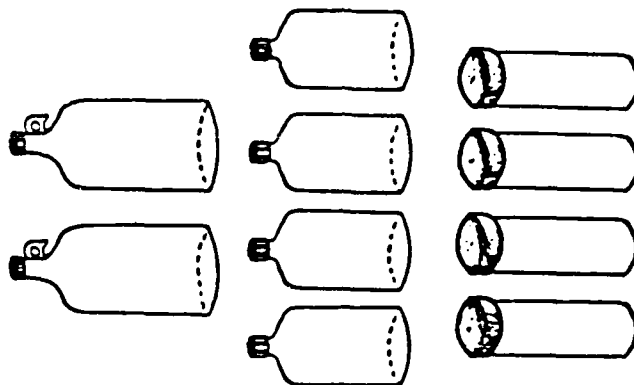

•ALL MEDIUM LEVEL SAMPLES TO BE SEALED IN METAL PAINT CAN FOR SHIPMENT

# ORGANIC SAMPLE COLLECTION REQUIREMENTS

<u>SOIL/SEDIMENT SAMPLES</u>	<u>REQUIRED VOLUME</u>	<u>CONTAINER TYPE</u>
EXTRACTABLE ANALYSIS (LOW OR MEDIUM LEVEL*)	6 OZ.	1 X 8-OZ. WIDE-MOUTH GLASS JAR
		OR
		2 X 4-OZ. WIDE-MOUTH GLASS JARS
VOLATILE ANALYSIS (LOW OR MEDIUM LEVEL*)	240 ML	2 X 120-ML WIDE-MOUTH GLASS VIALS

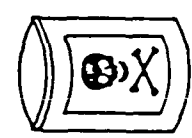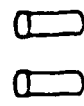

\*ALL MEDIUM LEVEL SAMPLES TO BE SEALED  
IN METAL PAINT CAN FOR SHIPMENT



# INORGANIC SAMPLE COLLECTION REQUIREMENTS

<u>WATER SAMPLES</u>	<u>REQUIRED VOLUME</u>	<u>CONTAINER TYPE</u>
METALS ANALYSIS (LOW LEVEL)	1 LITER	1 X 1-LITER POLYETHYLENE BOTTLE
METALS ANALYSIS (MEDIUM LEVEL*)	16 OZ.	1 X 16-OZ. WIDE-MOUTH GLASS JAR
CYANIDE (CN <sup>-</sup> ) ANALYSIS (LOW LEVEL)	1 LITER	1 X 1-LITER POLYETHYLENE BOTTLE
CYANIDE (CN <sup>-</sup> ) ANALYSIS (MEDIUM LEVEL*)	16 OZ.	1 X 16-OZ. WIDE-MOUTH GLASS JAR

\*ALL MEDIUM LEVEL SAMPLES TO BE SEALED IN METAL PAINT CAN FOR SHIPMENT



# INORGANIC SAMPLE COLLECTION REQUIREMENTS


<u>SOIL/SEDIMENT SAMPLES</u>	<u>REQUIRED VOLUME*</u>	<u>CONTAINER TYPE</u>
METALS AND CYANIDE (CN <sup>-</sup> ) ANALYSIS (LOW OR MEDIUM LEVEL*)	6 OZ.	1 X 8-OZ. WIDE-MOUTH GLASS JAR
		OR
		2 X 4-OZ. WIDE-MOUTH GLASS JARS

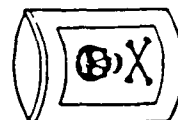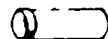

• ALL MEDIUM LEVEL SAMPLES TO BE SEALED  
IN METAL PAINT CAN FOR SHIPMENT

# HIGH HAZARD SAMPLE COLLECTION REQUIREMENTS

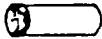
	REQUIRED VOLUME		CONTAINER TYPE
<u>LIQUID SAMPLES</u>			
ORGANIC AND INORGANIC ANALYSIS	6 OZ.		1 X 8-OZ. WIDE-MOUTH GLASS JAR
<u>SOLID SAMPLES</u>			
ORGANIC AND INORGANIC ANALYSIS	6 OZ.		1 X 8-OZ. WIDE-MOUTH GLASS JAR

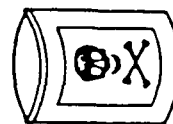

•ALL MEDIUM LEVEL SAMPLES TO BE SEALED  
IN METAL PAINT CAN FOR SHIPMENT

# DIOXIN SAMPLE COLLECTION REQUIREMENTS

<u>SOIL/SEDIMENT SAMPLES</u>	<u>REQUIRED VOLUME</u>	<u>CONTAINER TYPE</u>
2,3,7,8-TCDD (DIOXIN) ANALYSIS	4 OZ.	1 X 4-OZ. WIDE-MOUTH GLASS JAR
		OR
		1 X 8-OZ. WIDE-MOUTH GLASS JAR

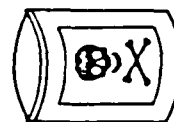

•ALL MEDIUM LEVEL SAMPLES TO BE SEALED  
IN METAL PAINT CAN FOR SHIPMENT

Required Containers, Preservation Techniques, and Holding Times

Measurement Table/Parameter	Container <sup>1</sup>	Preservative <sup>2,3</sup>	Maximum Holding Time For Water Samples <sup>4</sup>
<u>Inorganic Tests</u>			
IB 23-24, Cyanide, total and amenable to chlori nation	P, G	0.6g ascorbic acid <sup>5</sup> NaOH to pH >12 Cool, 4°C	14 days 6
<u>Metals<sup>7</sup></u>			
35, Mercury	P, G	HNO <sub>3</sub> to pH <2	26 days
IB 3, 5-8, 11, Metals, 13, 14, 19, except above 20, 22, 26 29, 30, 32- 34, 36, 37, 45, 47, 51 52, 58, 59, 60, 62, 63, 70-72, 74 75,	P, G	HNO <sub>3</sub> to pH <2	6 months

See following page for notes.



#### Notes

1. Polyethylene (P) or Glass (G).
2. Sample preservation should be performed immediately upon sample collection. For composite samples each aliquot should be preserved at the time of collection. When use of an automated sampler makes it impossible to preserve each aliquot, then samples may be preserved by maintaining at 4°C (+5°C) until compositing the sample splitting is completed.
3. When any sample is to be shipped by common carrier or sent through the United States Mails, it must comply with the Department of Transportation Hazardous Materials Regulations (49 CFR Part 172). The person offering such material for transportation is responsible for ensuring such compliance. For the preservation requirements of Table II, the Office of Hazardous Materials, Materials Transportation Bureau, Department of Transportation has determined that the Hazardous Materials Regulations do not apply to the following materials: Hydrochloric acid (HCL) in water solutions at concentrations of 0.04% by weight or less (pH about 1.96 or greater); Nitric acid (HNO<sub>3</sub>) in water solutions at concentrations of 0.15% by weight or less (pH about 1.62 or greater); Sulfuric acid (H<sub>2</sub>SO<sub>4</sub>) in water solutions at concentrations of 0.35% by weight or less (pH about 1.15 or greater; and Sodium hydroxide (NaOH) in water solution at concentration of 0.080% by weight or less (pH about 12.30 or less).
4. Samples should be analyzed as soon as possible after collection. The times listed are the maximum times that samples may be held before analysis and still considered valid.
5. Should only be used in the presence of residual chlorine.
6. Maximum recommended holding time is less when sulfide is present. Optionally, all samples may be tested with lead acetate paper before the pH adjustment in order to determine if sulfide is present. If sulfide is present, it can be removed by the addition of cadmium nitrate powder until a negative spot test is obtained. The sample is filtered and then NaOH is added to pH 12.
7. Samples should be filtered immediately on-site before adding preservative for dissolved metals.

## DEVELOPMENT OF THE QUALITY ASSURANCE TASK FORCE

Experts from many scientific and regulatory agencies met in Las Vegas, Nevada on February 18-20, 1987 to discuss Quality Assurance/Quality Control. The objectives of the group were to exchange ideas, share resources and technology, explore creative approaches, address key issues, and hopefully develop a unified plan for assuring quality data. A list of the attendees is included in Appendix A. The participants were divided into seven working groups so that the following topics could be discussed in greater detail:

Group 1 - Quality Assurance Management and Data Quality Objectives

Group 2 - On-Site Auditing, Data Review, and Evaluation

Group 3 - Performance Evaluation and Reference Material

Group 4 - Method Validation and Equivalency

Group 5 - Sample Management, Holding Times and Chain of Custody

Group 6 - Statistics and Chemometrics

Group 7 - Documentation and Data Communication

In order to facilitate the continuing exchange of ideas and resources, a proposal was submitted to the group for the formation of a Quality Assurance Task Force. The Quality Assurance Task Force would promote the continual development of a unified approach to QA/QC.

A summary of the findings of each of the working groups is presented below:

### GROUP 1 - QUALITY ASSURANCE MANAGEMENT AND DATA QUALITY OBJECTIVES

Group 1 emphasized that the Army and EPA should define what types of data are needed, note comparisons and differences in data packages, and address methods to meet data requirements. Follow-up meetings on these and other differences should be held. The decision-making personnel in each agency should be identified and be responsible for establishing equivalency. The results should then be communicated to all ten regions.

### GROUP 2 - ON-SITE AUDITING, DATA REVIEW, AND EVALUATION

Group 2 concluded that there was no consistency among the different agencies concerning precertification and certification. EPA's CLP program requires prospective labs to demonstrate, at their

own cost, administrative and technical capabilities before the contract is awarded. USATHAMA awards contracts, through a RFP process, on the basis of a written proposal and the history of the laboratory. Therefore, demonstration of technical proficiency is paid for by USATHAMA. USDA requires a performance evaluation sample for accreditation. Failure on the performance evaluation sample necessitates that the laboratory must wait for six months before reapplication. The EPA CLP considers the postaward performance evaluation samples to be a major topic for the on-site laboratory evaluation, unlike the USATHAMA, the Navy, or the USDA. The group agreed that the on-site evaluation checklist was fairly uniform, but that the frequency of the audits and the level of corrective action applied varied with the agency. The EPA CLP was the only program which looked for serious problems by reconstructing final results from the original raw data during the audit. Some members of the group expressed interest in on-site auditing of the field sampling process. Also, the group proposed that the issues and benefits of agencies sharing laboratory performance information should be addressed.

Group 2 found that even though the time frame allowed for the review of data ranged from one week to three months, the procedures were generally the same. All reviewers looked for outlying data, controls, suspicious calibrations, etc. The group did emphasize that they did not feel that data was over-reviewed. Even though data was reviewed in the same manner, the application of the data depended on the end user and could vary widely. Differences in reviews depended on the auditor's function in the overall project scheme. Any data that was involved in litigation and chain of custody would take longer to review. The group decided that the audit could be facilitated by computerized data scans, organized standard data packages, and easier access to lab personnel and sampling information.

#### GROUP 3 - PERFORMANCE EVALUATION AND REFERENCE MATERIAL

Group 3 began their working session by identifying several issues to discuss. The group felt that the level of effort required to "certify a material" needed to be defined and that increased traceability of materials to NBS would be desirable. Purity documentation of standard materials should include one identity and two purity analyses and only designated labs should be allowed to certify materials. Also reference materials should have a certificate of analysis with dates of preparation and expiration, methods used, and any other pertinent data. Interagency cooperation in making and using batches of reference materials would increase the amount of funding available to make and certify Standard Reference Materials from NBS. The group felt that an interagency work statement was needed to plan a feasibility study on soil reference material preparation and analysis.

Group 3 was then divided into three subgroups: definitions, performance evaluation materials, and reference materials. The definitions subgroup defined reference materials as a general term for characterized substances used for the following three functions: instrument calibration, intralaboratory QC, and interlaboratory performance evaluation.

1. The purpose of calibration is to establish the relationship between instrument response and concentration. Calibration is accomplished by using calibration standards which are well characterized as to purity, stability, and concentration.
2. The purpose of intralaboratory quality control is to provide an assessment of data quality within a laboratory. This information is developed in part by daily analysis of a laboratory check sample. These check samples can be prepared by the laboratory or obtained from an external source. In addition, the laboratory should periodically analyze externally supplied check materials of known composition, such as EPA's QC samples and the SRM's from NBS.
3. The purpose of performance evaluation is to provide an assessment of the comparability of analytical results between laboratories. Performance evaluation materials (PEM's) are periodically supplied to the laboratory as unknowns by the sponsoring agency.

PEM's may be derived from natural matrix materials or synthetically prepared. The evaluation of laboratories using PEM's may be based on comparison with known values or by comparison of individual results against the average performance.

The subgroup on performance evaluation materials discussed the need for solid organic performance evaluation samples, acknowledging that nonvolatile evaluation samples are relatively easy to prepare, but that studies were needed for preparation of volatiles in solids. A need for performance evaluation samples for explosives in soil was expressed. The subgroup decided that natural contaminated matrices were preferable to spiked matrices if possible. The frequency and character of the performance evaluation samples was discussed, as well as the selection of analytes, matrices, and concentrations. The necessity of keeping evaluation samples as blind as possible was recognized, with the minimum frequency being one blind performance evaluation sample per lot of samples. However, studies should be made of laboratory operations to determine the optimum frequency, with a combination of control charts and performance evaluation samples being the best approach. The subgroup decided that the performance evaluation samples should be at multiple levels of concentration

(taking into account the method), and should contain well characterized analytes of low intrinsic variability. Both easy and difficult analytes should be included, as well as those that are easily confused. Performance evaluation results should be purged of proprietary information and shared between the different agencies. Matrices needed to be representative of the sites under study. Selection of six to ten representative types of soil was suggested. The need for a catalog of sources of environmental performance evaluation samples was expressed. The subgroup also discussed the evaluation of the performance of laboratories from evaluation sample data. The use of the results and the formation of data quality objectives should be consistent with the end use of the data such as risk assessment and legal considerations. Education of the end users of the data was suggested because of misconceptions concerning the width of windows for evaluating the results. Participants also reaffirmed that performance evaluation sample results only comprised one portion of the laboratory evaluation and that results would be of marginal value if criteria limits were set near method detection limits. Comments were made that dilution errors were a problem with high level performance evaluation samples and outliers on the high and low sides needed to be handled differently.

#### GROUP 4 - METHOD VALIDATION AND EQUIVALENCY

The first concern of Group 4 was to define method validation in the following manner:

Method Validation is a process starting with definition of method analytes and sample matrices followed by identification of suitable existing methods for conducting the analysis. A method is selected, optimized and validated in a single laboratory study which must include determination of method detection limit, precision and bias in a range of matrices of interest and also include ruggedness testing and writing a complete method protocol. Following a rigorous single lab study, the method undergoes collaborative testing by a minimum of six to eight laboratories. The method is considered to be validated if the written protocol can be followed by the participating laboratories, the method is suitable for the tested matrices, and if the precision and accuracy of the collaborative results are within the limits set in the Data Quality Objectives.

The group also recognized that method validation is a separate process from certification of a laboratory to perform a method. In validating a method, the natural environmental matrix was preferred if the process of locating and characterizing the matrix did not exhaust the available resources. Fortified natural matrices could be used as well as a totally synthetic matrix. Problems associated with obtaining laboratories to participate in interlaboratory studies was discussed. The group felt that competition among the laboratories might make them commit to a validation effort, but there were no guarantees that the data would be delivered in a timely manner.

Guidelines for conducting dynamic validation were developed by the group. It should only be used when more than 20 laboratories are participating, the method must be based on a previously tested method for which there is a high degree of confidence that its performance will meet or exceed program requirements, and the method must include a strong quality control program.

Group 4 recommended that the EPA Superfund staff adopt a policy which would allow other federal agencies to demonstrate equivalency of their methods to Superfund methods. Problems mentioned which could result from adopting such a policy were different reporting requirements and proliferation of equivalent methods which would make data review more difficult.

The group recommended that other federal agencies be invited to attend future caucuses and that the QA workshop should be reconvened in one year or less.

#### GROUP 5 - SAMPLE MANAGEMENT, HOLDING TIMES AND CHAIN OF CUSTODY

Group 5 recognized that a uniform sample definition would be desirable, but may not be possible. However, sample terms must be defined and mutually understood by all agencies involved. QA procedures for field sampling, field logs, and chain-of-custody documentation should be uniform. Total compliance in maintaining chain-of-custody could be very difficult with more automated analyses. More stringent QA concerning field sampling is needed since this represents a huge source of error. Some estimate of this error would be desirable. The group felt that information on the validity of holding times would be desirable in unifying specifications among the various agencies. Requirements for reanalysis were different among the agencies if checking of data revealed that analyses were inappropriately performed. A centralized database as a means of obtaining summary information on laboratories such as current standing, date of most recent audit, and date of most recent performance evaluation sample analysis would be advantageous to all agencies concerned.

#### GROUP 6 - STATISTICS AND CHEMOMETRICS

Group 6 raised five issues for discussion. The first issue was detection limits. The group established that the detection limit represents a concentration where decisions about presence or absence are made and that the quantitation (reporting) limit is at some concentration above the detection limit. The detection limit is highly dependent on the individual characteristics of the apparatus, analyst, method, etc. Data was shown to suggest that EPA's MDL and THAMA's Hubaux and Vos estimates from some data sets show a maximum difference

in ratio of 1.5. The group concluded that the two procedures were not as different as thought at first, and expressed a need for more information on how to set limits for multi-analyte methods, using surrogates to evaluate matrix variations, more comparative evaluations of the different methods of estimating detection limits, more education on the variables which are included in the detection limit estimates, and a determination of the most effective way to specify concentration limits and evaluate inherent risks.

The second issue raised was chemometrics. During this discussion, the following needs were expressed: investigation of applications of composite sampling techniques to environmental monitoring, estimation of variability due to sampling, improved laboratory subsampling protocols, and more nested experimental programs to provide objective estimates based on real samples of the following: field sampling variability, lab subsampling/preparation, and analytical variability.

The third issue that the group discussed was the development of performance evaluation sample criteria. Double-blind performance evaluation sample submission was considered ideal when feasible at a frequency consistent with the needs of the program and cost benefit goals. Ultimately the group wanted to see capability limits established for the performance evaluation standards for various methods and for different types of evaluation materials. In establishing these criteria, the problem of editing data to exclude true "outliers" without unduly truncating data sets was recognized. The question was also raised about the effect of a large number of outliers in a data set upon future repeatability. A suggestion was made to use the Biweight Robust Estimation Procedure (JASA, June 1982) which provides the basis for USEPA performance evaluation criteria limits for water analysis. Some out-of-control data could be discarded if the frequency of performance evaluation checks is increased.

The group labeled the fourth issue as design/analysis comparability. Essentially, improved communication between lab personnel and field samplers was a necessity in order to carefully formulate all sampling protocols and analysis in advance so that all statistical computations would be compatible with the actual performance of the experiments. Any uncontrolled variables in the procedures also needed to be noted to aid in the design process as well as final use of the data (qualitative vs. quantitative).

Concerning the last issue, validation of reference materials, the group required detailed procedures to establish homogeneity and stability.

## GROUP 7 - DOCUMENTATION AND DATA COMMUNICATION

Group 7 agreed that the quality and quantity of the documentation required for a program varies depending on the end use of the data from the most intensive documentation necessary for litigation purposes to minimum documentation only needed for characterization/information gathering to make rapid decisions. The group stated that sampling and analytical contracts should require the level of documentation necessary for the program's purposes with the long term goal that the various agencies would reach some consensus on what documentation should be required. Meeting this goal would solve such problems as agreement of EPA and DOD on documentation when both agencies are involved in sites on DOD facilities and alleviate the frustration of contractor labs who are required to follow different documentation and reporting requirements depending on the agency. The group questioned whether software systems could be developed and implemented to process some of the quality control elements common to all the agencies and provide some documentation. USATHAMA followed up this discussion with a presentation of their computer software.



## DISTRIBUTION LIST

### Internal

1. Central Research Library, 4500-N, MS 286
2. Document Reference Section, 9711-1
- 3-12. M. R. Guerin, 4500-S, MS 120
13. Laboratory Records, 4500-N, MS 285
14. Mr. J. A. Lenhard, DOE/ORO
- 15-16. S. K. Holladay, 4500-S, MS 120
- 17-26. M. P. Maskarinec, 4500-S, MS 120
27. M. Miller, K-1004-C, MS 440
28. ORNL Patent Office, 4500-N, MS 258

### External

- 29-40. Defense Technical Information Center, ATTN: DTIC-DDA, Cameron Station, Alexandria, VA 22314.
- 41-42. Commander, U.S. Army Toxic and Hazardous Materials Agency, ATTN: AMXTH-CO-P, Aberdeen Proving Ground, MD 21010-5401.
- 43-44. Defense Logistics Studies Information Exchange, U.S. Army Logistics Management Center, Fort Lee, VA 23801.
- 45-49. Mary Ann Ryan, USATHAMA, AMXTH-TE-A, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
50. Nils Akerlind, U.S. Air Force, HQ AFESC/RDN, Tyndall AFB, FL 32403.
51. Andy Anderson, USATHAMA, AMXTH-IR, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
52. Stanley Blacker, USEPA/ORD, 401 M St., S.W., Washington, DC 20460.
53. Dave Bottrell, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
54. George Brilis, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
55. Paul Britton, USEPA EMSL/CIN, 26 W. St. Clair St., Cincinnati, OH 45268.
56. Ken Brown, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
57. Judy Burris, USAF, USAF OEHL/TS, Brooks Air Force Base, TX 78235.
58. Larry Butler, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
59. Joe Campana, Environmental Research Lab, UNLV, 4505 S Md. Parkway, Las Vegas, NV 89154.
60. Bob Clerman, MITRE Corp., 7525 Colshire Dr., McLean, VA 22102.
61. Wallace Eakes, NEESA, Code 112I, Port Hueneme, CA 93043.
62. Gene Easterly, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
63. Timothy Fisher, US Army Environmental Hygiene Agency, Aberdeen Proving Ground, MD 21010.

Distribution (Cont'd)

64. Jerry Fitzgerald, MITRE Corp., 7525 Colshire Dr., McLean, VA 22102.
65. George Flatman, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
66. Forest Garner, Lockheed-Las Vegas, P.O. Box 15027, Las Vegas, NV 89114.
67. Duane Geuder, USEPA/OSWER, WH548A, 401 M St., Washington, DC 20460.
68. C. L. Grant, University of New Hampshire, Parsons Hall, Chemistry Department, Durham, NH 03824.
69. Fred Haeberer, EPA/QAMS, WH548A, 401 M St., Washington, DC 20460.
70. Michael K. Hoffman, USDA/FSIS, 300 12th St. SW, Washington, DC 20250.
71. Mike Homsher, Lockheed EMSCO, P.O. Box 15027, Las Vegas, NV 89114.
72. Chuck Hoover, Lockheed EMSCO, P.O. Box 15027, Las Vegas, NV 89114.
73. William Horwitz, FDA, Center for Food, Safety and Applied Nutrition, HFF-7, Washington, DC 20204.
74. Peter Isaacson, Sample Management Office, P.O. Box 918, Alexandria, VA 22313.
75. Tom Jenkins, USACRREL, ATTN: CRREL-RC, Hanover, NH 03775.
76. Captain Philip Jung, USAF, OEHL/TF, Brooks Air Force Base, TX 78235-5000.
77. D. A. Kane, USAF, HQ USAF/LEEV, Bolling Air Force Base, Washington, DC 20332.
78. Ed Kantor, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
79. Suji Kumar, Lockheed EMSCO, P.O. Box 15027, Las Vegas, NV 89114.
80. Ken Lang, USATHAMA, AMXTH-TE-A, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
81. Kathy Laukner, UNLV-ERC, Las Vegas, NV 89114.
82. T. J. Lewis, USAF, HW USAF/SGPA Bolling Air Force Base, Washington, DC 20332.
83. Jim Longbottom, USEPA EMSL/CIN, 26 W. St. Clair St., Cincinnati, OH 45268.
84. Paul Marsden, S-Cubed, P.O. Box 1620, La Jolla, CA 92038.
85. Bill Maurits, Product Assurance Direct, AMSMC-QAO-C(A), Aberdeen Proving Ground, MD 21010.
86. Willie May, National Bureau of Standards, Gaithersburg, MD 20899.
87. Ron Mitchum, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
88. John Moore, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
89. J. Gareth Pearson, USEPA/LV, P.O. Box 15027, Las Vegas, NV 89114.

Distribution (Cont'd)

90. Florence Richardson, USEPA/OSWER, WH562B, 401 M St., Washington, DC 20460.
91. Pete Rissell, USATHAMA, AMXTH-TE-A, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
92. David Rozak, USAEHA, Bldg. E2100, Aberdeen Proving Ground, MD 21010-5401.
93. Amy Smiecinski, Environmental Research Lab, UNLV 4505 S. Md. Parkway, Las Vegas, NV 89154.
94. Robin Stein, USATHAMA, AMXTH-TE-A, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
95. Bud Sturtzer, NEESA, Code 112N, Port Hueneme, CA 93043.
96. Marty Stutz, USATHAMA, AMXTH-TE-A, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
97. J. Scott Warner, Battelle Columbus Division, 505 King Ave., Columbus, OH 43201.
98. C. Weber, USEPA EMSL/CIN, 26. W. St. Clair St., Cincinnati, OH 25268.
99. L. R. Williams, USEPA EMSL/LV, P.O. Box 15027, Las Vegas, NV 89114.
100. John Winter, USEPA EMSL/CIN, 25 W. St. Clair ST., Cincinnati, OH 25268.
101. Denny Wynne, USATHAMA, AMXTH-TE-A, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
102. Bob York, USATHAMA, AMXTH-TE-A, Bldg. E-4585, Aberdeen Proving Ground, MD 21010-5401.
- 103-207. Technical Information Center, U.S. Department of Energy, Oak Ridge, TN 37831. For DOE/TIC 4500 distribution under UC-4 category.

END

DATE

FILMD

3-88

DTIC